

MEDIT Link



Revision 4 (2021.10)

CONTENTS

1 Introduction and Overview	5
1.1 Cloud Service	6
1.2 Workflow Management Using Ordering System.....	6
1.3 Utilizing the Dashboard	6
1.4 Installation	7
1.4.1 System Requirements.....	7
1.4.2 Installation Guide	7
2 Registration	10
2.1 Landing Page.....	10
2.2 Registration.....	11
2.3 Login and Welcome Page	12
3 Medit Link Web.....	13
3.1 Dashboard	14
3.2 Case Box/Work Box	17
3.3 Case Talk	19
3.3.1 Creating a Case Talk.....	19
3.3.2 Web Viewer	24
3.3.3 Data Tree	27
3.4 Partners	29
3.4.1 Partner Search	29
3.4.2 Recommended Partners.....	29
3.4.3 My Partner List	30
3.4.4 Partner Invitation	31
3.5 Membership	32
3.5.1 Membership	33
3.5.2 Billing Info.....	40
3.5.3 Redeem.....	45
3.6 Payment.....	47
3.6.1 Payment Details (Clinics, Labs accounts).....	47
3.6.2 Payment account information (Lab account only)	47
3.7 My Information.....	48
3.7.1 My Account	48
3.7.2 My Organization (Admin account only)	48
3.7.3 My Price List (Lab: Admin account only)	49

3.8 Notice and Updates	50
3.9 Help Center	50
3.10 Download the App	50
4 Medit Link App Overview	51
4.1 Title Bar	52
4.1.1 Network Status	52
4.1.2 Case Synchronization Manager	53
4.1.3 User Information	55
4.1.4 Notifications	56
4.2 Boxes Overview	57
4.3 Settings	59
4.3.1 Patient Management	68
4.3.2 Notices and Updates	68
4.3.3 My Account	68
4.3.4 Help Center	68
4.4 Detailed Functions of Case Box/Work Box	69
4.4.1 Creating cases	69
4.4.2 Change viewing options (Clinic Account only)	69
5 Medit Link App – Box Detail.....	71
5.1 Case Box / Work Box.....	71
5.2 Order Box / In Box	73
5.3 Case Talk	74
5.4 App Box.....	76
5.5 Trash Box	77
6 Detailed Explanation of the Process	78
6.1 Managing Cases: Creating, Editing, Cloning, Deleting	78
6.1.1 Creating Cases	78
6.1.2 Editing cases	80
6.1.3 Cloning cases	82
6.1.4 Deleting cases (Lab and Clinic)	83
6.2 Inputting form information	84
6.2.1 Inputting Form Information in Case Box	84
6.2.2 Inputting Form Information in Work Box	85
6.2.3 Inputting Bridge Information	86
6.2.4 Setting Scanning Options.....	87
6.3 Copying Form Information	88

6.4 Deleting treatment information	89
6.5 Adding Product Information to “Favorite”	90
6.5.1 Adding to Favorites	90
6.5.2 Managing Favorite Templates	90
6.6 Copying form information between Scan and CAD tabs	91
6.6.1 From Scan tab to CAD tab	91
6.6.2 From CAD tab to Scan tab	92
6.7 Scanning process	93
6.8 CAD Process	93
6.9 Exporting Result Files.....	94
6.10 File Viewer	95
6.10.1 Data Tree	95
6.10.2 Model Control	96
6.10.3 Capture and Memo	97
6.10.4 Orthodontic Capture	97
6.10.5 Fill Holes	98
6.10.6 Trimming	99
6.10.7 Manage the Multi Occlusion Data in File Viewer	102
6.11 Order Process and Order Management	105
6.11.1 Ordering	105
6.11.2 Job Completed – Work Box	108
6.11.3 Checking Order information.....	111
6.11.4 Managing Expected/Requested Delivery Dates	112
6.12 Patient Management.....	116
6.12.1 Adding Patients	116
6.12.2 Deleting Patients	117
6.12.3 Merging Patients	118
6.13 Case filtering.....	119
6.13.1 Filtering by Status.....	119
6.13.2 Filtering Using Period Settings	120
6.13.3 Grouping Cases.....	121
6.13.4 Filtering with Search Capabilities	122

1 Introduction and Overview

Medit Link is a collaborative tool designed to facilitate communication between clinics and labs to support patient-centered workflows. It includes tools according to your business type (Clinic or Lab) and specific features can be used according to the account type. Medit Link helps users manage a series of digital dentistry workflows including scanning with Medit's intraoral and table top scanners, prosthesis order and manufacturing process for prosthesis through a user-friendly system.

Medit Link consists of a web-service and a computer application.

Medit Link Overview

		Web	App
		Register, search and manage partners, share cases via Web Viewer with a sharable link; access cases made through the App.	Scan, send/receive orders to/from the partners found through Web, enter form information; do CAD and manage patients.
Clinic	Admin	Search and manage partners, as well as your organization. ● Dashboard ● Case Box ● Case Talk ● Partners ● Membership ● Payment ● My Info	Scan, do CAD, send orders. ● Dashboard ● Case Box ● Order Box ● Case Talk ● App Box ● Trash Box ● Patient Management
	Member	● Case Box ● Case Talk ● My Info	● Case Box ● Case Talk ● App Box ● Trash Box
Lab	Admin	Search and manage partners, as well as your organization. ● Dashboard ● Work Box ● Case Talk ● Partners ● Membership ● Payment ● My Info	Scan, do CAD, receive orders. ● Dashboard ● In Box ● Work Box ● Case Talk ● App Box ● Trash Box
	Member	● Work Box ● Case Talk ● My Info	● Work Box ● Case Talk ● App Box ● Trash Box

1.1 Cloud Service

Medit Link supports automatic upload/download of database and files through AWS (Amazon Web Services). The data is automatically synchronized with computers of multiple users in the organization so that they can access this data from any computer at any time. Being connected to the Cloud helps with data backup.

The Cloud is also useful when it comes to quickly and easily sharing data between the dentist and the labs, creating a working environment that minimizes the hassle of sharing data using email. All data patient-related information is securely encrypted, stored, and transmitted, and complies with medical law, HIPAA regulations, and GDPR regulations.

1.2 Workflow Management Using Ordering System

Medit Link provides order and payment system.


You can quickly and easily share your data by signing a contract between the clinic and the lab.

You can view/manage the status of your orders.

Credit card-backed system works in conjunction with online payment platform companies Stripe (www.stripe.com) and I'mPort (<http://www.iampart.kr>) payment modules, making payment services easy and convenient.

1.3 Utilizing the Dashboard

You can view the statistics for the entire work done by Medit Link in a visualized graphs and chart – both in App and Web versions of Medit Link.

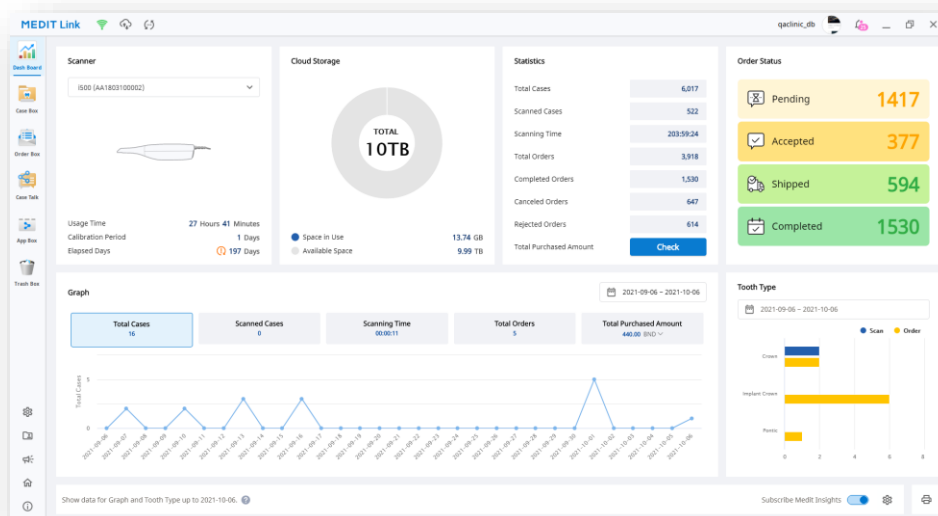

1.4 Installation

1.4.1 System Requirements

	Laptop	Desktop
CPU	Intel Core i5 2.6 GHz or higher	
RAM	16 GB or higher	
Graphics	NVIDIA GeForce GTX 760 (2 GB) or higher / or equivalent AMD video card	
OS	Windows 8 64 Bit (unavailable in 32 Bit) or higher	


1.4.2 Installation Guide

Medit Link is installed as a package with Medit Scan for Clinics and Medit Scan for Labs.


- ① Run Medit_LinkSetup.exe.
- ② Select the language for installation.




- ③ Select the folder for installation of Medit Link. Read “License terms & conditions” carefully and check the box “I agree to the License terms & Conditions.” Click “Install” to proceed to the installation.

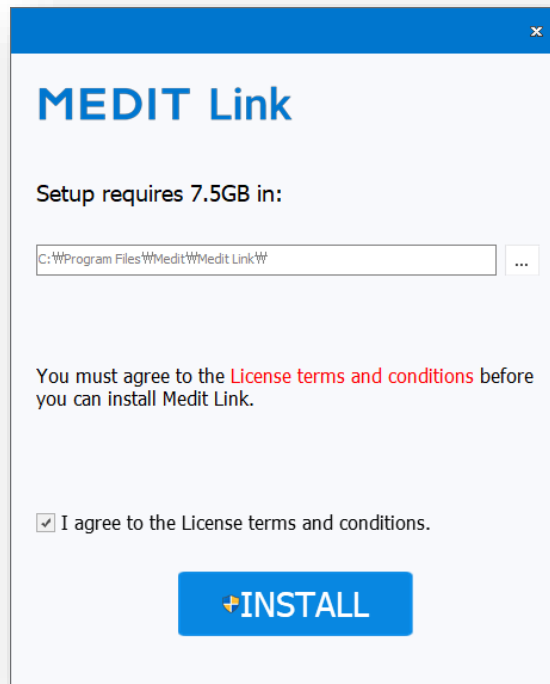

- ④ As Medit Link is installed with Medit Scan Application, the installation will not proceed if the scanner is connected to the computer.
- ⑤ Make sure that the scanner cable is disconnected from the computer before proceeding to the next step.

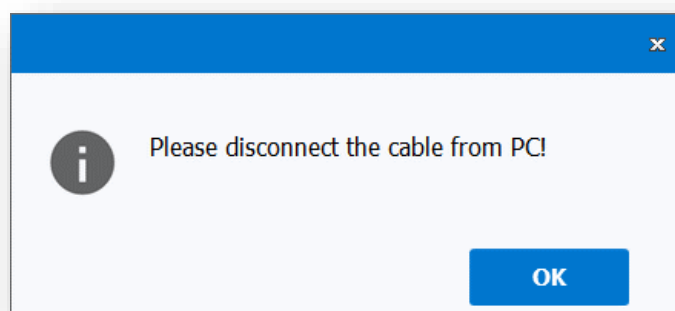

- ⑥ From this step, the installation will continue automatically and may take a few minutes. Do not turn off the computer until the installation is complete.

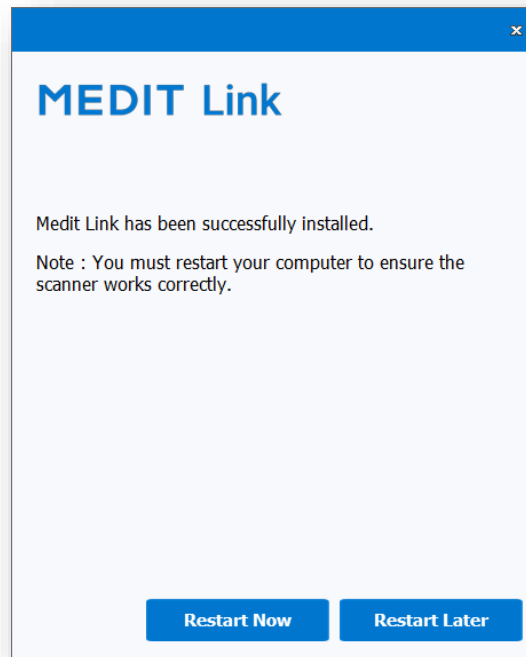

- ⑦ After the installation is complete, restart the computer to use Medit Scan Application with Medit Link. You can skip this step If you do not have a scanner.




The scanners may not work without restarting the computer.

2 Registration

2.1 Landing Page

Register at www.meditlink.com homepage.

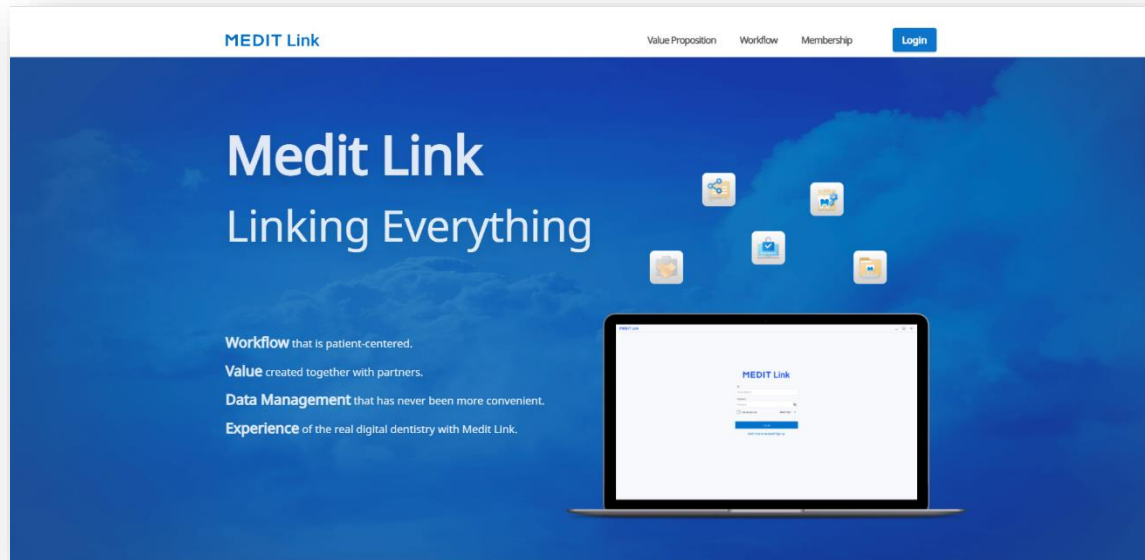

2.2 Registration

The registration interface is divided into two main columns: **Clinic** and **Lab**. Each column contains two registration options: **Administrator** and **Member**. Each option includes a list of available functions and a **Register** button.

Category	Role	Available Functions
Clinic	Administrator Dentist (Manager)	✓ Order / Payment ✓ Member Management
	Member Dental Hygienist and Dental Staff	✓ Oral Scanning ✓ Ordering of Dental Prostheses
Lab	Administrator Representative / Leader (Manager)	✓ Price Management for Dental Prostheses ✓ Order Management ✓ Member Management
	Member Laboratory Employees	✓ Receiving Data for Making Dental Prostheses ✓ Sharing Work Files and Production Processes

Medit Link requires an initial registration as Clinic (administrator & member) or Lab (administrator & member).


Administrator is usually a representative of your organization and has unlimited access to all the services, whereas members have a limited access to Medit Link services.

- ① After specifying the name of your clinic or lab, click “Check Availability.”
 - If you are registering as a member, your organization should be activated prior to your registration.
- ② Once you complete all the fields and click “Register,” a verification email will be sent to you.
- ③ Verify your email address through the link in the email to complete the registration process.

2.3 Login and Welcome Page

On the login page, enter your email address and password and click Login.

- ① To use Medit Link services, information about your organization should be entered using the administrator account.
- ② The information entered on the **Welcome Page** (My Organization etc.) can be edited anytime on “My Organization” page.
- ③ You can also register and manage intraoral and table top scanners and milling machines.
- ④ If you agreed to receive marketing information during registration, you will receive useful information for the product and service from the manufacturer.

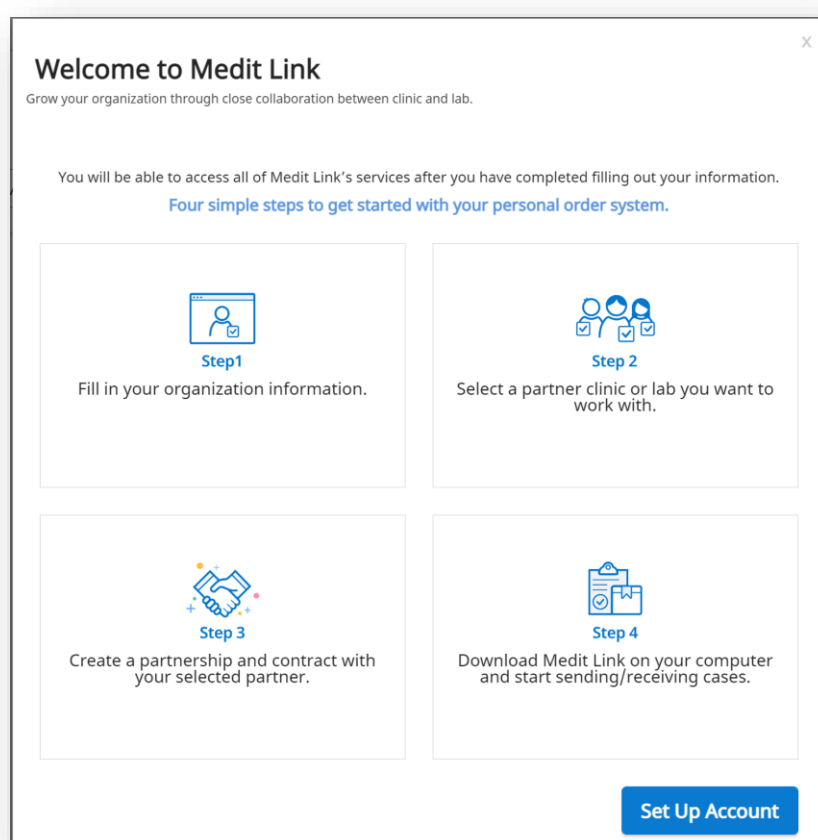


3 Medit Link Web

The screenshot shows the 'Case Box' interface of the Medit Link web application. The top navigation bar includes a sidebar with icons for Dashboard, Case Box, Case Talk, Partners, Membership, Payment, and My Info. The main header area contains filters for status (All, Form Only, Need Processing, Scan Completed, CAD, Ordered, Case Completed), a date range (2018-01-01 - 2021-10-07), and a search bar for Case Name. The table below lists individual cases with columns for Status, Case Name, Patient Name, Form Information, Last Modified Date, Scan Date, Order Date, and Lab Name. Each row includes a 'Share' button.

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scan Date	Order Date	Lab Name
Scan Completed	BS_4_Scan for Clinics's Case	BS_4_Scan for Clinics	Maxilla, Mandible	2021-10-06 14:07	2020-09-25 14:12	-	-
Ordered	Clinical Sample 19	QA	44, 45	2021-10-01 16:49	2021-06-11 10:53	2021-10-01 16:48	qalab_db
Scan Completed	Clinical Sample 19	QA	44, 45	2021-10-01 16:36	-	-	-
Scan Completed	Clinical Sample 19	QA	44, 45	2021-10-01 16:35	-	-	-
Scan Completed	Clinical Sample 19	QA	44, 45	2021-10-01 15:50	-	-	-
Scan Completed	Clinical Sample 19	QA	44, 45	2021-10-01 15:47	-	-	-
Ordered	임플란트 46, 47's Case	임플란트 46, 47	46, 47	2021-09-16 16:37	2020-10-20 11:53	2021-09-16 16:07	qalab_db
Case Completed	임플란트 46, 47's Case - Clone	임플란트 46, 47	46, 47	2021-09-16 16:29	2020-10-20 11:53	2021-09-16 16:21	qalab_db
Ordered	재강 확인 크라운 브릿지 44-46's Case	크라운 브릿지 44-46	44, 46, 45	2021-09-16 16:06	2020-10-20 11:46	2021-09-16 16:06	qalab_db
Ordered	임플란트 46, 47's Case	임플란트 46, 47	46, 47	2021-09-16 15:54	2020-10-20 11:53	2021-09-16 15:54	qalab_db

3.1 Dashboard

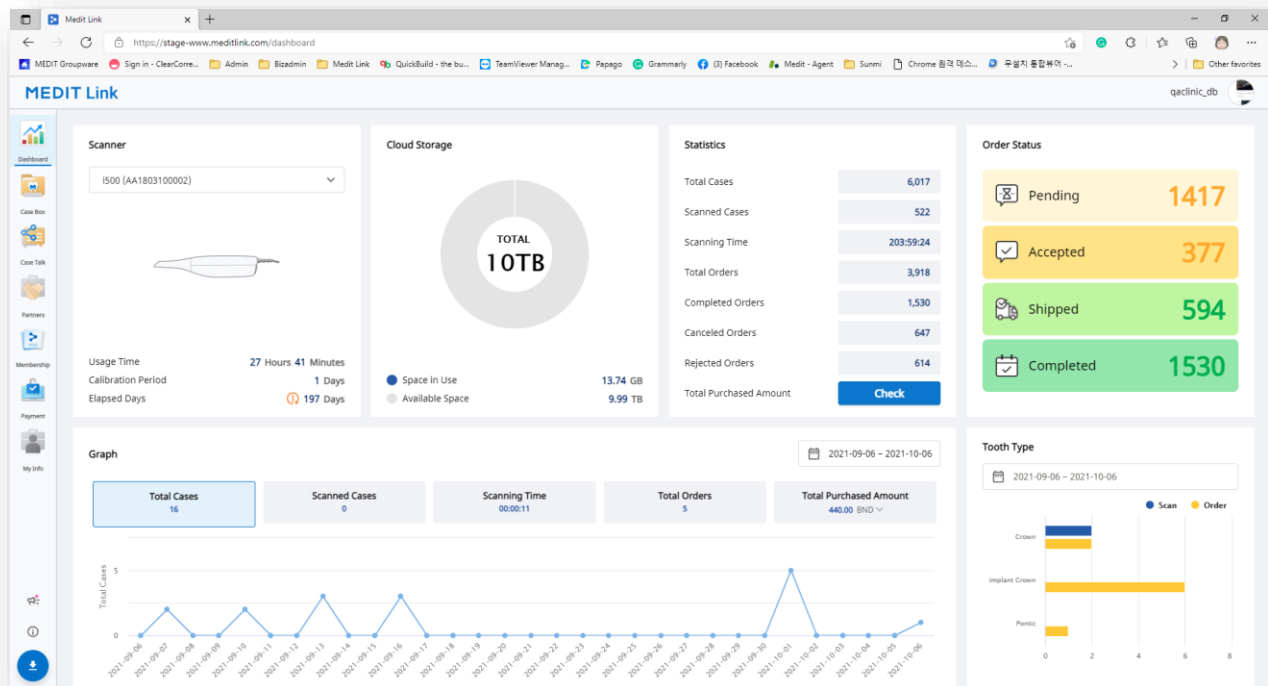

Graphs and charts show statistics for the entire work done in Medit Link.


Scanner	Provides information about the active scanners, including scanner usage time, calibration interval, and calibration elapsed time.
Cloud Storage	Displays the total capacity of the cloud and the amount in use.
Statistics	Provides statistics on cases based on job status.
Order Status	Provides statistics on cases based on order status.
Graph	Graph job information by period.
Tooth Type	Provides statistics on the number of scans and orders by treatment information.

Medit Insights

Medit provides the reporting service and helps you analyze your usage statistics. You can receive a report email on the first day of every month. No additional sign-up process is required. Just enable “Subscribe Medit Insights.”



Medit Insights service is available only for the Clinic admin account.

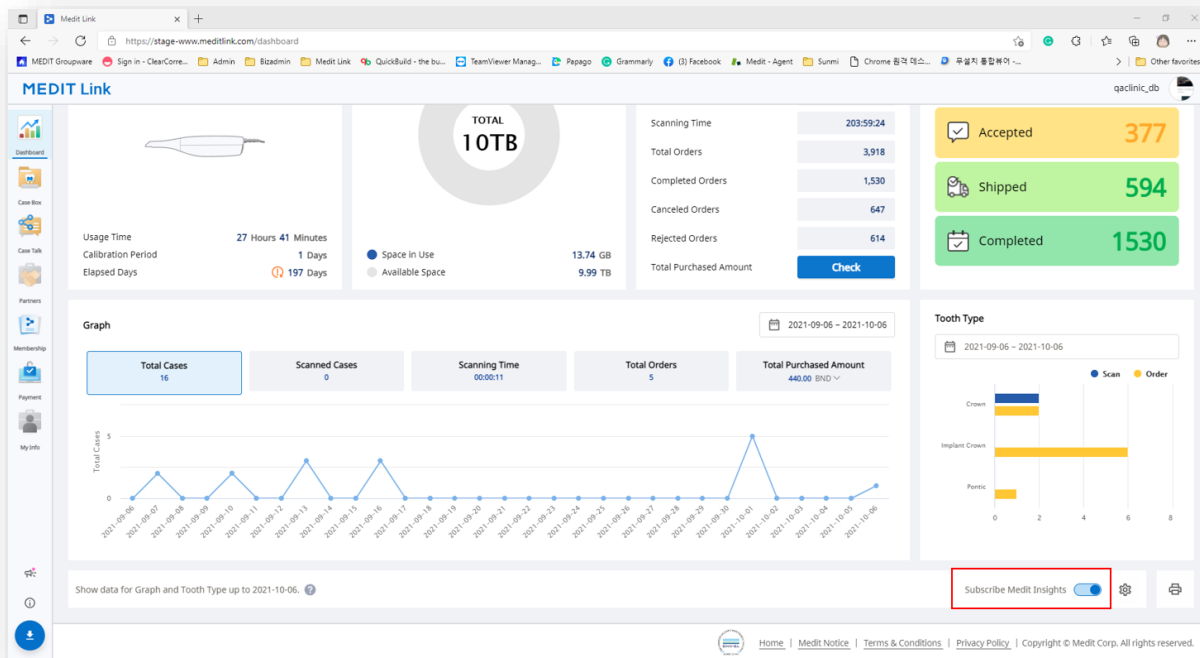

If you don't use the email address you used for the Medit Link account, you can set another email address to receive Medit Insights.

Settings

Receive the report via

☒ Medit Link account


☐ Another email address

Please enter the email address.

The report will be sent on the 1st day of every month.

Cancel OK

You will receive the statistics email like below on the first day of every month.

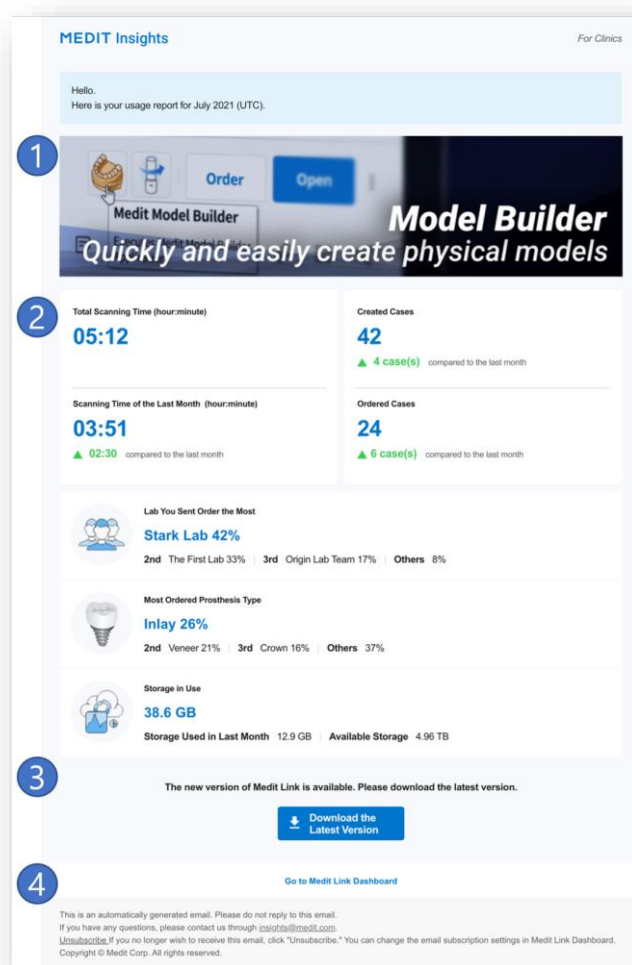

① Banner

Medit provides useful tips and information here. If you don't want to see the banner, uncheck the "Marketing and Promotion Agreement" check box in My Info > My Account.

② Statistics Information

The statistics information is displayed, such as scanning time, created cases, most ordered lab, etc.

③ Medit Link version

If you are not using the latest Medit Link version, we provide the download link for the latest version.

④ Dashboard


All statistics information is based on Dashboard. You can see more details in Dashboard.

How to unsubscribe the email

- ① Disable “Subscribe Medit Insights” in Dashboard.
- ② Click “Unsubscribe” in the received email.

3.2 Case Box/Work Box

Shows all cases which the user has created on the Medit Link application. See files via a web browser without installing a separate program.

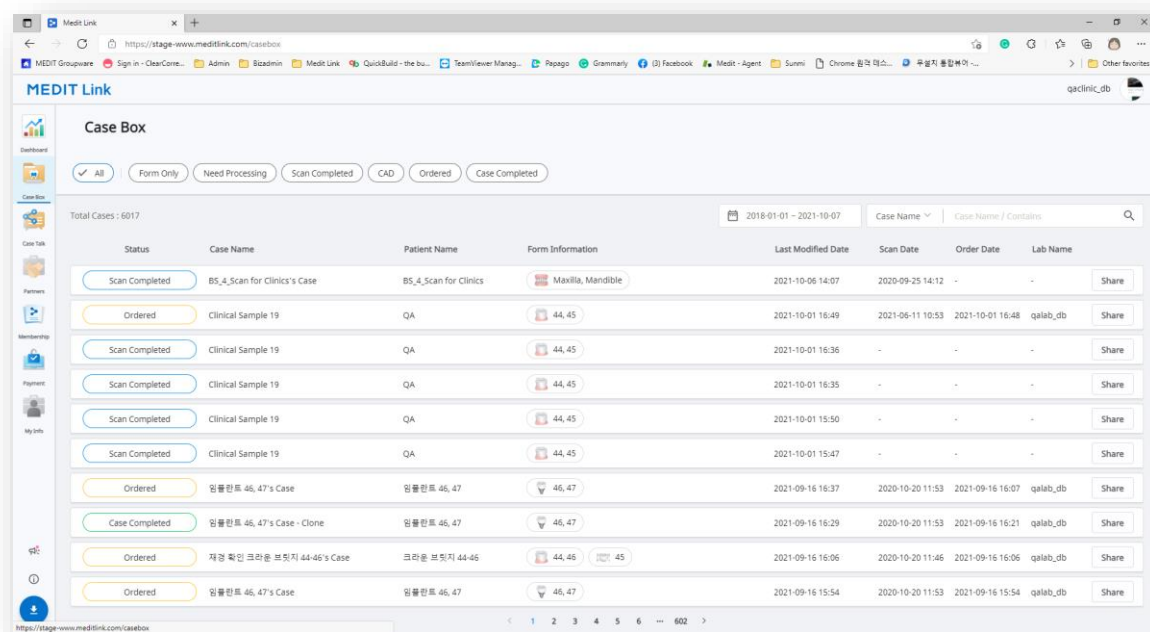


The screenshot shows the 'Case Box' interface in the Medit Link application. At the top, there are filter buttons: 'All', 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'. Below the filters, it says 'Total Cases : 6017'. There is a date range selector set to '2018-01-01 - 2021-10-07' and a search bar. The main table displays the following columns: Status, Case Name, Patient Name, Form Information, Last Modified Date, Scan Date, Order Date, and Lab Name. The table contains several rows of case data, including 'B5_A_Scan for Clinics's Case', 'Clinical Sample 19', and '임플란트 46, 47's Case'.

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scan Date	Order Date	Lab Name
Scan Completed	B5_A_Scan for Clinics's Case	B5_A_Scan for Clinics	Maxilla, Mandible	2021-10-06 14:07	2020-09-25 14:12	-	-
Ordered	Clinical Sample 19	QA	44, 45	2021-10-01 16:49	2021-06-11 10:53	2021-10-01 16:48	qalab_db
Scan Completed	Clinical Sample 19	QA	44, 45	2021-10-01 16:36	-	-	-
Scan Completed	Clinical Sample 19	QA	44, 45	2021-10-01 16:35	-	-	-
Scan Completed	Clinical Sample 19	QA	44, 45	2021-10-01 15:50	-	-	-
Scan Completed	Clinical Sample 19	QA	44, 45	2021-10-01 15:47	-	-	-
Ordered	임플란트 46, 47's Case	임플란트 46, 47	46, 47	2021-09-16 16:37	2020-10-20 11:53	2021-09-16 16:07	qalab_db
Case Completed	임플란트 46, 47's Case - Clone	임플란트 46, 47	46, 47	2021-09-16 16:29	2020-10-20 11:53	2021-09-16 16:21	qalab_db
Ordered	재강 확인 크라운 브릿지 44-46's Case	크라운 브릿지 44-46	44, 45	2021-09-16 16:06	2020-10-20 11:46	2021-09-16 16:06	qalab_db
Ordered	임플란트 46, 47's Case	임플란트 46, 47	46, 47	2021-09-16 15:54	2020-10-20 11:53	2021-09-16 15:54	qalab_db

Filtering by status

Provides the ability to filter cases according to status.

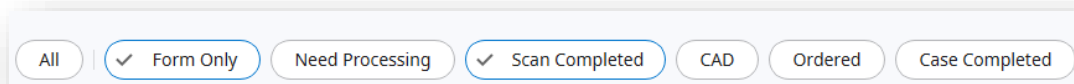

- ① The user can click the status they want to view. The browser will display only cases in the chosen status.
- ② The user can select multiple statuses at once.


Case Box – Clinic	All/Form Only/Need Processing/Scan Completed/CAD/Ordered/Case Completed
Work Box – Lab	All/Form Only/Need Processing/Scan Completed/CAD/Done

Filtering using period settings

Filters and displays cases within the period of time you set.

① Click  2020-08-10 ~ 2020-09-09 .

② Select the range.

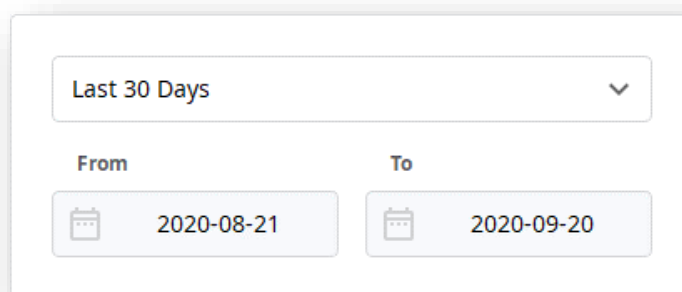


- All: Shows all cases.
- Custom Range: Specifies start and end dates to display cases within this period.
- Today: Shows cases for today's date.
- Yesterday: Shows cases with yesterday's date.
- Past 7 days: Shows cases before 7 days except today.
- Past 30 days: Shows cases before 30 days except today.
- This Month: Shows this month's cases.
- Previous Month: Shows previous month's cases.

Case corresponding with the date range chosen will be updated in real time.

Filtering with search capabilities

Leverage search capabilities in order to quickly display the information.

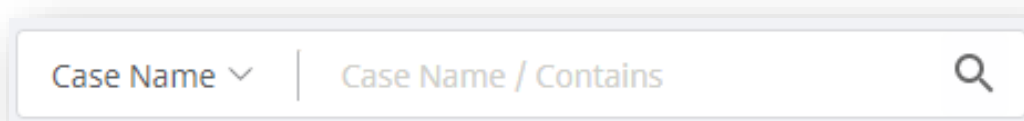


Each box provides the following search options:

Case Box – Clinic	Case Name, Patient Name, Patient ID, Lab Name
Work Box – Lab	Case Name, Patient Name, Clinic Name, Order ID

3.3 Case Talk

Cases created in the Medit Link can be shared with restrictions, without restrictions, with specific people, and registered members of Medit Link. The user can see all shared lists through “Case Talk” in the Medit Link application and Medit Link web service.

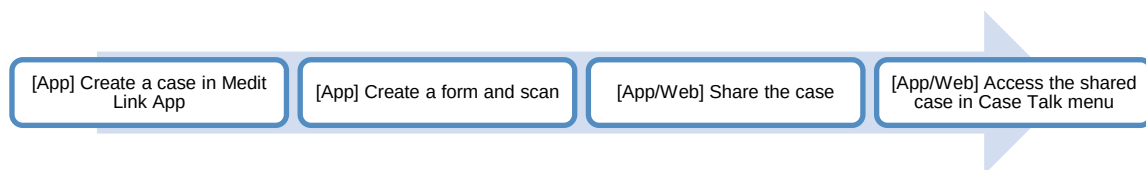


Shared cases can be viewed via a web browser without installing a separate program, enabling easy communication over the 3D data on various devices.


3.3.1 Creating a Case Talk

Create a shared link from a case created in Case Box/Work Box.

- ① Select the case you want to share.

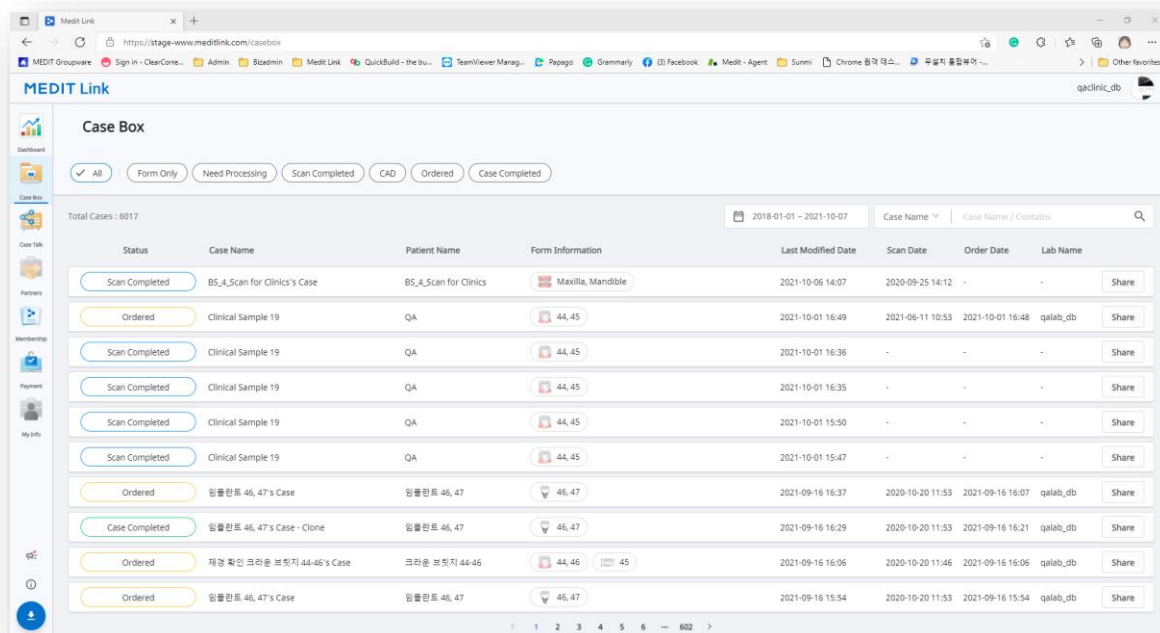

- ② Click the “Share” button.
- ③ Select the files you want to share.

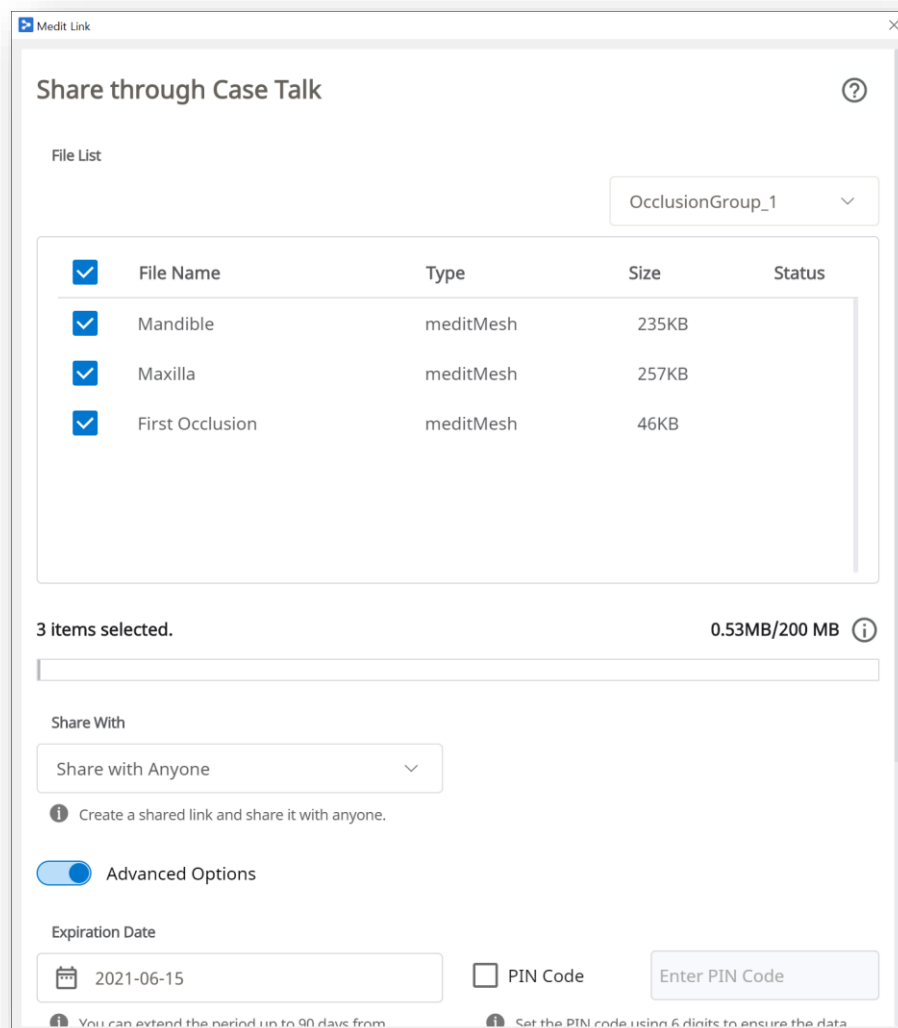

- ④ Set share scope and access rights.
- ⑤ Click “Get Link” to create a shared link.

Set the share scope and access rights

Medit Link

1 items selected. 0.19MB/200 MB

Share With: Share with specific people

Share via: [Email icon] [Input field]

☒ Advanced Options

Expiration Date: 2021-06-15

☐ PIN Code [Enter PIN Code]

☒ Share the patient's name

☒ Allow comments

By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

Close Get Link

Share With

- Share with anyone
 - Creates a shared link. You can share the files with anyone.
- Share with a Medit Link Organization
 - Creates a shared link so you can send the link to your partners.
- Share with specific people
 - Creates a shared link so you can send the link to specific people.

Share via

If you select "Share with a Medit Link Organization" or "Share with specific people," you can add the specific people using email address, phone number, Medit Link partner. Then you can share the link with them.


Advanced Options

- Expiration Date
 - Set the expired date of the link. After the link expired, the people you shared the link with will no longer be able to access it. You can extend the expiration period up to 90 days from today.
- PIN Code
 - Set a 6-digit access code to ensure data security.
- Share the patient name
 - If you check this option, all people who are connected to link will see the patient name.
- Allow comments
 - If you check this option, you will be able to communicate through the shared link.
 - If you uncheck the option, the recipient can only view your comments, but will not be able to add the ones of their own.


Case Talk Details

You can see files in the shared case in detail. In case detail, the below functions are provided.

In the Medit Link Web

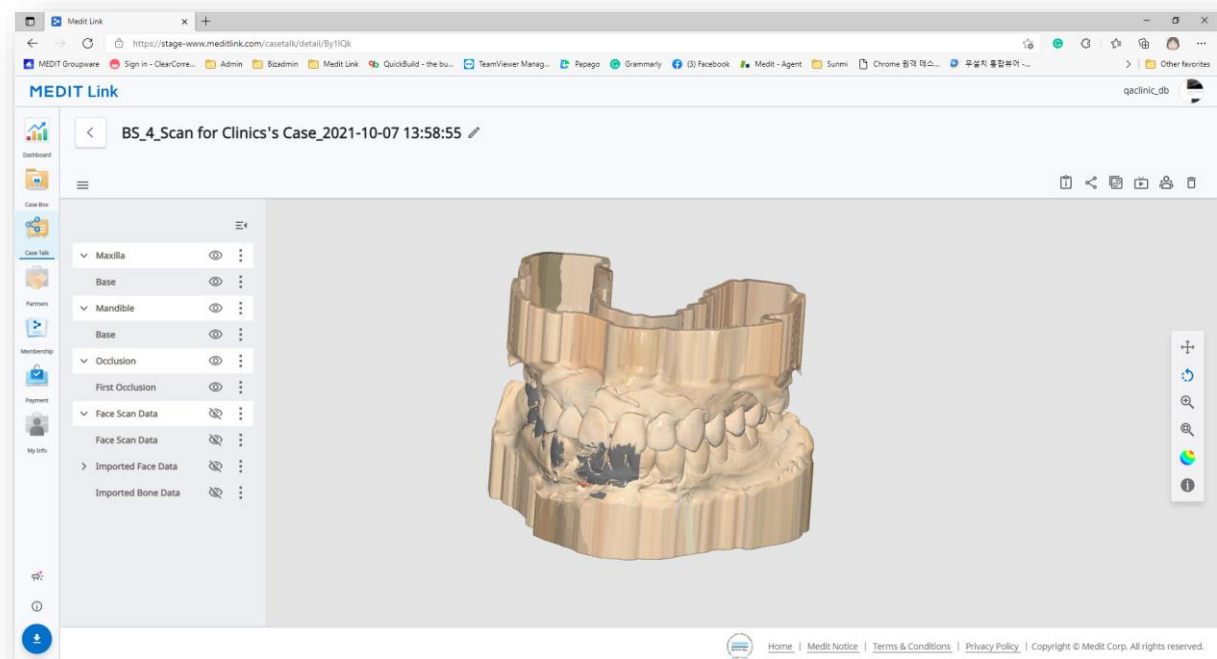

In the Medit Link application

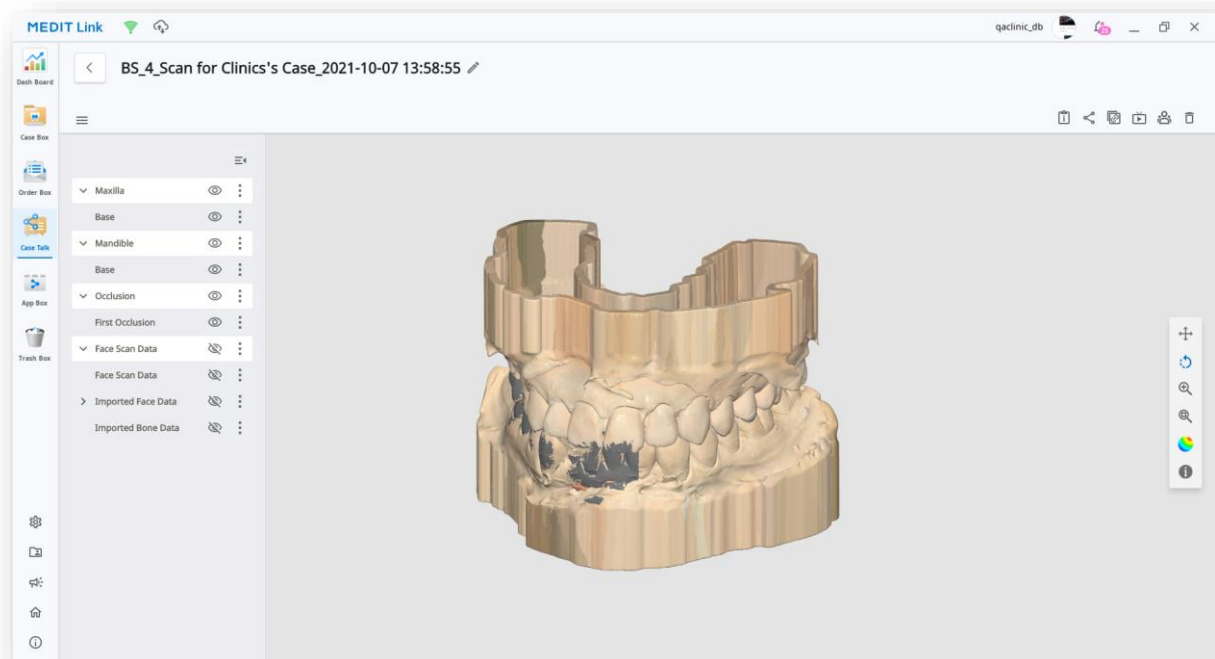


	Case Information	Shows detail information of the case.
	Share	Modifies share setting.
	Copy link	Copies the shared link to the clip board. You can see the shared item through the web browser using the copied link.
	Live	Share your screen with the connected users. The screen sharing is synchronized with the connected users.
	Attendees	Shows all users who are connected to the shared link.
	Delete	Deletes shared link. All shared items will be deleted and the people who try to access to the link will no longer be able to access it.

3.3.2 Web Viewer


Use a browser to view shared data in Case Talk.

It is possible to communicate with the people who have access to the link using annotations and comments.

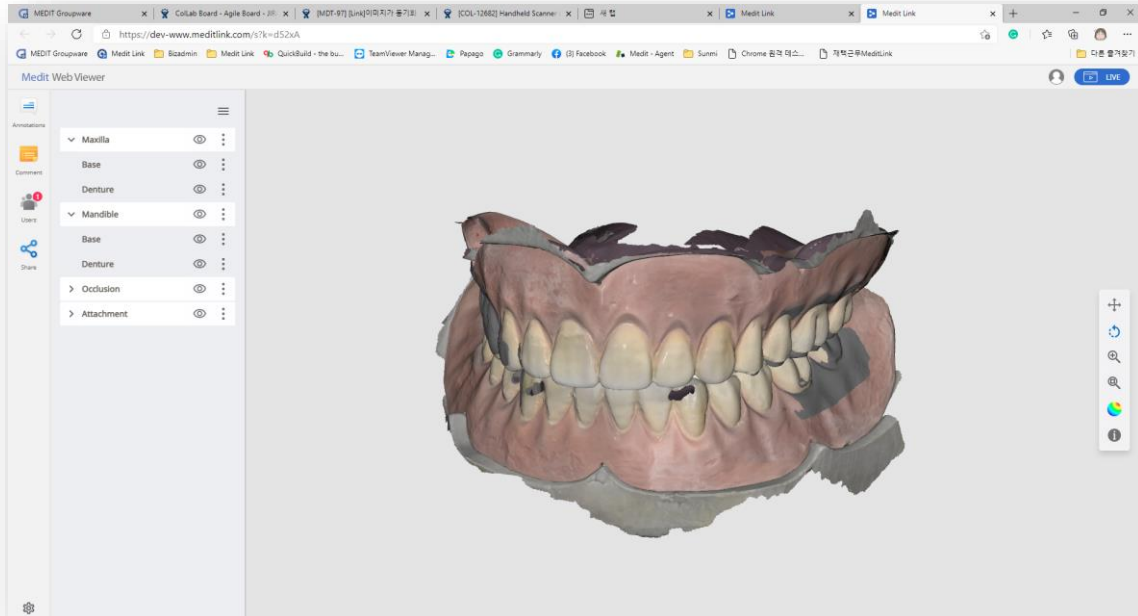


Annotation

You can add annotations to a specified location on the 3D model. The annotations you add can be viewed by anyone with the link.

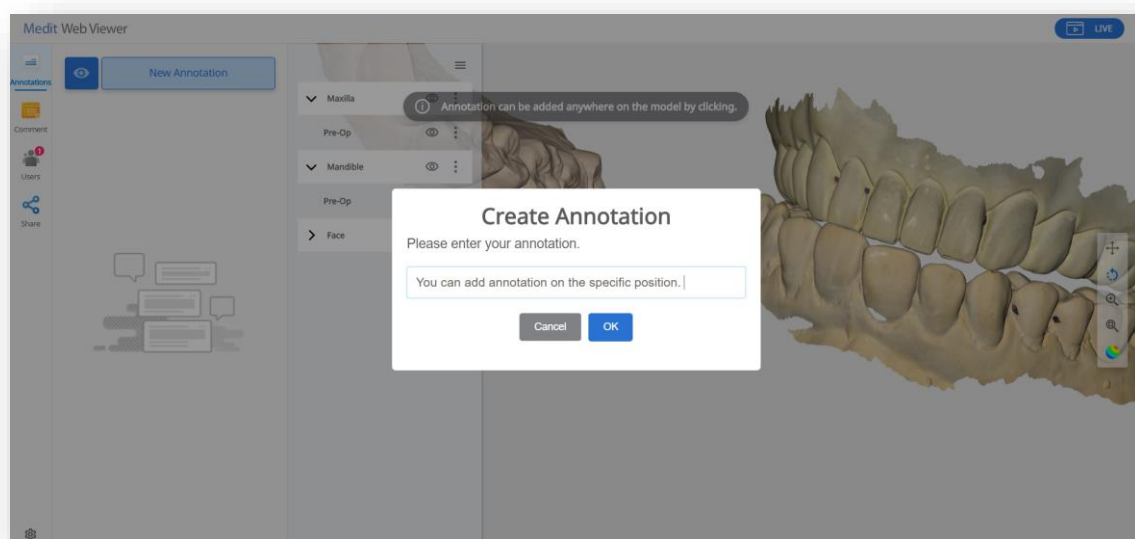


- ① Click on “New Annotation” on an any specific spot.

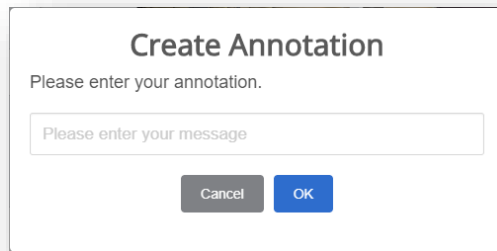

A dialog box titled "Create Annotation". It contains the text "Please enter your annotation." followed by a text input field with the placeholder "Please enter your message". At the bottom, there are two buttons: "Cancel" and "OK".

- ② Type what you want to deliver, then click “OK” to create an annotation.

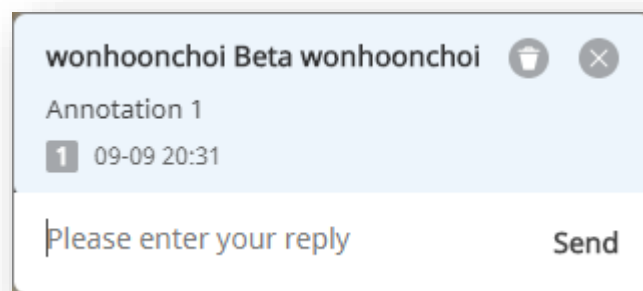

An annotation card with a light blue header. The header contains the text "wonhoonchoi Beta wonhoonchoi" and two icons (trash and close). Below the header, it says "Annotation 1" and "1 09-09 20:31". The main body of the card has a text input field with the placeholder "Please enter your reply" and a "Send" button.

- ③ Enable the conversation through annotation.

Comments

Communicate with users who have access to the link.

Users

Displays the currently connected users.


Settings

Change sharing settings.

Real-time screen sharing

Share your screen in real time.


- ① Click the “LIVE”  to share your Web Viewer screen with the connected users.

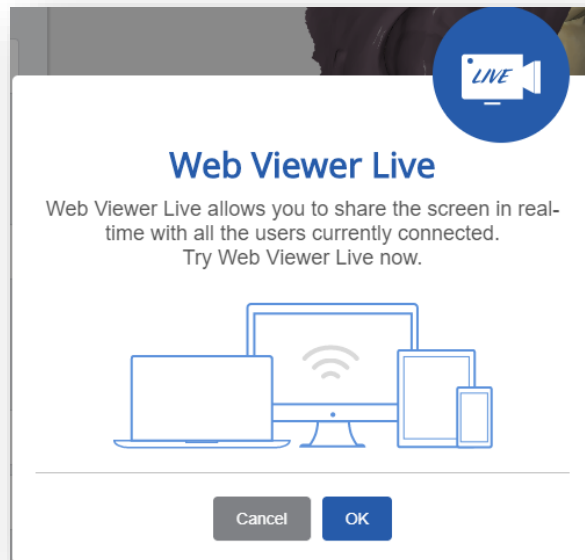

- ② Click “OK” to start sharing your screen.
- ③ Other people who have accessed the link accept screen sharing.
- ④ The screen of the user who requested screen sharing is synchronized with the screen of the person who received the request.

3.3.3 Data Tree

Provides the ability to group and display data according to attributes, view/hide data, adjust transparency, and so on.

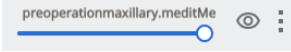Show/Hide	Click on the  icon to adjust the display status of your data.
Opacity adjustment	When you mouse over an element, a slider bar is displayed to adjust the transparency. Use the slider bar to adjust the transparency of your data. 
Extended menu	When you mouse over an element, the Extensions menu is displayed.  ● Show This Only: Hides all other but the current one. ● Hide: Hides the current data.

3D Model Control

	Pan	Moves the model.
	Rotate	Rotates the model.
	Zoom in/out	Zooms in/out the model.
	Zoom Fit	Fits the model to the screen.
	Change Color Display Mode	Changes the color display mode for the model.

3.3.3.1 Recommended browsers

PC	Chrome, Firefox, Safari, Edge (Chromium)
Android	Chrome, Firefox
iOS	Safari, Chrome, Firefox

3.3.3.2 Minimum specifications for mobile devices

Device	Android Same or above the specs of Galaxy S7	iOS Same or above the specs of iPhone 7
OS	Android 7.1 Nougat	iOS 11.0
RAM	Above 3 GB	Above 2 GB
AP	Above Snapdragon 710	A10 Fusion


3.4 Partners

You can search for a partner to request a contract, and see which partners are currently contracted.

3.4.1 Partner Search


You can search for the Medit Link users based on their business name and/or address.

In the search results, you can request Partnership from an organization you would like to work with.

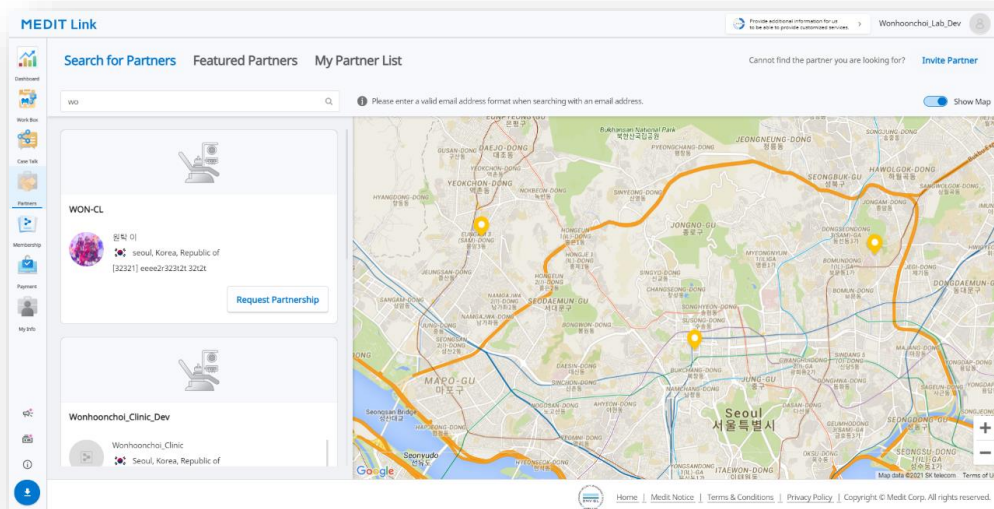

3.4.2 Recommended Partners

Medit Link displays its users located close to you based on your address information.

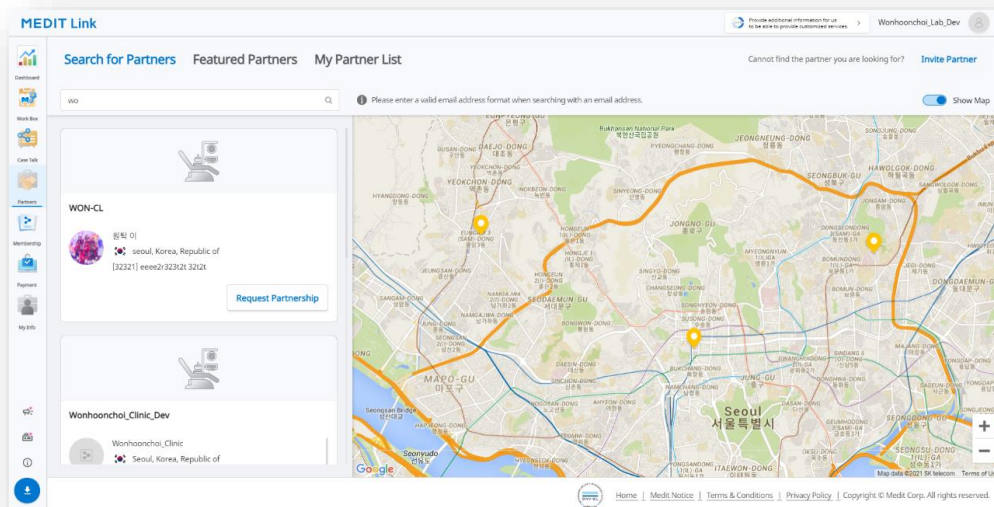


3.4.3 My Partner List

My Partner List displays a list of the companies that you currently have a contract with.



On Partnership vs. Contract


- If you are in a partnership with a business, you can request a contract with the organization.
- Once you request a contract, you can proceed with payment for the placed orders.

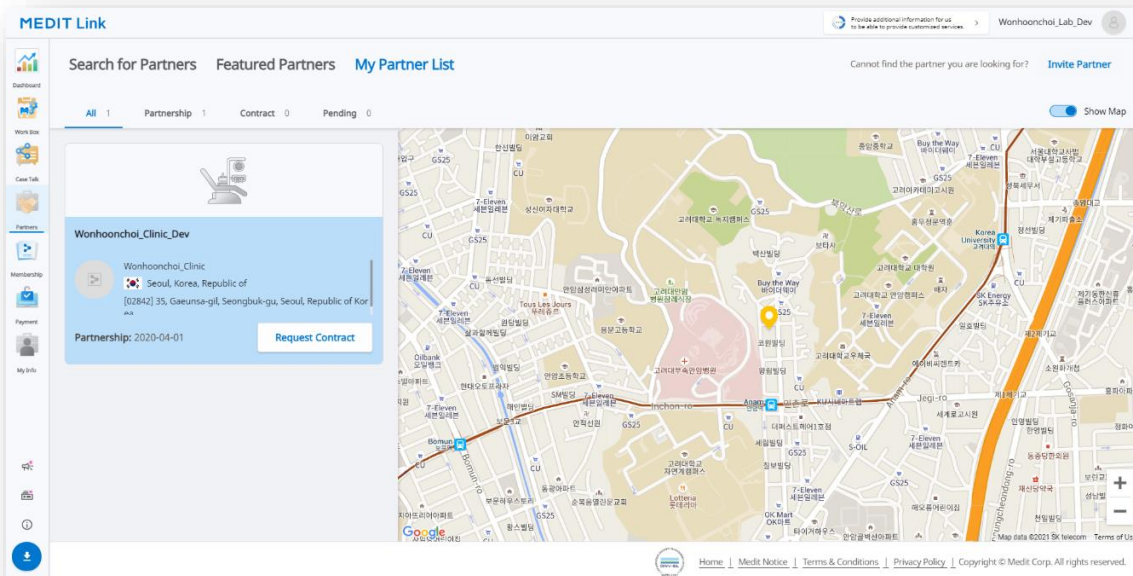

3.4.4 Partner Invitation

You can proceed with a partnership request with a subscription request to a user who is not a member of Medit Link.

- ① Click “Invite Partner” to display the following screen:

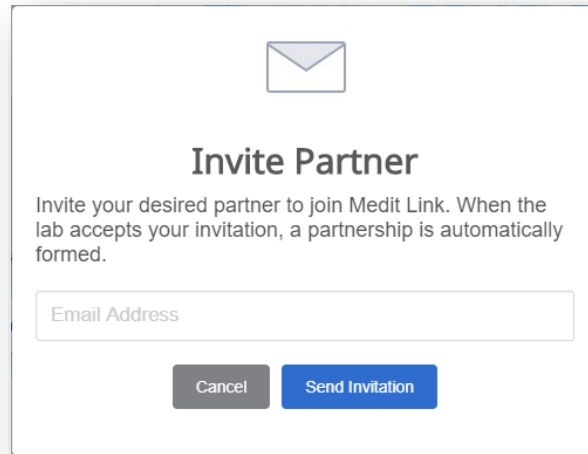A screenshot of a web-based dialog box titled "Invite Partner". At the top center is an envelope icon. Below the icon, the title "Invite Partner" is displayed in a bold font. Underneath the title is a short instruction: "Invite your desired partner to join Medit Link. When the lab accepts your invitation, a partnership is automatically formed." Below this text is a text input field with the placeholder text "Email Address". At the bottom of the dialog box are two buttons: a grey "Cancel" button on the left and a blue "Send Invitation" button on the right.

- ② Enter the email address of the vendor you want to request a partnership from and send the invitation mail to proceed with the invitation.
- ③ If the invited company confirms the membership and proceeds with the membership, the partnership will be established, and the Medit Link service based on this will be available.

3.5 Membership


Medit Link's membership policy aims to establish a more extensive digital environment, gradually expanding to provide various services based on the cloud platform. We aim to provide our customers a more consistent and personalized experience with our new membership policy.

The available membership plans are

- Standard
 - A free plan for every member
 - 1GB cloud storage
 - Efficient for new users
- Premium
 - A paid plan
 - 10 TB of cloud storage
 - Efficient for users with many cases and various data
 - One-month free trial included


	Standard	Premium
Cloud Storage*	~ 1GB	~ 10TB
Price (Monthly)	Free	\$0.99/Month
*Raw data will be excluded from the cloud storage.		

3.5.1 Membership

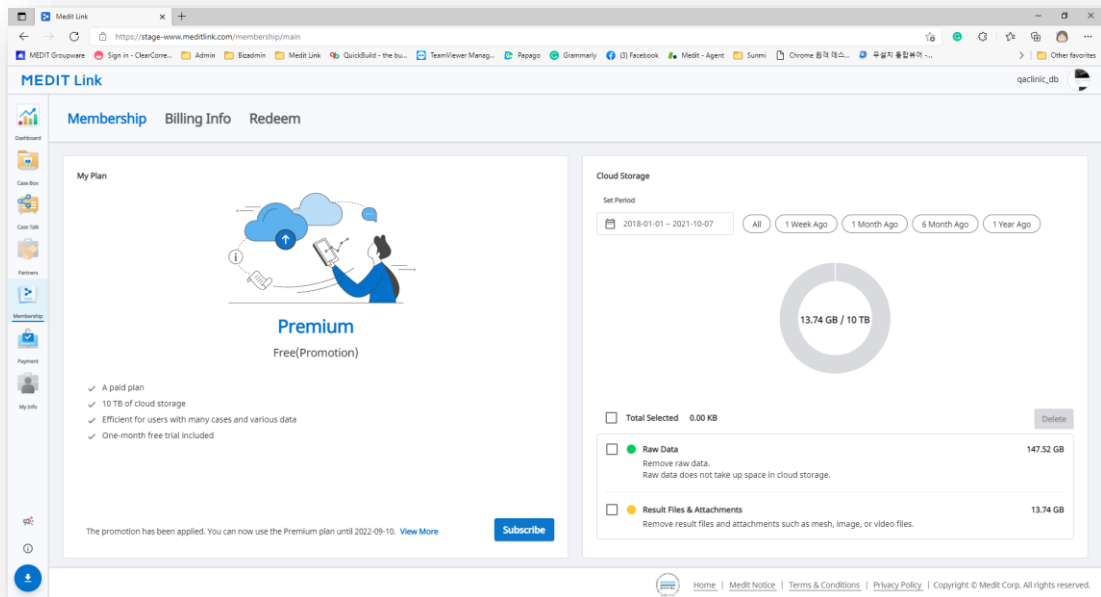


My Plan

You can check the membership plan you are using now. You can also subscribe to or unsubscribe from the membership when you want.

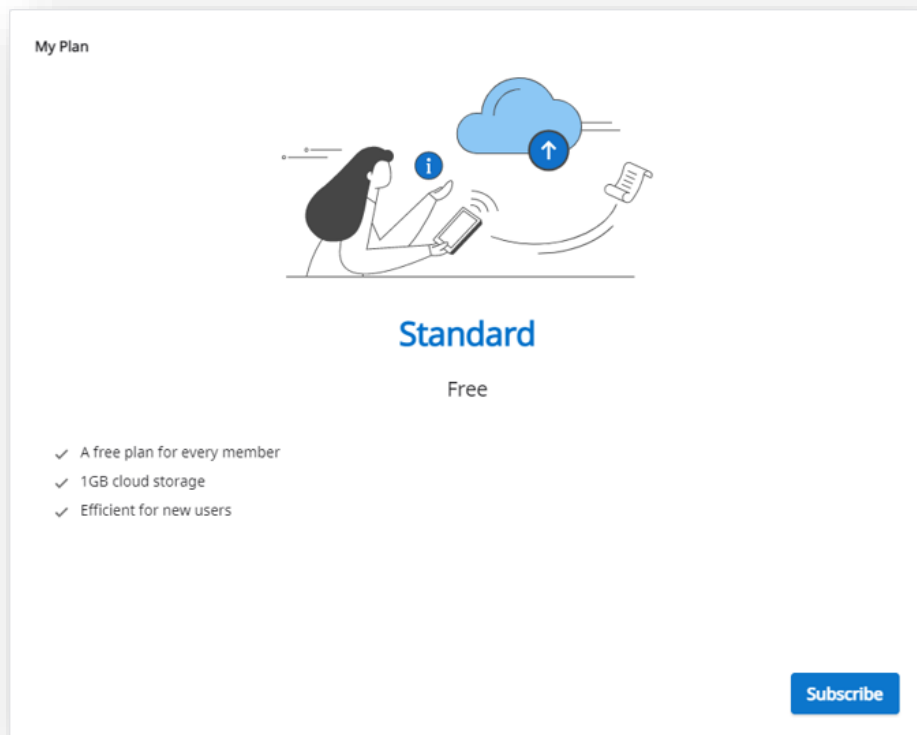

Cloud Storage

You can manage your cloud storage.

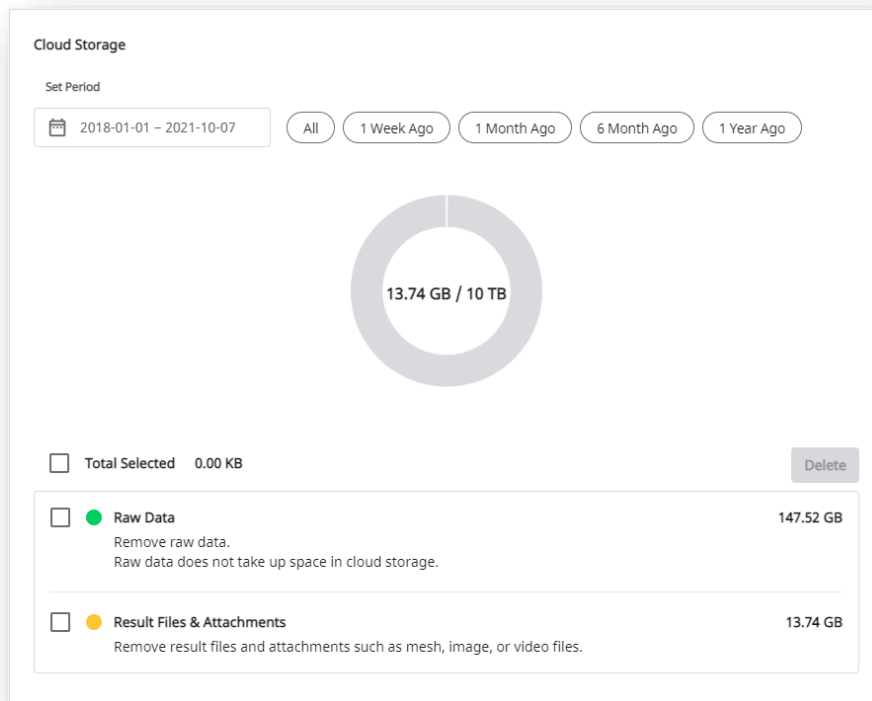


Set Period	Set the period. You can select the files to remove in the period you set. ● All: See all files. ● 1 Week Ago: See the files that have passed a week after creating. ● 1 Month Ago: See the files that have passed a month after creating. ● 6 Months Ago: See the files that have passed six months after creating. ● 1 Year Ago: See the files that have passed a year after creating.
Raw Data	The raw data does not take up the cloud storage, so it will not be in the pi-chart. You can manage the raw data as well for a specific purpose, such as data security.
Result & Attachment	The result and attachment files take up the cloud storage. If the cloud storage is full, you can secure the storage space by removing unnecessary files.
Buy Storage	If you use Standard membership, you can subscribe to the Premium plan to expand to more cloud storage.
Delete	Remove selected files from the cloud storage to free up the space.


Subscribe Membership

If you are using the Standard plan or using the Premium plan with a redeem code, you can subscribe to the Premium membership. The payment system is based on Stripe. (www.stripe.com)

- ① Click the “Subscribe” button.

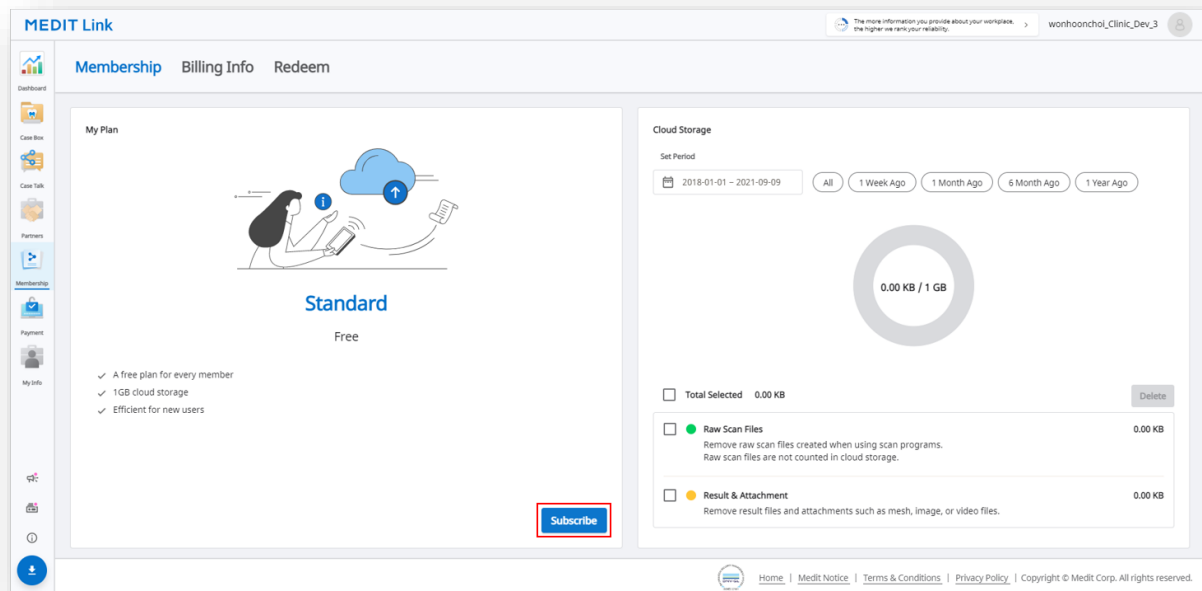

- ② Check the description for the plan and click “Subscribe.”

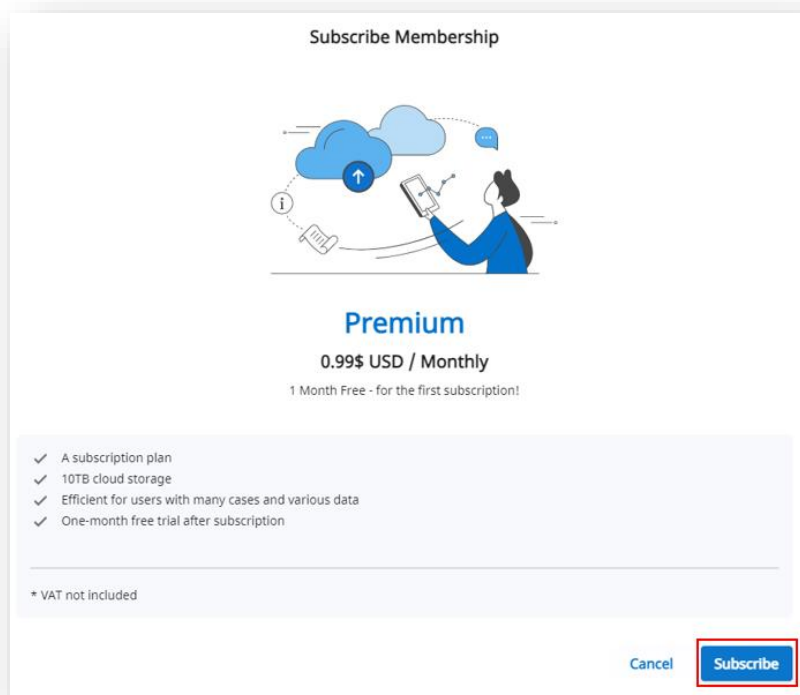

- ③ Check the price and conditions and input necessary information for the payment. Then, click “Subscribe.”

← MEDIT Link TEST MODE

Subscribe to Premium membership for Clinic

\$0.00

Then \$0.99 per day

Medit Link Membership

Premium membership for Clinic

Billed daily

\$0.99

Subtotal

\$0.99

👉 1 MONTH FREE – for first subscription

- \$0.99

100% off for a day

Total due today

\$0.00

Powered by stripe | Terms Privacy

Pay with card

Email

Card information

1234 1234 1234 1234

VISA

MasterCard

Amex

Discover

MM / YYCVC

Name on card

Country or region


South Korea

☐ Save information to pay faster next time ⓘ

Subscribe

By confirming your subscription, you allow MEDIT SG PTE. LTD. to charge your card for this payment and future payments in accordance with their terms.


- ④ The payment will take several seconds. When it's done, you will see the below message.

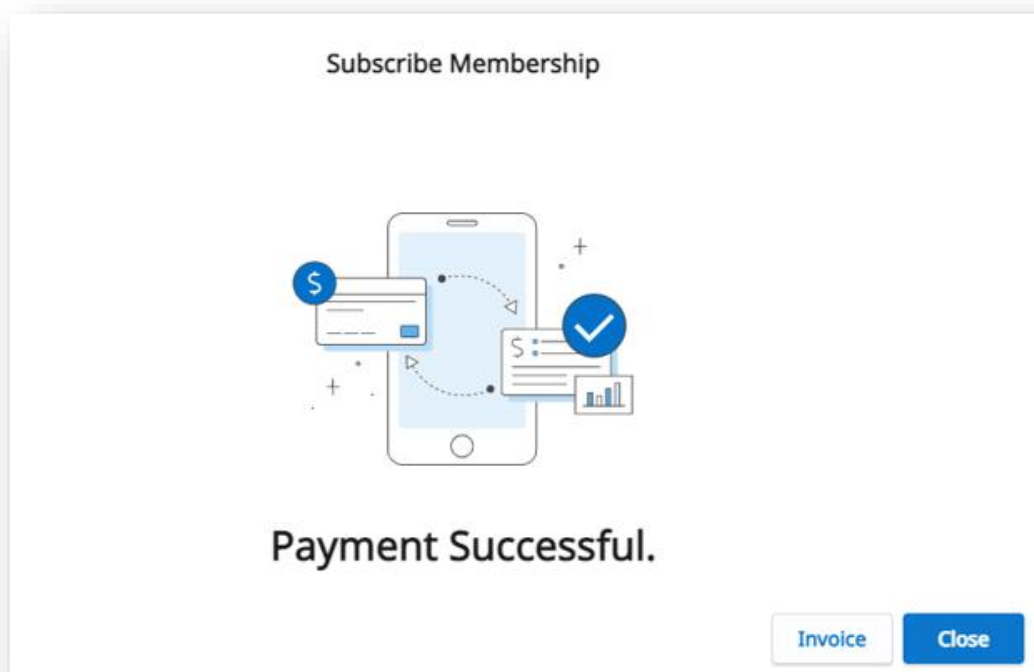


Unsubscribe Membership

You can unsubscribe from the membership whenever you want. After unsubscribing, you can still use the Premium membership for the remaining days. The plan will be changed to the Standard plan from the next payment day, and the cloud storage will be limited to 1 GB. If the cloud storage capacity is exceeded, the files will not be backed up, and the order and share functions will not be available.


① Click “Unsubscribe.”

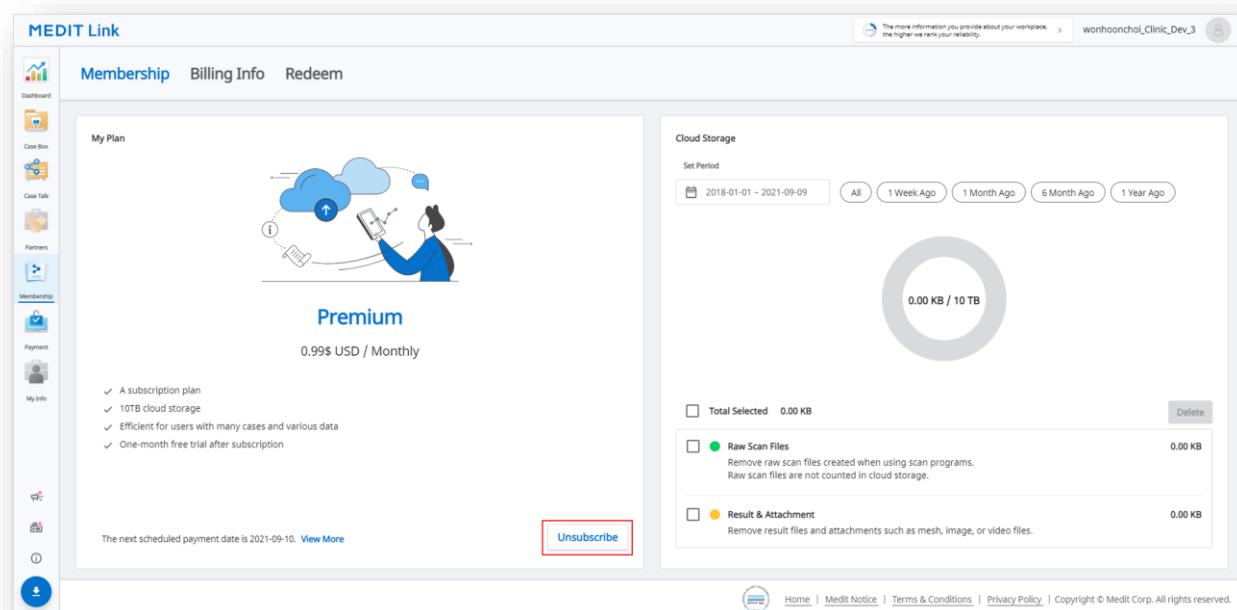

② Click “Cancel Subscription.”

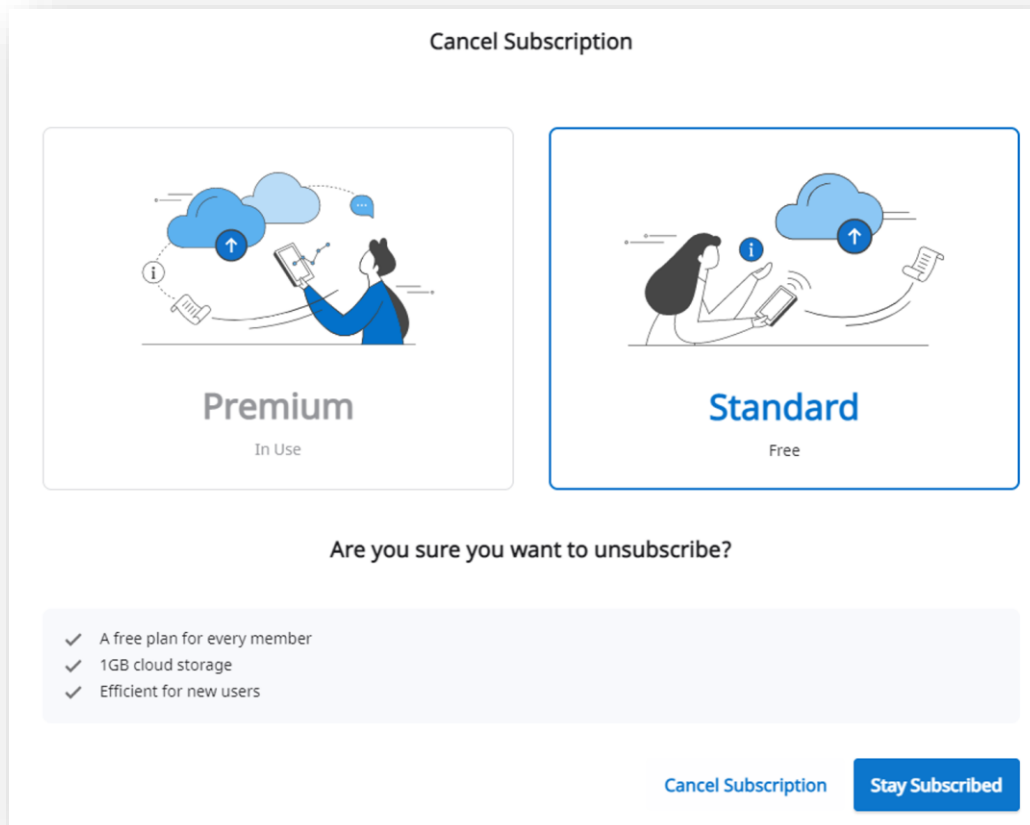

③ Click “Cancel Subscription.”

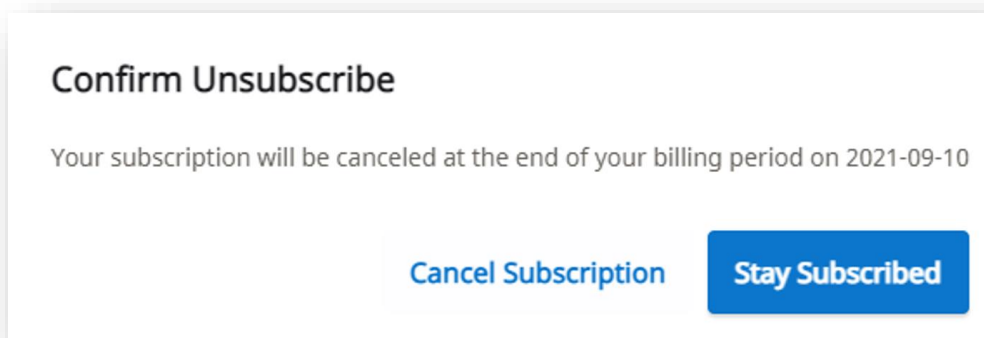

- ④ The subscription has been canceled.

The plan will be changed to the Standard plan from the next payment day, and the cloud storage will be limited to 1 GB. If the cloud storage capacity is exceeded, the files will not be backed up, and the order and share functions will not be available.

Confirm Unsubscribe

Your subscription has been canceled.

Your current plan can use up to 2021-09-10. After that, the cloud capacity is limited to 1GB.

If the cloud capacity is exceeded, the files will not be backed up, and the order and share functions will not be available.

Close

3.5.2 Billing Info

You can check the payment history for membership subscriptions, including an invoice for each payment. You can also add or change the payment method whenever you want. The new payment method will be applied from the next payment.

The screenshot displays the 'Billing Info' section of the MEDIT Link user interface. The top navigation bar includes 'Membership', 'Billing Info' (selected), and 'Redeem'. The left sidebar contains icons for Dashboard, Case Box, Case Talk, Partners, Membership, Payment, and My Info. The main content area is titled 'My Plan and Billing Information' and includes the following details:

- Plan: Premium
- Amount: 0.99\$ USD / Monthly
- Next Payment Date: 2021-09-09
- Payment Method: Visa *****4242 (11/21) | Stripe [Edit Payment Method](#)


Below this information, a message states: 'If you haven't subscribed yet? Experience a differentiated premium.' The 'Billing History' section contains a table with the following data:

Payment Date	Membership	From	To	Payment Amount	Status
2021-09-09	Premium	2021-09-09	2021-09-10	0.99 \$ (USD)	Payment Completed Invoice
2021-09-09	Premium	2021-09-09	2021-09-09	0 \$ (USD)	Payment Completed Invoice


The bottom of the page features a footer with links: Home, Medit Notice, Terms & Conditions, Privacy Policy, and Copyright © Medit Corp. All rights reserved.

Add Payment Method

- ① Click “Edit Payment Method.”

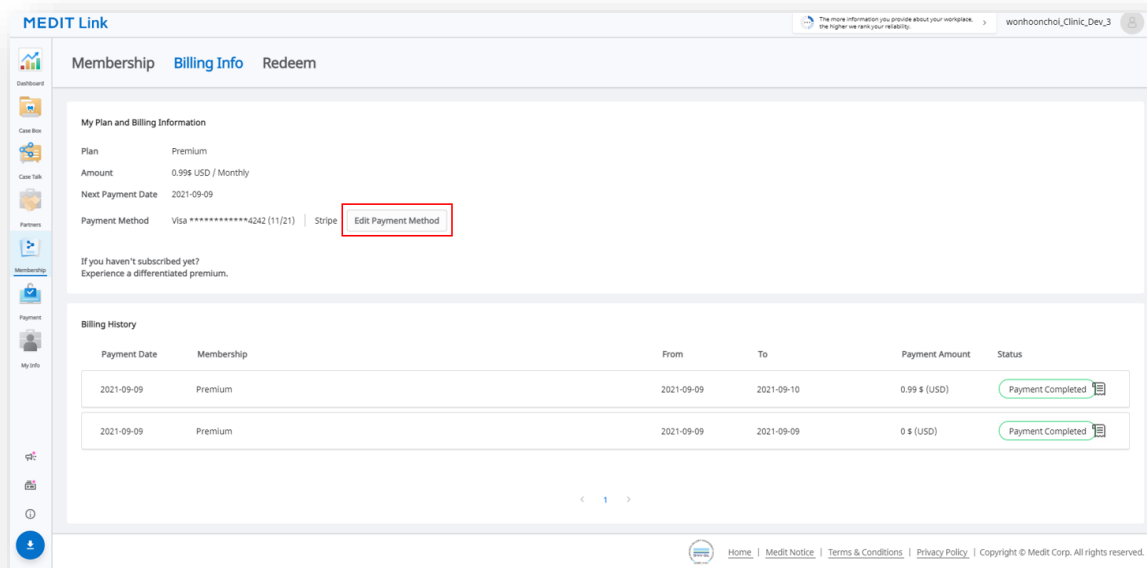

- ② Click “Add payment method.”

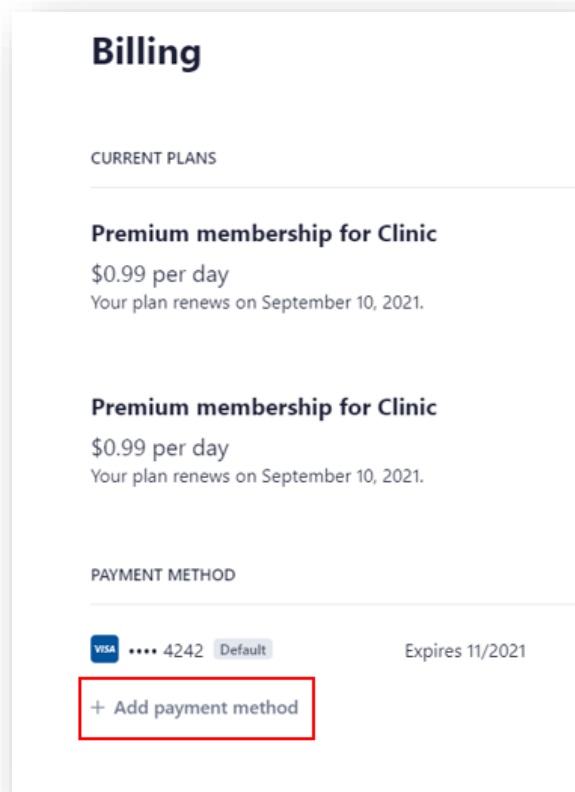

- ③ Input the required information and click “Add.”

Add payment method

Card information

 Card number

MM / YY

☒ Use as default payment method

Add

Go back

By adding a payment method, you agree to MEDIT SG PTE. LTD.'s Terms of Service and Privacy Policy.

- ④ The newly added payment method will be set as default and used from the next payment date.

Billing

CURRENT PLANS

Premium membership for Clinic

\$0.99 per day

Your plan renews on September 10, 2021.

Premium membership for Clinic

\$0.99 per day

Your plan renews on September 10, 2021.

PAYMENT METHODS

 **** 5556 Default

Expires 11/2021

×

 **** 4242

Expires 11/2021

...

+ Add payment method

Change Payment Method

- ① Click "...".

Billing

CURRENT PLANS

Premium membership for Clinic

\$0.99 per day

Your plan renews on September 10, 2021.

Premium membership for Clinic

\$0.99 per day

Your plan renews on September 10, 2021.

PAYMENT METHODS

 **** 5556	Default	Expires 11/2021	×
 **** 4242		Expires 11/2021	...

+ Add payment method

- ② Click “Make default” to change the default payment method.

Billing

CURRENT PLANS

Premium membership for Clinic
\$0.99 per day
Your plan renews on September 10, 2021.

Premium membership for Clinic
\$0.99 per day
Your plan renews on September 10, 2021.

PAYMENT METHODS

VISA

**** 5556

Default

Expires 11/2021

×

VISA

**** 4242

Expires 11/2021

...

+ Add payment method

Make default

Delete

BILLING INFORMATION

3.5.3 Redeem

Register the redeem code provided as a part of the promotion. You can use premium membership with redeem code.

Register and Activate Redeem Code

- ① Input redeem code and click “Register.”

The screenshot shows the 'MEDIT Link' web application. The top navigation bar includes 'Membership', 'Billing Info', and 'Redeem'. The 'Redeem' tab is active. Below the navigation bar, there is a 'Register Redeem Code' section with a text input field for the 'Redeem Code' and a blue 'Register' button. Below this is a 'My Redeem Code' section with a table. The table has columns: Promotion, Code Number, Dealer Name, Registration Date, From, To, and Status. The table is currently empty, displaying 'No search results were found.' at the bottom. The footer of the application includes links for Home, Medit Notice, Terms & Conditions, Privacy Policy, and a copyright notice for Medit Corp.

- ② Click “Register” or “Apply.”

The modal dialog titled 'Register Redeem Code' provides information about the registration process. It states: 'You can register for the promotion or apply immediately. Register allows you to apply and activate when you need it, If applied immediately, it will start after the expiration of the active promotion. Periodic payments are stopped automatically.' At the bottom of the dialog, there are three buttons: 'Close', 'Register', and 'Apply'.

Register: Register the redeem code. You can activate and use the redeem code you registered.

Apply: Register and activate the redeem code. The promotion starts immediately.

Activate Redeem Code

You can have multiple redeem codes as “Registered” status. In this case, you can activate the code you want.

① Click “Activate.”

If there isn't any activated code, the one you choose will be started immediately. If there is an activated code, the other code you choose will be in “Waiting” status, and it will be activated automatically right after the current one expires.



Only one code can be in “Waiting” status.


The screenshot shows the MEDIT Link web application interface. The top navigation bar includes 'Membership', 'Billing Info', and 'Redeem'. The 'Redeem' section is active, displaying a 'Register Redeem Code' form and a table of 'My Redeem Code'.

Promotion	Code Number	Dealer Name	Registration Date	From	To	Status
Test Promotion - 1 month	[Redacted]	Medit Link	2021-09-09	-	-	Registered Activate
Test Promotion - 1 month	[Redacted]	Yoshida	2021-09-09	-	-	Registered Activate
Test Promotion - 1 month	[Redacted]	Medit Link	2021-09-09	-	-	Registered Activate
Test Promotion - 1 month	[Redacted]	Medit Link	2021-09-09	-	-	Registered Activate


3.6 Payment

3.6.1 Payment Details (Clinics, Labs accounts)

You can check all payment history and transactions through the Payment Detail page.

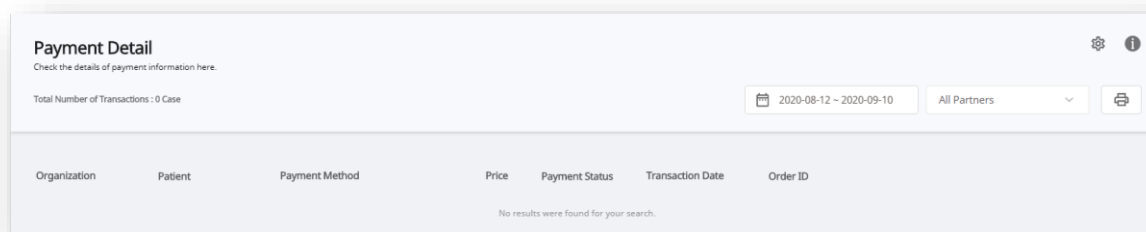

3.6.2 Payment account information (Lab account only)

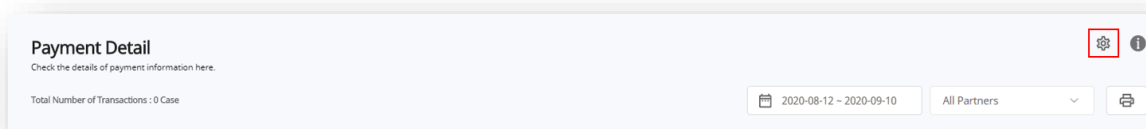

This feature is only available to the administrator of Lab account.

- You must enter your payment account information so that your account can use the online payment system when ordering.



If you select a country that can support Stripe, you will be notified of the Stripe connection, and you can create a Stripe account by pressing the appropriate button.

3.7 My Information

3.7.1 My Account

You can view/modify the information you entered during the registration and check your membership.

The screenshot shows the 'My Account' page with a navigation bar at the top containing 'My Account', 'My Organization', and 'My Price List'. The main content is divided into two columns. The left column, titled 'Account Information', contains a profile picture placeholder, an 'Email Address' field with the value 'coldgreentea@naver.com' and a 'Change Email Address' button, a 'Password' field with a 'Change Password' button, a 'Name' field with the value 'WonhoonChoi_Lab_Stage', a 'Language' dropdown menu set to 'English(United States)', and a 'Country' dropdown menu set to 'United States'. Below these fields is an 'Organization Information' section with an 'Edit Information' button. The right column, titled 'Group Information', shows a list of groups with one entry: 'WonhoonChoi_Lab_Stage' with a user icon, the email 'coldgreentea@naver.com', and the registration date '2020-07-03'.

3.7.2 My Organization (Admin account only)

View/modify information of your organization.

The screenshot shows the 'My Organization' page with a navigation bar at the top containing 'My Account', 'My Organization', and 'My Price List'. Below the navigation bar is a sub-navigation bar with four tabs: 'Address Information' (active), 'Organization Information', 'Additional Information', and 'Scanner Information'. The main content area is titled 'Address Information' and is marked as 'Mandatory *'. It contains several form fields: 'Address1 *' with the value '23 Incheon-ro 22-gil, Seongbuk-gu, Seoul, Korea' and a 'View Location On Map' button; 'Address2' with the value 'Apartment, Suite, Unit, Building, Floor, etc.'; 'Zip/Postal Code *' with the value '02835'; 'City *' with the value 'Seoul'; 'Country *' with a dropdown menu set to 'Korea, Republic of'; and 'State/Province *' with a dropdown menu set to 'Seoul'. There are also informational icons and text: 'Street Address, P.O. Box, Company Name, C/O' and 'Please enter the correct address in order for labs to check your company's location on the map.' at the top, and 'Apartment, Suite, Unit, Building, Floor, etc.' and 'Please select your country information carefully because it is important for completion of your transactions.' at the bottom.

3.7.3 My Price List (Lab: Admin account only)

Enter the products offered by your lab and the price information for each product.



You can make a transaction based on the product list and price information if you have a contract relationship with the client.



Medit Link does not support price-based transactions for partnerships.

3.8 Notice and Updates

Displays notifications about using and updating the Medit Link service.

3.9 Help Center

Once you press “Help Center” you get transferred to Medit’s Help Centre which contains information about the company’s products and their functions.

3.10 Download the App

Download the installer for the latest version of the program.

- You can install the program on another computer if you have it already installed after the purchase of the product.



Program updates are automatic when the Internet is connected.

4 Medit Link App Overview

MEDIT Link Meditcompany

Case Box

Filters: All (selected), Form Only, Need Processing, Scan Completed, CAD, Ordered, Case Completed

Total Count: 462 | Date Range: 2018-01-01 ~ 2021-06-14 | Search / Contains

Status	Case Name	Patient Name	Birth Date	Form Information	Date Modified	Date Scanned	Date Ordered	Lab Name
Scan Completed	孙天柱案例 - Clone	孙天柱	-	12-Crown/11-Pontic/21-Pontic/22-Crown	2021-06-11 오후 2:38	2021-05-17 오후 4:47	-	-
Ordered	Chair Side CAD TEST's ...	Chair Side CAD TEST	-	11-Crown/21-Crown/22-Crown	2021-06-08 오전 11:58	2021-05-25 오후 12:21	2021-06-08 오전 11:58	meditlab
Form Only	webAPI test의 케이스	webAPI test	-	11-Crown/21-Crown	2021-06-03 오전 10:07	-	-	-
Form Only	addCaseItem의 케이스	addCaseItem	-	38-Inlay/37-Inlay	2021-06-02 오후 6:55	-	-	-
Scan Completed	孙天柱案例	孙天柱	-	12-Crown/11-Pontic/21-Pontic/22-Crown	2021-05-26 오후 6:19	2021-05-17 오후 4:47	-	-
Scan Completed	test's Casetest	test	-	11-Crown/21-Crown/22-Crown/23-Crown/...	2021-05-25 오후 7:23	-	-	-
Form Only	cha의 케이스	cha	-	37-Custom Abutment/36-Custom Abutme...	2021-05-25 오후 3:23	-	-	-
Need Processing	오토 알로드 온	스캔 알로드 테스트	-	26-Crown	2021-05-25 오후 3:22	2021-05-03 오후 3:28	-	-
Scan Completed	스캔 알로드 테스트와 ...	스캔 알로드 테스트	-	26-Inlay	2021-05-25 오후 12:03	2021-05-03 오후 3:28	-	-
Ordered	cha의 케이스	cha	-	11-Crown	2021-05-21 오후 6:10	2019-02-11 오후 5:47	2021-05-21 오후 6:10	jaeyeonlab
Ordered	orderTEST222의 케이스	orderTEST222	-	26-Crown/27-Crown	2021-05-21 오후 5:34	-	2021-05-21 오후 5:34	meditlab
Ordered	orderTEST의 케이스	orderTEST	-	27-Implant Crown/28-Implant Crown	2021-05-21 오후 5:29	-	2021-05-21 오후 5:29	메디트링크_알_테스트
Ordered	test's Case	test	0001-01-01	11-Crown/21-Crown/22-Crown	2021-05-21 오후 5:27	-	2021-05-21 오후 5:27	메디트링크_알_테스트

1. Title Bar

2. Boxes

3. Detailed Functions

4. Case List

4.1 Title Bar

The title bar displays network status, cloud synchronization, membership grade, and notifications.

4.1.1 Network Status

Medit Link enables you to use all the features of the service seamlessly when you are connected to a network.

The Title Bar displays the network status for the smooth use of the service.

If the network connection is not good, the following restrictions occur:

- Some functions, such as ordering, creating and sharing Case Talk, and cloud synchronization are restricted.
- The overall performance of the service can be reduced due to the impact of the network environment.

In this case, you can enter the Offline Mode to perform basic tasks, and then use the features that require internet once you are connected.




The network speed symbol indicates the response speed between the current area and AWS(Amazon Web Service). It does not indicate your network speed. Sometimes it may show the status as poor or bad, even though your network speed is fast enough.

If the other process uses the same network, the speed may show poor or bad, and the overall performance of Medit Link could be bad. In this case, please do not share the network with other processes. For example, the “Upload Raw Data” feature takes up a huge part in the network process. If the symbol shows poor or bad status, go to Settings > Upload Raw Data and turn “Upload Automatically” off.

4.1.2 Case Synchronization Manager

You can upload and download the cases in Medit Link. You can also upload multiple cases to the cloud or download them to your local PC.

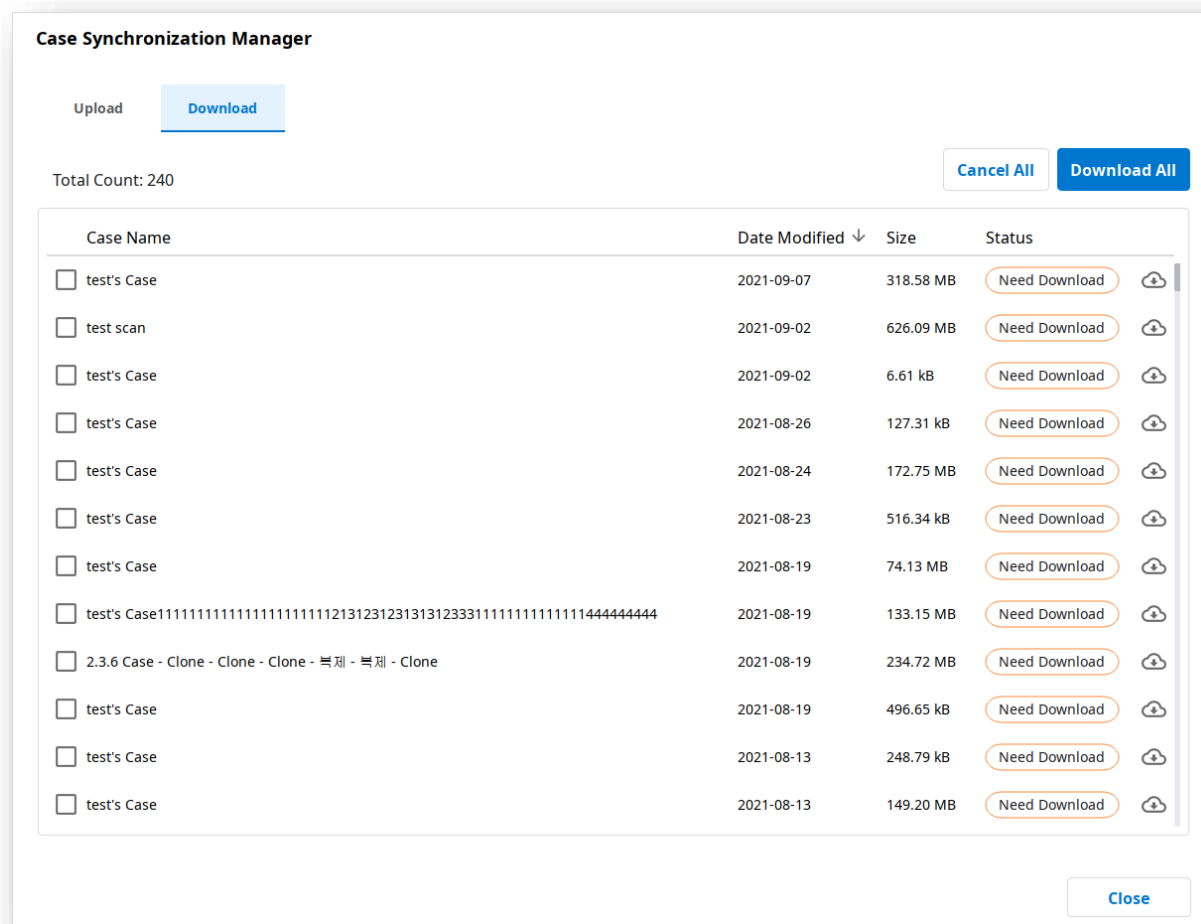

Upload Case


Total Count / Item Selected	Check the number of all/selected cases that need to be uploaded.
Upload All	Upload all cases to the cloud.
Cancel All	Cancel all requests.
	Upload selected cases to the cloud.
	Cancel uploading of the selected cases.
	Retry uploading of the selected cases.

Download Case

Total Count / Item Selected	Check the number of all/selected cases that need to be downloaded.
Download All	Download all cases to the PC.
Cancel All	Cancel all requests.
	Download selected cases to PC.
	Cancel downloading of the selected cases.
	Retry downloading of the selected cases.

4.1.3 User Information

User information displays basic member information such as profile image, organization name, email address, and membership name. It also provides logout option and home page shortcuts.

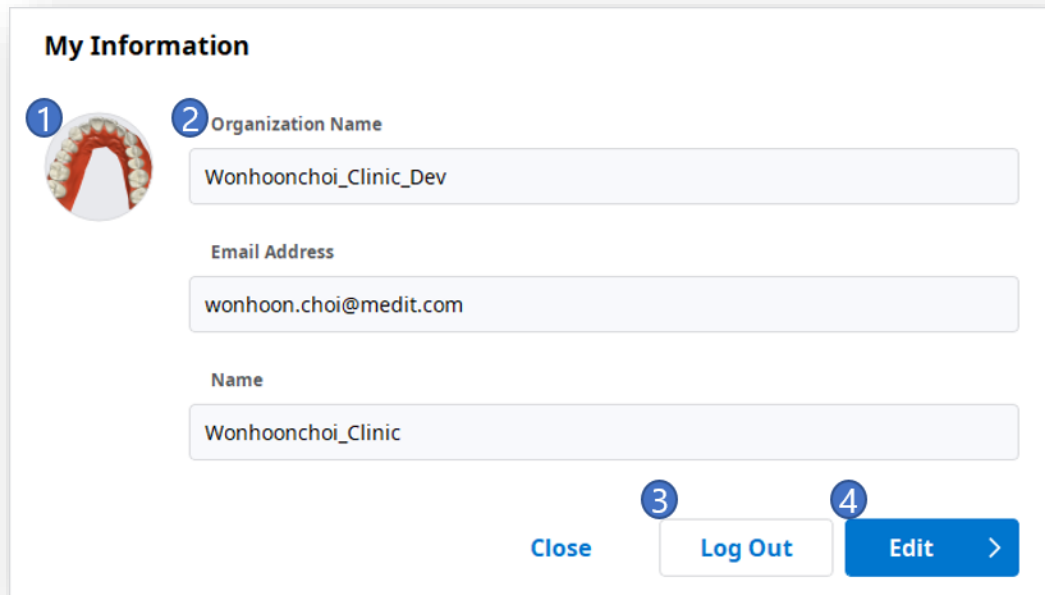

My Information

1 

2 Organization Name
Wonhoonchoi_Clinic_Dev

Email Address
wonhoon.choi@medit.com


Name
Wonhoonchoi_Clinic

Close Log Out Edit >

-
- | | |
|---|-----------------------------|
| 1 | Displays the profile image. |
|---|-----------------------------|
-
- | | |
|---|---|
| 2 | Displays the organization name, email address, and membership name. |
|---|---|
-
- | | |
|---|--|
| 3 | Logs out of the currently connected account. |
|---|--|
-
- | | |
|---|--|
| 4 | Logs in with the currently accessed account information and access the homepage.
Check, modify, and manage registered information |
|---|--|
-

4.1.4 Notifications


Provides notifications related to contracts/orders and work processes.

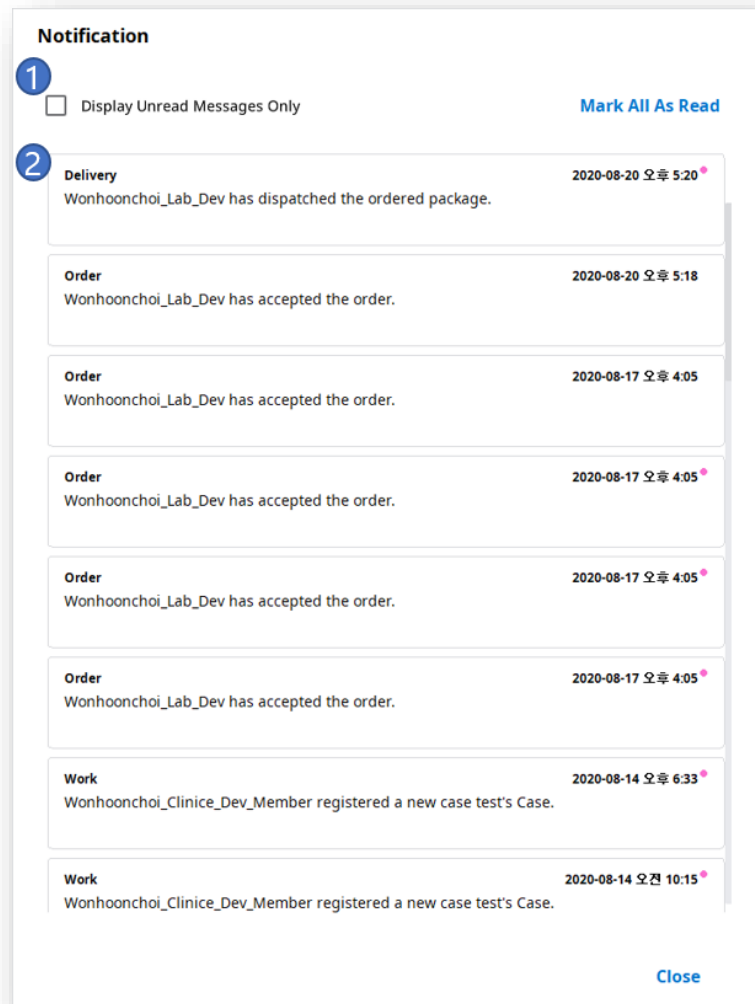

1 Shows only the unread messages.

2 Quickly navigate to the specific case or webpage by selecting the corresponding notification.

4.2 Boxes Overview

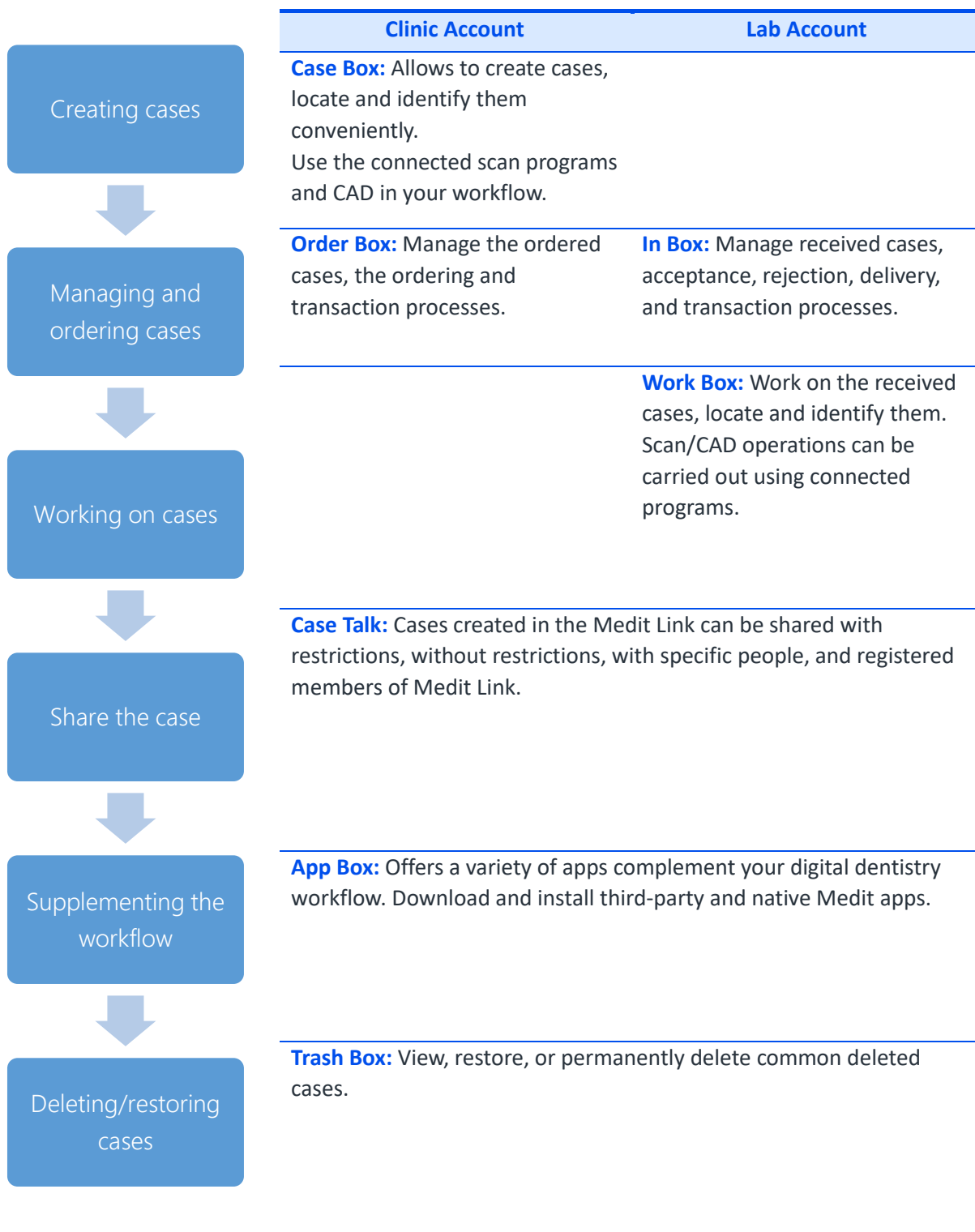


Boxes available in Lab and Clinic accounts

		Case Box	Order Box	In Box	Work Box	Case Talk	Trash Box	App Box
Clinic Account	Admin	O	O	X	X	O	O	O
	Member	O	X	X	X	O	O (Restoring cases)	O
Lab Account	Admin	X	X	O	O	O	O	O
	Member	X	X	X	O	O	O (Restoring cases)	O


4.3 Settings

General Settings: Provides features for overall system setup.

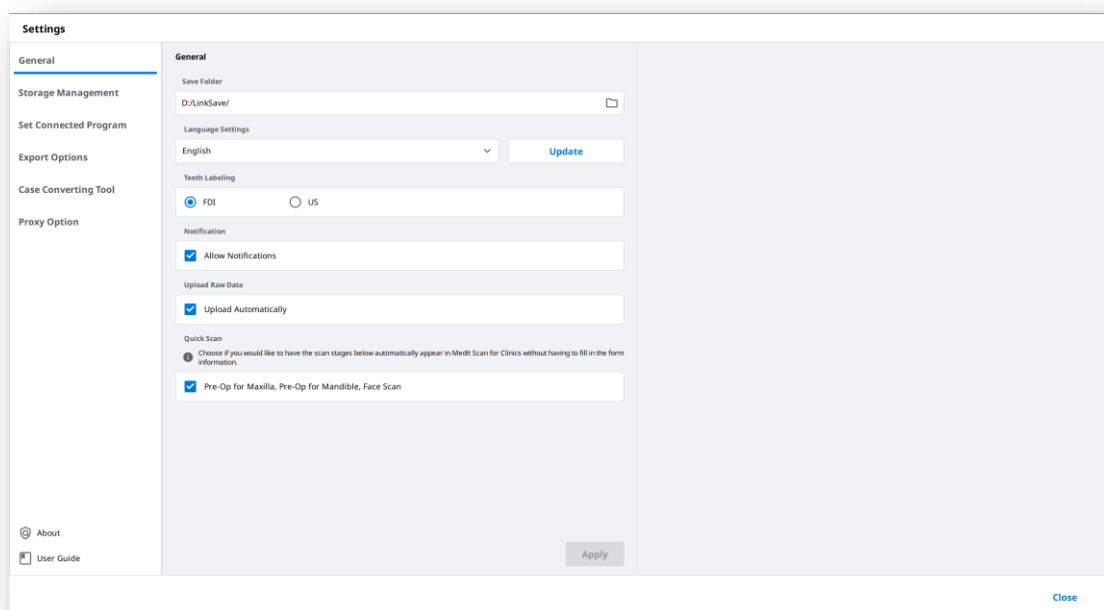

Save Folder	Specify a folder to store the results generated from Medit Link and connected programs.
Languages Settings	Select the program language. Improve your current language and update in real time when a new language is added.
Teeth Labeling	Select a labeling system.
Notifications	Set whether to receive notifications for program updates.
Upload Raw Data	Determines whether to automatically upload raw scan data to the cloud. It helps to back up the original data and perform seamless scanning when using Medit Link with the same account on another PC. The original data cloud backup only applies to data acquired by Medit Scan for Clinics.
Quick Scan	Choose if you would like to have the scan stages below to automatically appear in Medit Scan for Clinics without having to fill in the form information. ● Pre-op for Maxilla ● Pre-op for Mandible ● Face Scan

Storage Management: You can check and manage the local and cloud storage space.

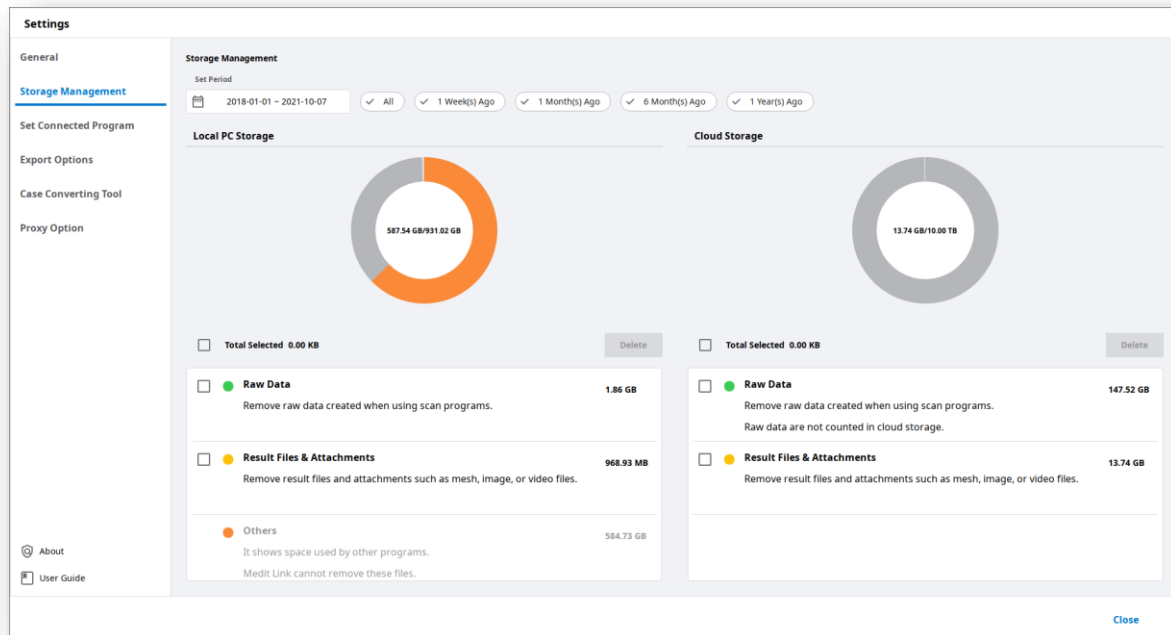

Set Period

Set the period. You can select the files to remove in the period you set.

- All: See all files.
- 1 Week(s) Ago: See the files that have passed a week after creating.
- 1 Month(s) Ago: See the files that have passed a month after creating.
- 6 Month(s) Ago: See the files that have passed six months after creating.
- 1 Year(s) Ago: See the files that have passed a year after creating.

Local Storage

Manage the PC storage space.


Raw Data	Remove the selected raw data from the PC. If the raw data is already uploaded to the cloud, you can download and use it whenever you want.
Result & Attachment	Remove the result and attachment files such as mesh, image, video from the PC. If these files are already uploaded to the cloud, you can download and use them whenever you want.
Others	See the files used to operate the system and other programs. It shows the storage of your PC. You cannot remove these files directly in Medit Link.
Delete	Remove selected files from the PC.

Cloud Storage

Manage the cloud storage space.

Raw Data	The raw data does not take up the cloud storage, so it will not be in the pi-chart. You can manage the raw data as well for a specific purpose, such as data security.
Result & Attachment	The result and attachment files take up the cloud storage. If the cloud storage is full, you can secure the storage space by removing unnecessary files.
Buy Storage	If you use Standard membership, you can subscribe to the Premium plan to expand to more cloud storage.
Delete	Remove selected files from the cloud storage to free up the space.

Set Connected Program: Set/manage programs that work with Medit Link.

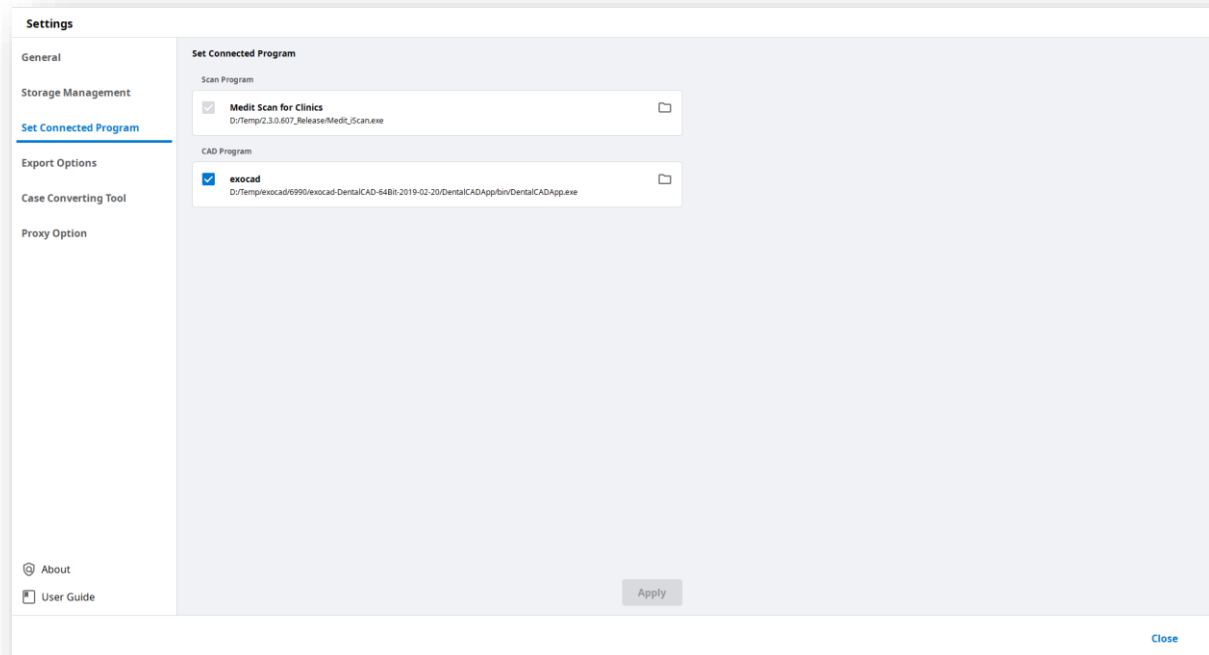

You can connect the following programs:

	Clinic Account	Lab Account
Scanning S/W	Medit Scan for Clinics	Medit Scan for Clinics Medit Scan for Labs
CAD	exocad	exocad

Export Option: Set how to export files from Medit Link.

The screenshot shows the 'Settings' window with the 'Export Options' tab selected. The left sidebar contains links for General, Storage Management, Set Connected Program, Export Options (active), Case Converting Tool, and Proxy Option. The main area is divided into several sections:


- CAD Type:** A dropdown menu currently set to 'exocad'.
- Export Files To:** A text field showing the path 'C:/Users/Medit-Wonhoonchoi/Desktop/' with a folder icon to the right.
- Mesh Format:** Radio buttons for 'Keep Original Mesh Format' and 'Mesh Format' (selected). Below 'Mesh Format' are checkboxes for 'OBJ' (checked), 'PLY', and 'STL'.
- Image Format:** Radio buttons for 'Keep Original Image Format' and 'Image Format' (selected). Below 'Image Format' are radio buttons for 'JPG' (selected), 'JPEG', 'PNG', and 'BMP'.
- Compress & Combine:** Checkboxes for 'Compress Files' and 'Combine Individual Mesh'.
- Options when same folder or file already exists:** Radio buttons for 'Always ask when exporting' (selected), 'Export a new file and overwrite the existing file. The file currently stored will be lost', and 'Create a folder and export the files'.
- Rule for Folder Name:** A section with checkboxes for 'Date Created', 'Date Scanned', 'Date Ordered', 'Patient Name', 'Case Name', 'Partner Name', 'Patient ID', and 'Order ID'. Below these is a text field for 'Folder Name' containing 'Date Created - Patient Name'.
- Rule for File Name:** A section with the same checkboxes as above. Below these is a text field for 'File Name' containing 'Date Created - Patient Name'.

At the bottom right, there is an 'Apply' button and a 'Close' button.

CAD Type	Select the CAD program you want to use for your design. The program will adjust the files based on this selection, allowing you to proceed with the design without any additional axis alignment.
Export Files To	Specify the path where you want to store the results.
Mesh Format	Specify the format of the final mesh file.
Image Format	Specifies the format of the exported image files.
Compress & Combine	Compress Files: Compresses all results into a single file. Combine Individual Mesh: Combines the result of maxilla or mandible with individual tooth. Useful when importing several files into a program is not possible.
Options when same folder or file already exists	Select an action when same folder or file name exists. Always ask when exporting. Export a new file and overwrite the existing files. The file currently stored will be lost. Create a folder and export the files.


Rule for Folder Name	Define the rules for the folder name. Creates a new folder with a defined name and saving the result files in it. Below terms is available for the rule. Clinic account Date Created, Scanning Date, Order Date, Patient Name, Case Name, Partner Name, Patient ID, Order ID Lab account Date Created, Scanning Date, Order Date, Patient Name, Case Name, Partner Name, Order ID
Rule for File Name	Define the rules for the output file name. Below terms is available for the rule. Clinic account Date Created, Scanning Date, Order Date, Patient Name, Case Name, Partner Name, Patient ID, Order ID Lab account Date Created, Scanning Date, Order Date, Patient Name, Case Name, Partner Name, Order ID

Case Converting Tool: Used to move cases between different accounts.

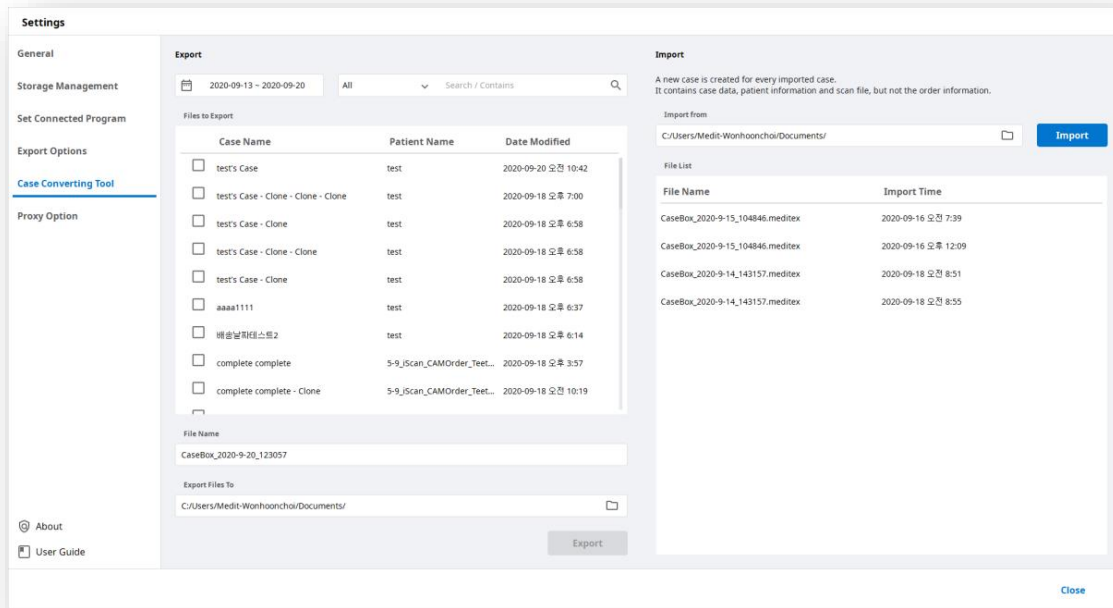


How to use

- ① Open Medit Link and press the “Settings” button located in the left bottom corner of the window.

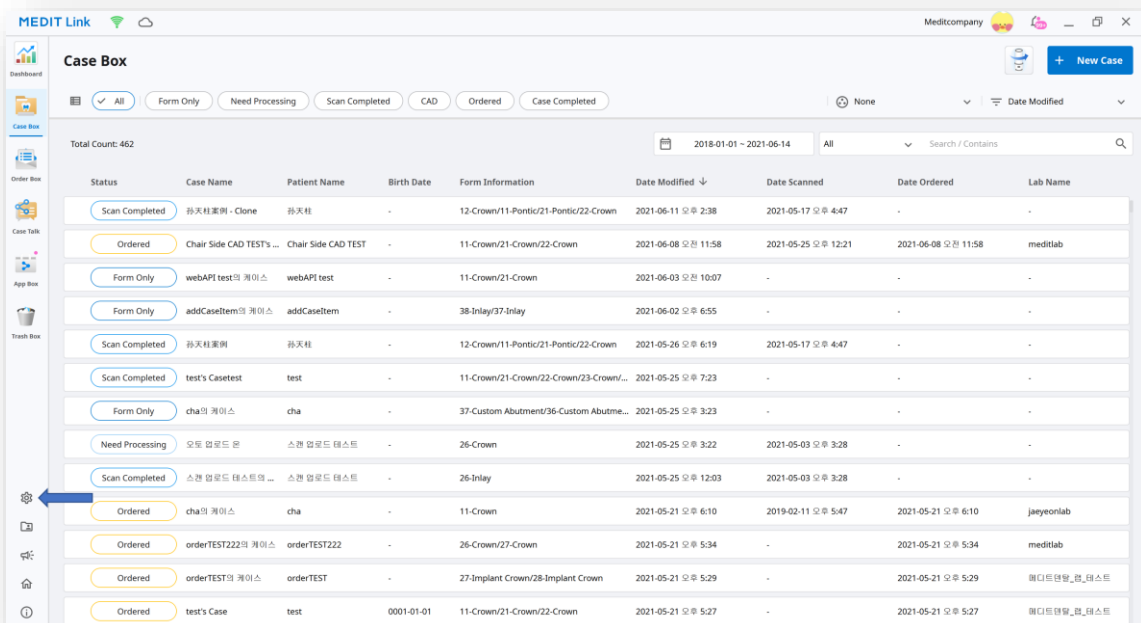


- ② Choose the “Case Converting Tool” option from the left side menu as shown on the picture below:

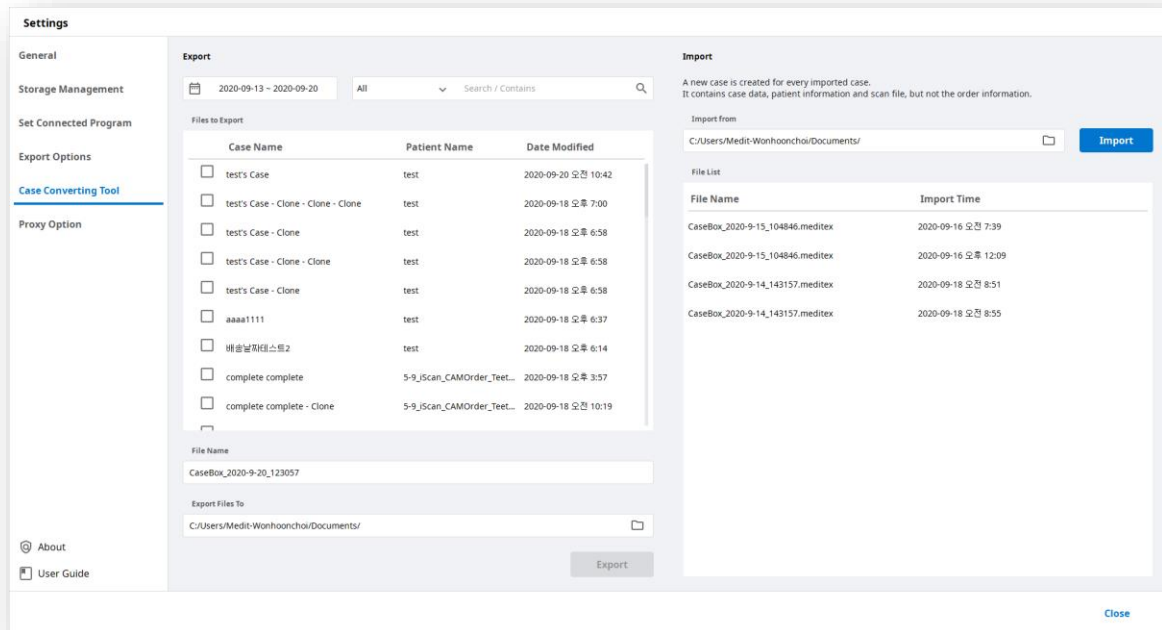


- ③ Click on the  icon to select the file to import from your PC.

- ④ Press “Import” button.

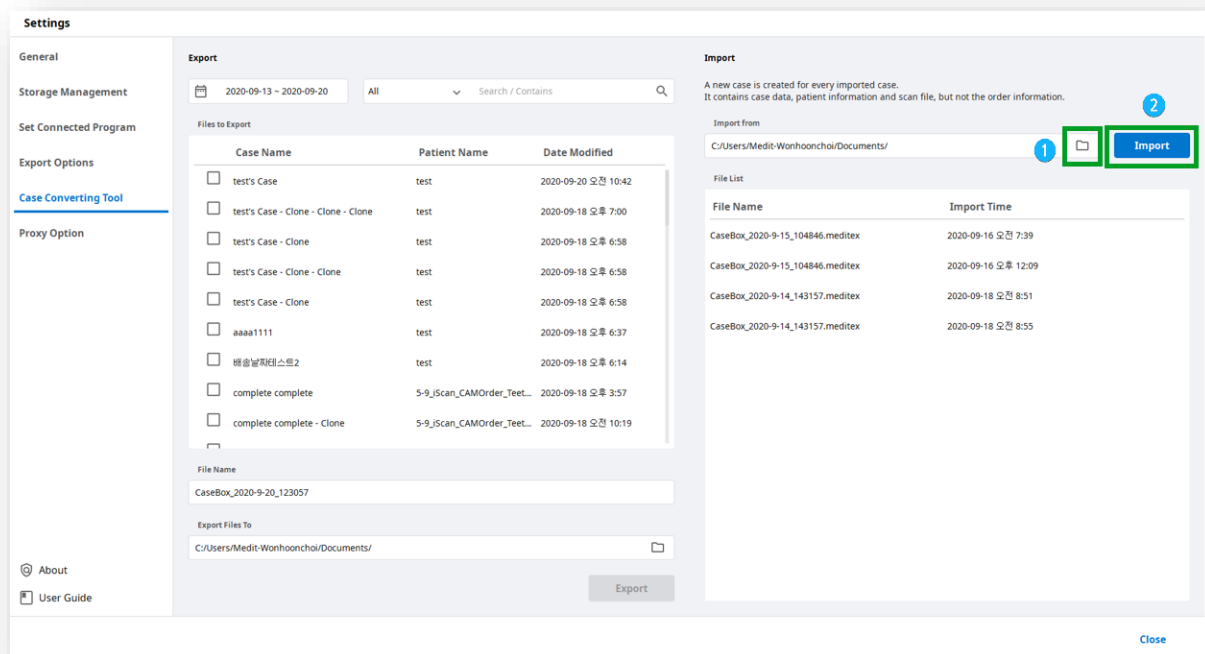


- ⑤ The case will appear in your Case Box (Clinic account) or Work Box (Lab account).

Proxy Options: In some limited network environments, you can change network settings.

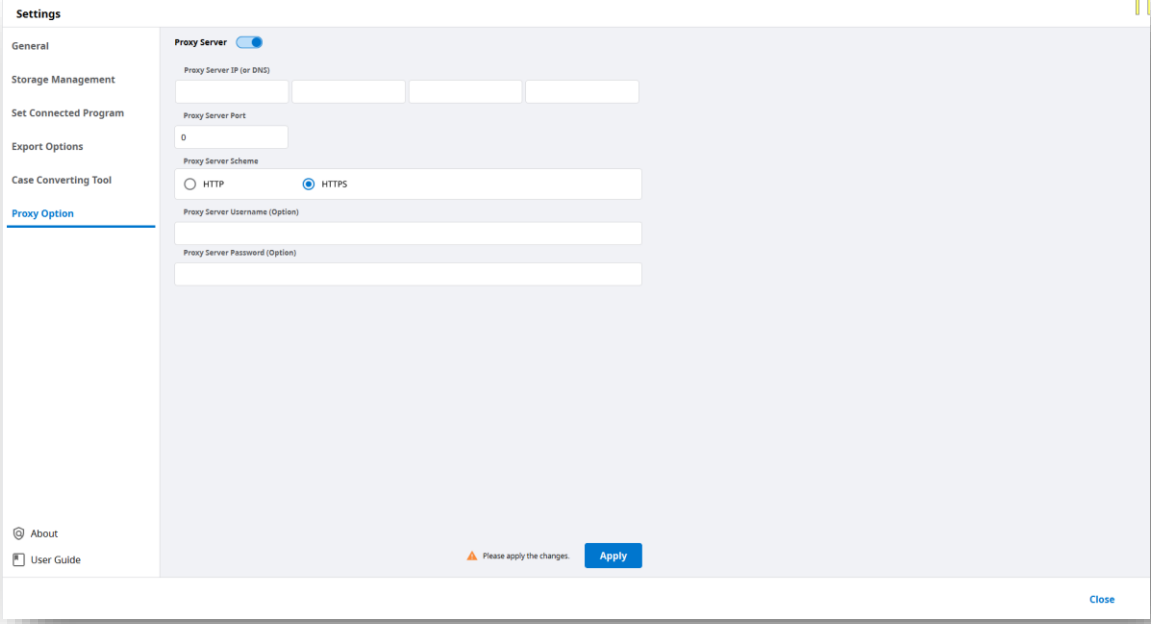

The screenshot shows the 'Settings' window with the 'Proxy Option' tab selected. The left sidebar lists various settings categories: General, Storage Management, Set Connected Program, Export Options, Case Converting Tool, and Proxy Option (which is highlighted). The main area is titled 'Proxy Server' and has a toggle switch turned on. Below the toggle are several input fields: 'Proxy Server IP (or DNS)' with three separate boxes, 'Proxy Server Port' with a box containing '0', 'Proxy Server Scheme' with radio buttons for 'HTTP' and 'HTTPS' (where 'HTTPS' is selected), 'Proxy Server Username (Option)' with an empty box, and 'Proxy Server Password (Option)' with an empty box. At the bottom of the main area, there is a warning icon and the text 'Please apply the changes.' followed by an 'Apply' button. The bottom of the window has a 'Close' button on the right and 'About' and 'User Guide' links on the left.

About: Displays the details of Medit Link.

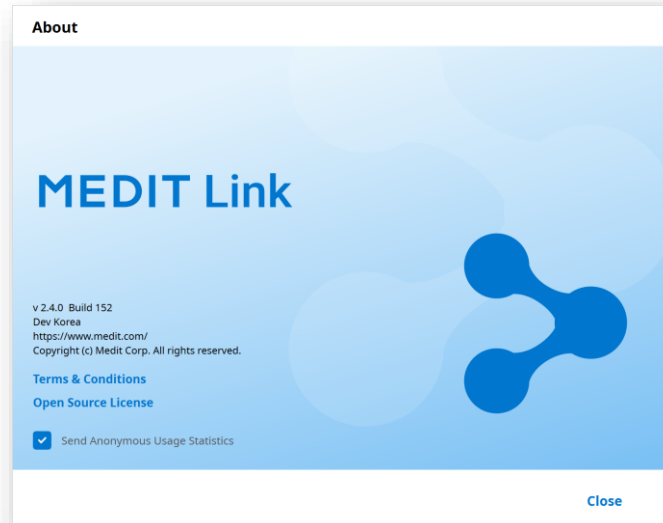

User Guide: Gets the user guide for Medit Link. Requires an Internet connection.

4.3.1 Patient Management

Provides the ability to add/edit/delete and merge the existing patients.

4.3.2 Notices and Updates

Displays program updates and critical alerts.

4.3.3 My Account


Transfers you to Medit Link homepage under the account that you are currently logged in with.

4.3.4 Help Center

Goes to the Help Center page. You can get various information about using Medit Link.

4.4 Detailed Functions of Case Box/Work Box

Case Box/Work Box provides the function to create cases and filter them using various options.

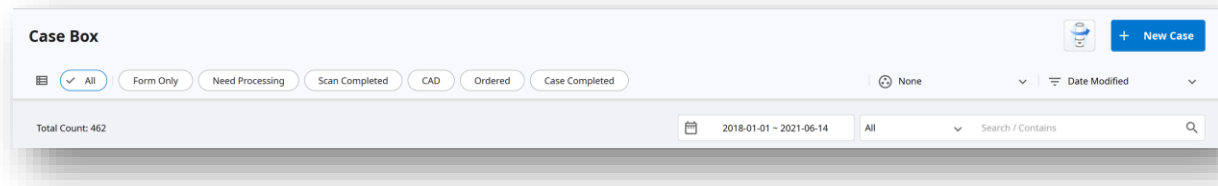


4.4.1 Creating cases

Create a case by pressing the “New Case” button in the top right corner and entering information such as patient and case name.



4.4.2 Change viewing options (Clinic Account only)

Change how you view cases between gallery and list forms.

Gallery

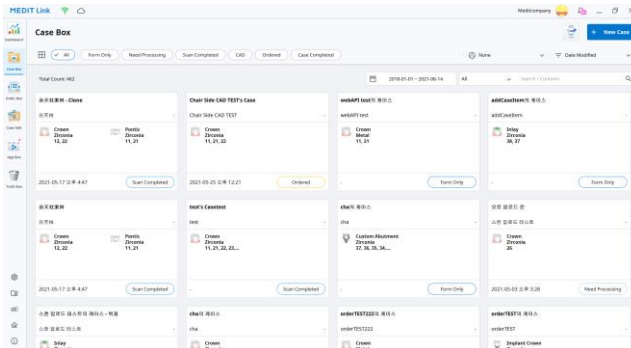

List

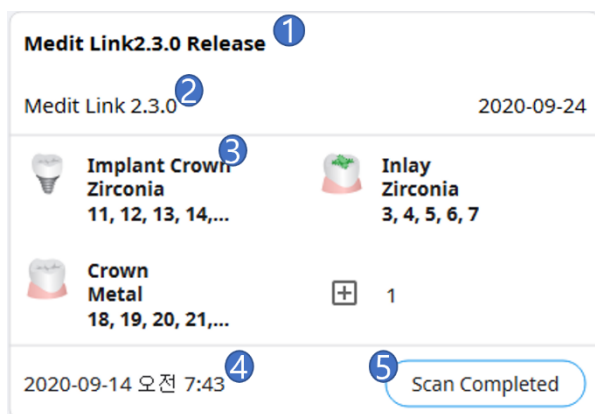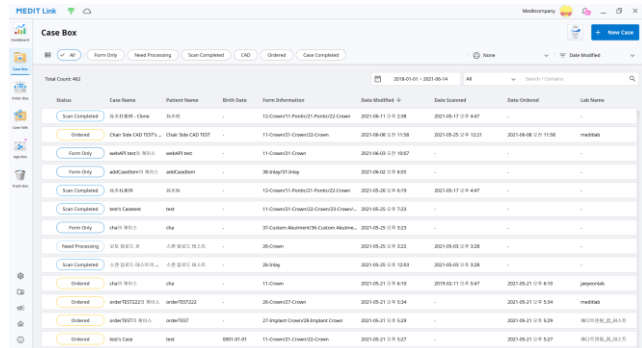


- 1 Case name
- 2 Patient name and birth date
- 3 Form information
- 4 Scanning Date
- 5 Case status

Filtering Options of Case List

Filtering using status	Quickly view/show cases of the selected status.
Filtering by groups	View cases by groups.
Filtering using period settings	Displays cases within a specified time period.
Filtering by search	Use a variety of search criteria to quickly locate a case.

5 Medit Link App – Box Detail

5.1 Case Box / Work Box

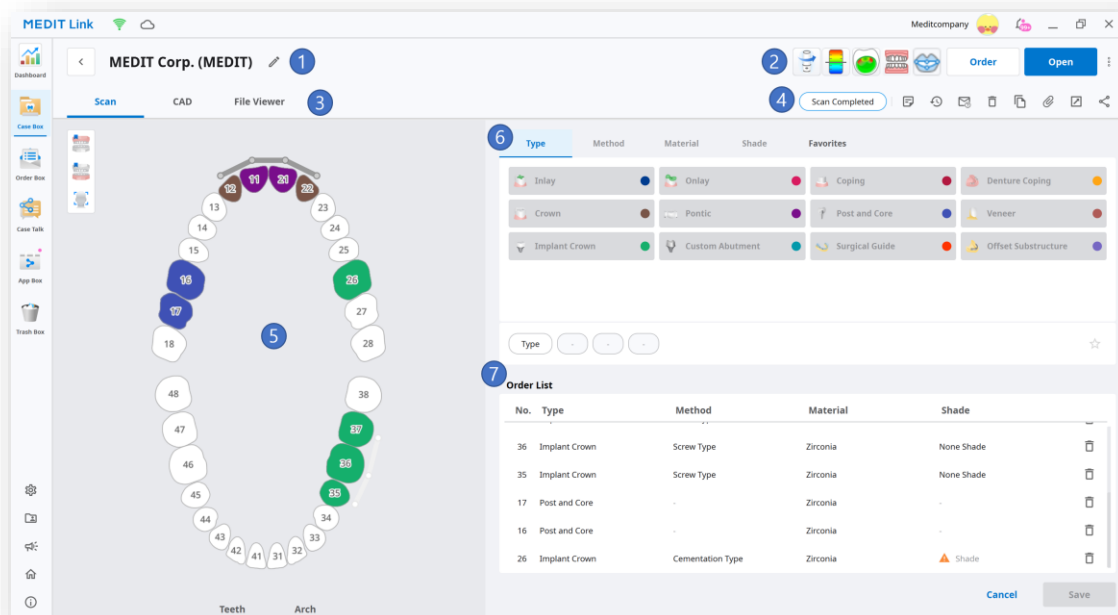

- 1 Displays the case name.
- 2 Displays the key available features: scanning, CAD, ordering, Apps.
- 3 Tabs for each of the tasks.
- 4 Functions for managing cases and files (case deletion, sharing, file attachment, export).
- 5 Displays corresponding tab for inputting dental form that can be used in selected scan or CAD program.
- 6 Allows to select type/method/material/shade information for form based on the CAD settings.
- 7 Displays the type/method/material/shade information applied to the selected teeth.



About Tabs

Scan	CAD	File Viewer
Enter information for basic scan tasks (form and scanning options). The information is passed to the connected scan program and used as the order information.	Set CAD options. The information you set up is passed to the CAD program.	View result files generated by scanning programs and CAD. Provides the ability to effectively display and manage 3D and 2D data, both acquired and attached.

About detailed functions of the case

	Status	Displays the status of the case.
	Go to Detail	Provides a shortcut for the same case in another box (Case Box ↔ Order Box; In Box ↔ Work Box).
	Memo	Shows the notes for the case. The memos will be included during ordering and can be used to communicate more detailed information about the case.
	Case History	Provides detailed information about the case. History Date Modified 2020-09-21 오후 5:29 Date Scanned Date Created 2020-09-21 오후 3:02 Case Size OMB Close
	Support Request	Asks the Medit Technical Support Center for technical support.
	Delete	Deletes the case (only for the cases created in Online Mode). You cannot delete cases that have been ordered.
	Clone	Clones a case to the Case Box. Useful for replicating treatment information for the same patient, or for additional orders after replicating all information/files in the same case.
	Attach	Attaches files in various formats.
	Export	Exports scan/CAD result files in the specified file format. Exported files can be used for future processes such as CAD, CAM, etc.
	Share	Create a link for sharing. View and check the data through a browser, without installing additional programs. Share the link freely or to specific people exclusively. Set a PIN code to access the link through the Medit Link website.

5.2 Order Box / In Box

2.3.6 Medit Scan Clinics_Basic (2.3.6 Data)

Order | **File Viewer**

Order ID: 3706292

Patient Name: 2.3.6 Data | **Partner Name:** Wonhoonchoi_Lab_Dev

Date Scanned: 2020-10-27 02:52 | **Date Ordered:** 2021-03-11 오전 8:39

Requested Delivery Date: 2021-03-18 | **Expected Delivery Date:** 2021-03-18

Memo:

```
<Scanbody Library>
36 : nt trading Preform L-Series 1.3 - RC - default
16 : nt trading Preform L-Series 1.3 - RC - default
14 : DESS - 3i OS TI Base - Desktop Scan Abutment - RP - 1.5mm - Engaging
```

No.	Type	Method	Material	Shade	Price
16	Implant Crown	Screw Type	Zirconia	A1	0.00
14	Custom Abutment	-	Zirconia	-	0.00
12	Veneer	Laminate	Press Ceramic	A1	0.00
11	Veneer	Laminate	Press Ceramic	A1	0.00
23	Custom Abutment	-	Zirconia	-	0.00
25	Coping	PFZ	Zirconia	A1	0.00
26	Pontic	Reduced	Zirconia	A1	0.00
27	Coping	PFZ	Zirconia	A1	0.00
36	Implant Crown	Screw Type	Zirconia	A1	0.00
35	Custom Abutment	-	Zirconia	-	0.00
46	Inlay	-	Zirconia	A1	0.00

Total USD 0.00

- 1 Displays the case name and patient name.
- 2 Provides functions for order management (accepting/cancelling; delivering/receiving).
- 3 Tabs for each of the tasks.
- 4 Shows case status and shortcuts.
 - Go to Detail: A shortcut to check the details of ordered cases. Click the button to go to the Case Box or Work Box case.
 - Status: Displays the job status of the case.
- 5 Displays order number, patient name, account name, scanning Date, order date, memo, etc.
- 6 Displays the product information and delivery information.



About Tabs

Order

Check the order number, patient name, account name, scanning date, order date, memo, and other detailed information related to the order.


File Viewer

View files generated by scan/CAD operations. It effectively displays and manages 3D models and 2D images, both acquired or attached.

5.3 Case Talk

Cases created in the Medit Link can be shared with restrictions, without restrictions, with specific people, and registered members of Medit Link. You can see all shared lists through “Case Talk” in the Medit Link application and Medit Link Web viewer.


“Case Talk” service requires an internet connection.

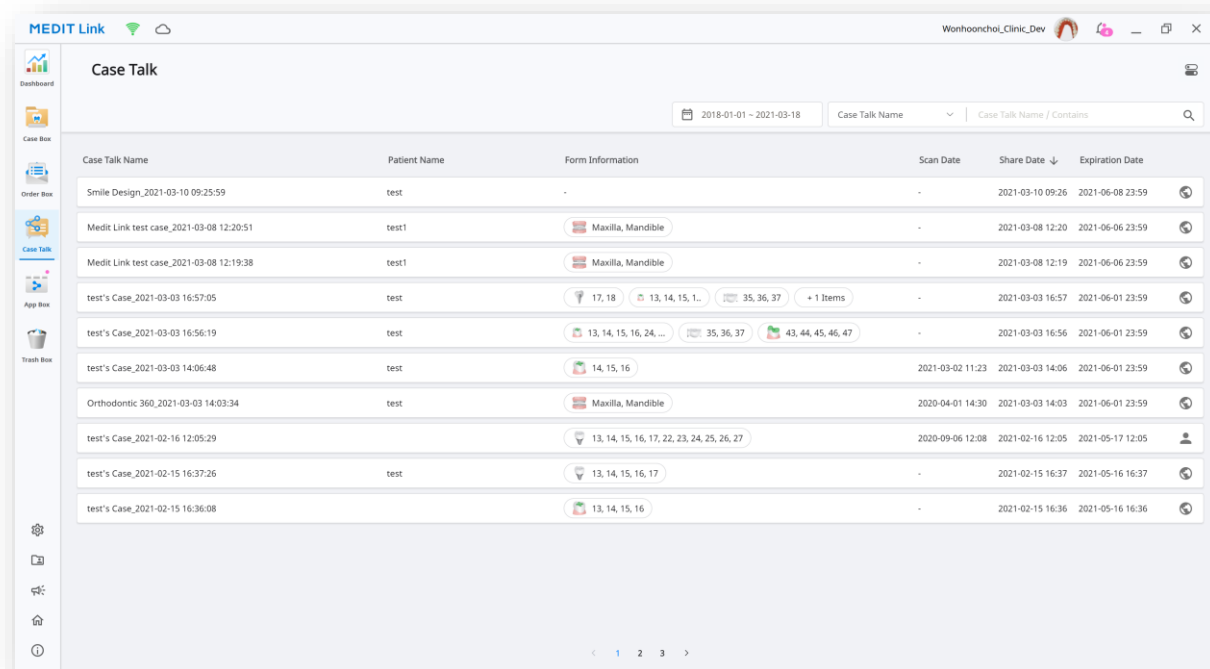

The screenshot displays the 'Case Talk' section of the Medit Link application. The interface includes a sidebar with navigation icons for Dashboard, Case Box, Order Box, Case Talk (active), App Box, and Trash Box. The main area shows a table of shared cases with columns for Case Talk Name, Patient Name, Form Information, Scan Date, Share Date, and Expiration Date. A search bar and date range filter (2018-01-01 ~ 2021-03-18) are at the top right. The table lists 11 cases, including 'Smile Design', 'Medit Link test case', and 'test's Case', with various form information and dates.

Case Talk Name	Patient Name	Form Information	Scan Date	Share Date ↓	Expiration Date
Smile Design_2021-03-10 09:25:59	test	-	-	2021-03-10 09:26	2021-06-08 23:59
Medit Link test case_2021-03-08 12:20:51	test1	Maxilla, Mandible	-	2021-03-08 12:20	2021-06-06 23:59
Medit Link test case_2021-03-08 12:19:38	test1	Maxilla, Mandible	-	2021-03-08 12:19	2021-06-06 23:59
test's Case_2021-03-03 16:57:05	test	17, 18, 13, 14, 15, 1, 35, 36, 37, + 1 Items	-	2021-03-03 16:57	2021-06-01 23:59
test's Case_2021-03-03 16:56:19	test	13, 14, 15, 16, 24, 35, 36, 37, 43, 44, 45, 46, 47	-	2021-03-03 16:56	2021-06-01 23:59
test's Case_2021-03-03 14:06:48	test	14, 15, 16	2021-03-02 11:23	2021-03-03 14:06	2021-06-01 23:59
Orthodontic 360_2021-03-03 14:03:34	test	Maxilla, Mandible	2020-04-01 14:30	2021-03-03 14:03	2021-06-01 23:59
test's Case_2021-02-16 12:05:29	test	13, 14, 15, 16, 17, 22, 23, 24, 25, 26, 27	2020-09-06 12:08	2021-02-16 12:05	2021-05-17 12:05
test's Case_2021-02-15 16:37:26	test	13, 14, 15, 16, 17	-	2021-02-15 16:37	2021-05-16 16:37
test's Case_2021-02-15 16:36:08	test	13, 14, 15, 16	-	2021-02-15 16:36	2021-05-16 16:36

Detail information in Case Talk

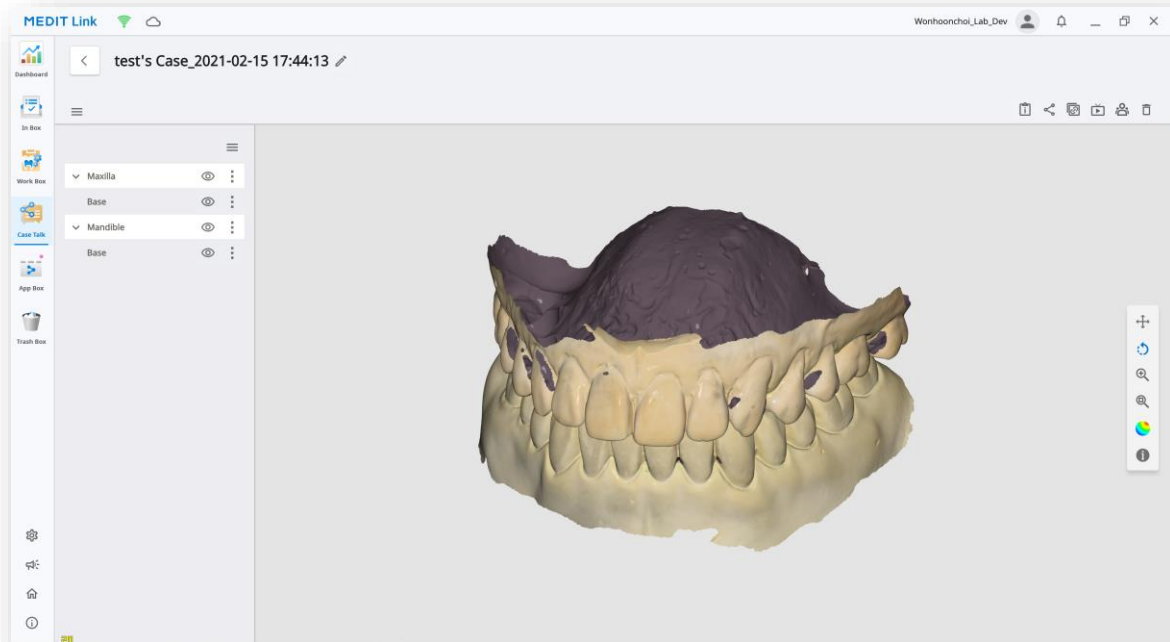

Case Information Shows detailed information of the case.



Share Modifies share settings.



Copy link Copies the shared link to the clip board. You can see the shared item through the web browser using the copied link.



Live Shares your screen with connected users. Screen sharing is synchronized with the connected users.



Attendees Shows all users who are connected to this item.




Delete Deletes shared link. All shared items will be deleted and people with access to the link will no longer be able to view it.

5.4 App Box

In addition to the built-in programs, we offer a variety of additional ones.

Download the programs of your liking and use them freely.

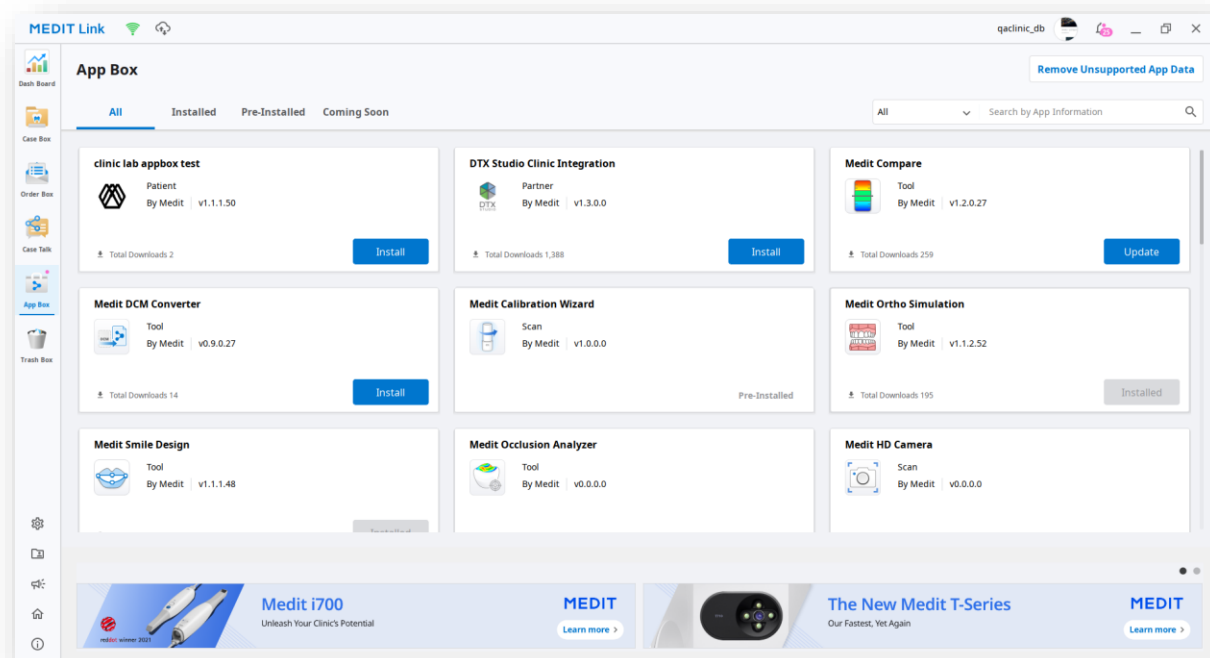

5.5 Trash Box

Trash Box is a space for keeping deleted cases before they are completely removed or recovered.



Data deleted from the Trash Box is deleted from the PC and server and cannot be recovered.

Status	Case Name	Patient Name	Lab Name	Date Deleted	Deleted By	Size
☐ Form Only	Multi Occlusion	test		2021-03-05 2:51	Wonhoonchoi_Clinic	0 KB
☐ Form Only	test's Case	test		2021-03-05 2:51	Wonhoonchoi_Clinic	0 KB
☐ Form Only	test's Case - Clone	test		2021-03-05 2:51	Wonhoonchoi_Clinic	0 KB
☐ Scan Completed	test's Case - Clone	test		2021-03-05 2:51	Wonhoonchoi_Clinic	21,325 KB
☐ Scan Completed	Multi Occlusion Sample - Clone	test		2021-03-05 2:51	Wonhoonchoi_Clinic	1,059,768 KB
☐ Scan Completed	test's Case	test		2021-03-04 8:24	Wonhoonchoi_Clinic	50,862 KB
☐ Scan Completed	test's Case	test		2021-03-04 8:24	Wonhoonchoi_Clinic	13,411 KB
☐ Scan Completed	test's Case	test		2021-03-04 8:24	Wonhoonchoi_Clinic	14,069 KB
☐ Scan Completed	test's Case	test		2021-03-04 8:24	Wonhoonchoi_Clinic	30,765 KB
☐ Scan Completed	test's Case	test		2021-03-04 8:24	Wonhoonchoi_Clinic	13,061 KB
☐ Scan Completed	test's Case	test		2021-03-04 8:24	Wonhoonchoi_Clinic	30,133 KB
☐ Form Only	test's Case	test		2021-03-02 1:52	Wonhoonchoi_Clinic	0 KB


6 Detailed Explanation of the Process

6.1 Managing Cases: Creating, Editing, Cloning, Deleting

6.1.1 Creating Cases

Creating cases in Case Box (Clinic)

- ① Click “New Case” button in the upper right corner.

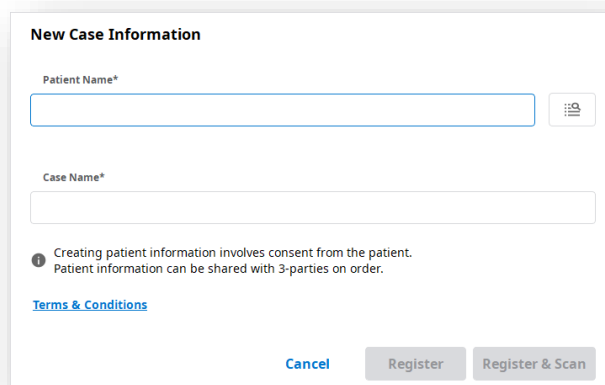


The 'New Case Information' form contains two text input fields: 'Patient Name*' and 'Case Name*'. To the right of the 'Patient Name*' field is a button with a magnifying glass icon. Below the 'Case Name*' field is a small information icon followed by the text: 'Creating patient information involves consent from the patient. Patient information can be shared with 3-parties on order.' Below this is a link labeled 'Terms & Conditions'. At the bottom of the form are three buttons: 'Cancel', 'Register', and 'Register & Scan'.

- ② Enter basic information such as patient name.

- Entering a name that is not in the patient list creates a new patient. You can enter a name in the patient list to select from the list.

- ③ In case you are registering a case for an already existing patients, click the  button.

- You can select a new patient after registering or editing information for an existing patient.

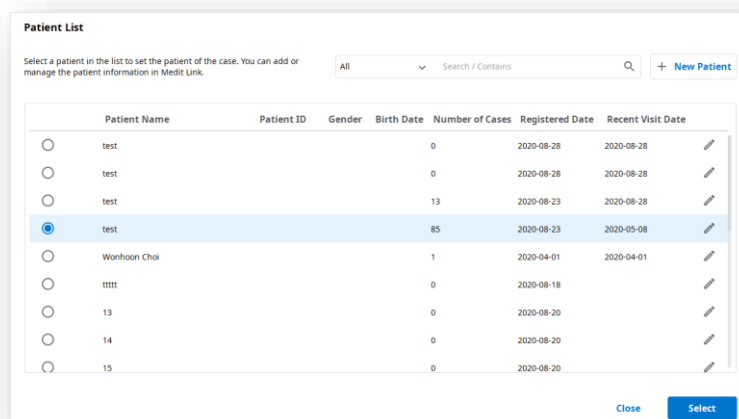

The 'Patient List' table displays a list of patients with columns for Patient Name, Patient ID, Gender, Birth Date, Number of Cases, Registered Date, and Recent Visit Date. Each row has a radio button on the left and an edit icon on the right. The fourth row, with Patient Name 'test' and Number of Cases '85', is selected. Below the table are 'Close' and 'Select' buttons.


	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Registered Date	Recent Visit Date	
☐	test				0	2020-08-28	2020-08-28	
☐	test				0	2020-08-28	2020-08-28	
☐	test				13	2020-08-23	2020-08-28	
☒	test				85	2020-08-23	2020-05-08	
☐	Wonhoon Choi				1	2020-04-01	2020-04-01	
☐	ttttt				0	2020-08-18		
☐	13				0	2020-08-20		
☐	14				0	2020-08-20		
☐	15				0	2020-08-20		

- Case name will be filled in automatically, however you can still change it.

- ④ Press “Register & Scan” to complete the case registration and proceed to scanning immediately.
- Otherwise, press “Register” which is useful when registering multiple cases at the same time.
 - You will be redirected to the main page of the Case Box.

Creating cases in Work Box (Lab)

- ① Click “New Case” button in the upper right corner.

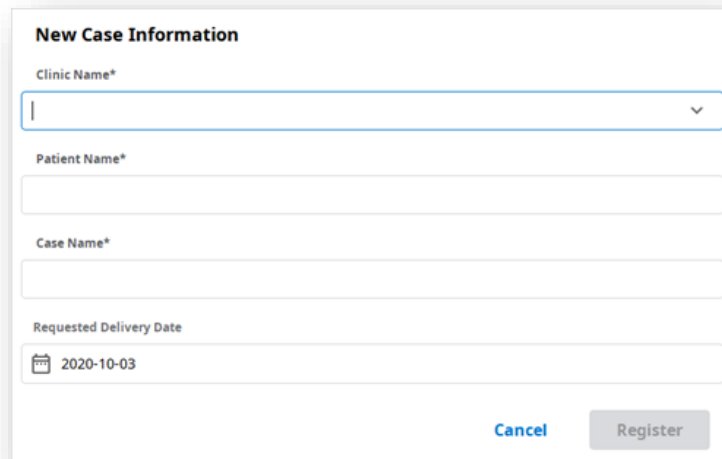


The screenshot shows a modal form titled "New Case Information". It contains four input fields: "Clinic Name*" (a dropdown menu), "Patient Name*" (a text box), "Case Name*" (a text box), and "Requested Delivery Date" (a date picker showing "2020-10-03"). At the bottom right, there are two buttons: "Cancel" and "Register".

- ② Enter or select the clinic name.
- ③ Enter the patient name and case name.
- Case name will be filled in automatically, however you can still change it.
- ④ Enter the expected delivery date considering the completion date of the work.
- ⑤ After registering the case, go to the job information entry screen.

6.1.2 Editing cases

Editing cases in Case Box (Clinic)

- ① Click on the  icon next to the case name.

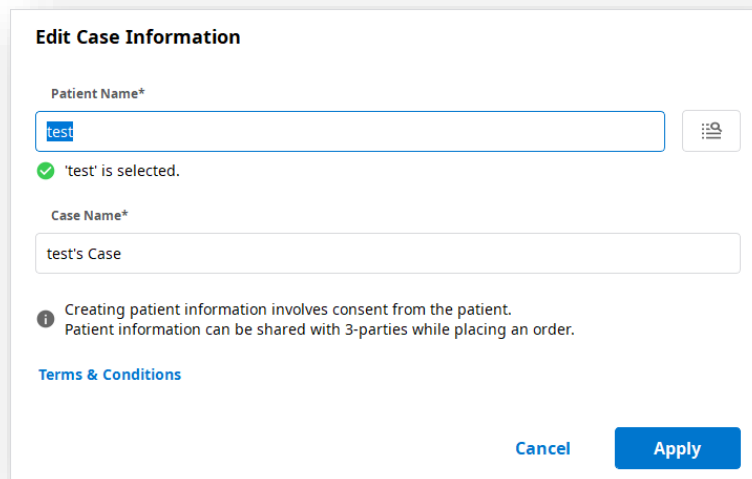

Edit Case Information

Patient Name*

test

✔ 'test' is selected.

Case Name*


test's Case

Creating patient information involves consent from the patient.
Patient information can be shared with 3-parties while placing an order.

[Terms & Conditions](#)

Cancel Apply

- ② Click on the  icon to change the patient name. Note that it cannot be edited directly.; instead, you can find a patient name you would like to replace the current one with and can edit the patient information from the list.

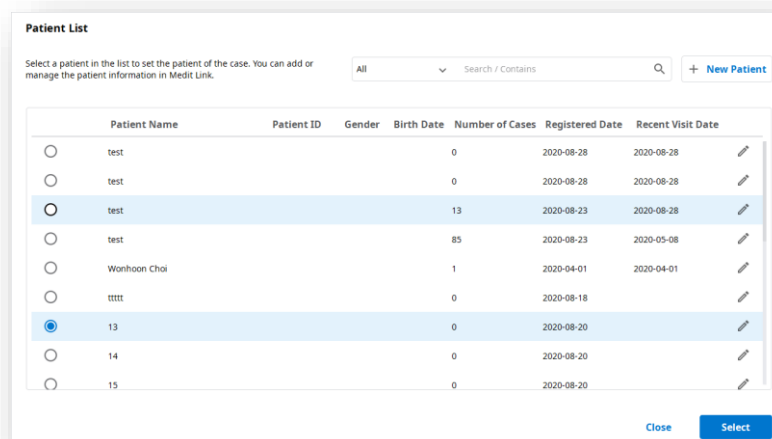

Patient List

Select a patient in the list to set the patient of the case. You can add or manage the patient information in Medit Link.

All Search / Contains + New Patient


	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Registered Date	Recent Visit Date	
☐	test				0	2020-08-28	2020-08-28	
☐	test				0	2020-08-28	2020-08-28	
☒	test				13	2020-08-23	2020-08-28	
☐	test				85	2020-08-23	2020-05-08	
☐	Wonhoon Choi				1	2020-04-01	2020-04-01	
☐	ttttt				0	2020-08-18		
☒	13				0	2020-08-20		
☐	14				0	2020-08-20		
☐	15				0	2020-08-20		

Close Select

- ③ Check the patient name you would like to replace the current one with.
- ④ To create a new patient, press the “New Patient” button.
- ⑤ Press  to edit the registered information about the patient.
- ⑥ Press “Apply” to finish the process.

Editing cases in Work Box (Lab)

- ① Click on the  icon next to the case name.

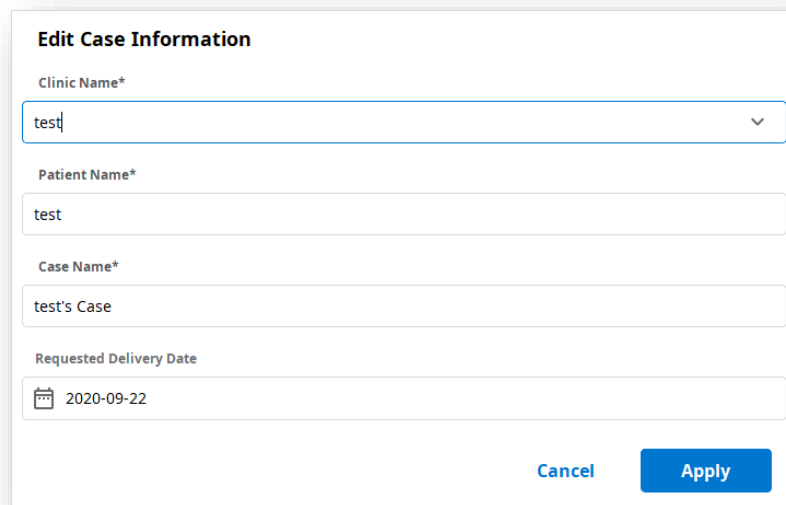

The dialog box is titled "Edit Case Information". It contains four input fields and two buttons at the bottom right. The first field, "Clinic Name*", is a dropdown menu with "test" selected. The second field, "Patient Name*", is a text box with "test". The third field, "Case Name*", is a text box with "test's Case". The fourth field, "Requested Delivery Date", is a date picker showing "2020-09-22". The buttons are "Cancel" and "Apply".

Field Label	Value
Clinic Name*	test
Patient Name*	test
Case Name*	test's Case
Requested Delivery Date	2020-09-22

- ② Press "Apply" to after the editing is done.




You cannot edit the information of the cases received through the order in In Box.

6.1.3 Cloning cases

You can use the Clone feature to clone a case to rescan with the same treatment information. You can duplicate the entire data to order a case that has already been ordered or choose other cloning options.

Cloning cases in Case Box (Clinic)

- ① Click on the  icon.
- ② You will be asked what kind of information you would like to clone to the new case.

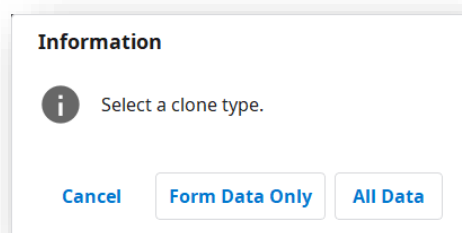

- ③ Choose whether to replicate only treatment information or data all together.
- ④ When cloning is complete, a case with the name "[Case name] – Clone" will be created.


Cloning cases in Work Box (Lab)

- ① Click on the  icon.
- ② Duplicate all data to create a new case.
- ③ When cloning is complete, a case with the name "[Case name] – Clone" will be created.

6.1.4 Deleting cases (Lab and Clinic)

You can delete cases.

- ① Click on the  icon.
- ② Click “OK” to delete the case.

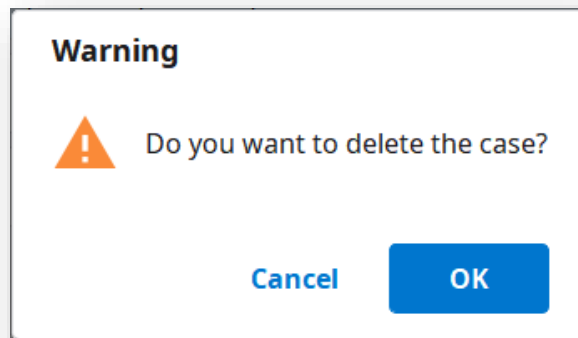

- ③ Deleted cases are moved to the Trash Box and can be recovered if necessary.

6.2 Inputting form information

You can enter treatment information based on the connected program of your selection.

The information you enter will be used as basic information for scanning, CAD, and ordering.

6.2.1 Inputting Form Information in Case Box

Entering dental information in the Scan tab

- ① Select whether you are going to input treatment information by teeth or by arch.
- ② Select Type – Method – Material – Shade.
 - If you entered information incorrectly, you go back to select it again.
 - The selected information is automatically registered to the order list.
- ③ Click “Save” to save your registration.

Entering dental information in the CAD tab

- ① Select one of the following options that will define your scanning strategy: Teeth/Flexible Multi-Die/Multi-Die.
- ② Select the tooth number/arch/die number.
- ③ Define the settings for detailed options and CAD operations.
- ④ Click “OK” to save the settings.
 - The selected information is automatically registered to the order list.
- ⑤ Click “Save” to save your registration.

6.2.2 Inputting Form Information in Work Box

Entering dental information on the Scan tab

- ① Select the scan program you want to use.
- ② Select one of the following options that will define your scanning strategy: Teeth/Flexible Multi-Die/Multi-Die.
- ③ Methods may differ depending on the selected scan program.

	Medit Scan for Clinics	Medit Scan for Labs
Teeth	O	O
Arch	O	O
Flexible Multi-Die	X	O
Multi-Die	X	O

- ④ Select teeth, arch, or die numbers for the scanning process.
- ⑤ Select Type – Method – Material – Shade.
 - If you entered information incorrectly, you go back to select it again.
 - The selected information is automatically registered to the list.
- ⑥ Click “Save” to save your registration.


Entering dental information in the CAD tab

- ① Select one of the following options that will define your scanning strategy: Teeth / Flexible Multi-Die / Multi-Die.
- ② Select the tooth number/arch/die number.
- ③ Define the settings for detailed options and CAD operations.
- ④ Click “OK” to save the settings.
 - The selected information is automatically registered to the list.
- ⑤ Click “Save.”


6.2.3 Inputting Bridge Information

You can enter information for a bridge in case same information is applied to all selected teeth when selecting several teeth of the same type a bridge icon will appear by their side. You can turn the bridge on/off by clicking the icon.

Before setting a bridge

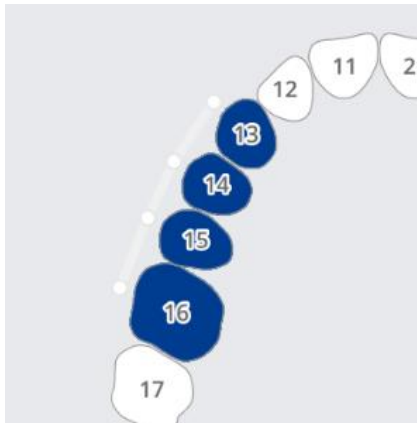

After setting a bridge

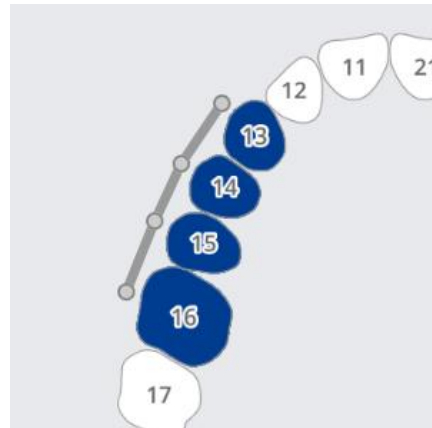

6.2.4 Setting Scanning Options

Depending on the program you are working with, various scan options are provided.

The selected scan options are passed to each scan program.

Medit Scan for Clinics

	Pre-op for Maxilla	When checking this option, you will be able to scan maxillary diagnostic (pre-operation) model.
	Pre-op for Mandible	When checking this option, you will be able to scan mandibular diagnostic (pre-operation) model.
	Face Scan	When checking this option, you will be able to acquire face scan data such as teeth, lips, and nose with the intraoral scanner.

Medit Scan for Labs

General options

	Adjacent teeth	Mark the prepared and adjacent teeth and pass this information to the scanning program.
	Antagonist	Automatically marks the antagonist and passes this information to the scanning program.
	Articulator	Displays options for selecting the articulator in the scanning program.
	Movement Marker	Allows to use the Jaw Motion scan strategy in the scanning program.

Options by teeth

	Scan body	Allows to set the strategy for scanning scan body in the scanning program.
	Pre-op Scan	Allows to set the strategy for scanning pre-op model in the scanning program.
	Wax-up Scan	Allows to set the strategy for scanning wax-ups in the scanning program.

6.3 Copying Form Information

You can copy the treatment information

Copy dental information using context menu

- ① Select the tooth for which the tooth information is set.
- ② Right-click and select “Copy” from the context menu.
- ③ Select the target tooth.
- ④ Click  to complete the copy.

Copy dental information using shortcut keys

- ① Select the tooth for which the tooth information is set.
- ② Press and hold CTRL and select the target tooth.

6.4 Deleting treatment information

You can delete the dental information by clicking  in the order list.

No.	Type	Method	Material	Shade	
16	Inlay	-	Zirconia	A4	
25	Inlay	-	Gold	-	
26	Onlay	-	Gold	-	

Delete dental information using context menu

- ① Select the tooth for which the tooth information is set.
- ② Right-click to select “Delete” from the context menu.
- ③ Select the target tooth.
- ④ Click  to complete the delete.

Delete dental information using shortcut keys

- ① Select the tooth for which the tooth information is set.
- ② You can delete the selected information by pressing the keyboard “Del” button.

6.5 Adding Product Information to “Favorite”

You can quickly select a product by registering it as a favorite.

The list of registered favorites is managed as a template and can be shared through the export and import features.

6.5.1 Adding to Favorites

- ① Select the tooth number.
- ② Select Type – Method – Material – Shade.
 - Type – Method – Material – Shades can be registered as favorite even if not all of them are selected.
- ③ Click  .
- ④ A combination of the selected parameters, excluding tooth number, is listed on your favorite search list.

6.5.2 Managing Favorite Templates

You can register, manage, and share different types of templates.


Add	Adds a new template.
Rename	Renames the current template.
Delete	Deletes the current template.
Export	Exports the current template to a file.
Import	Imports a template file.

6.6 Copying form information between Scan and CAD tabs

You can copy the treatment information entered between Scan tab and CAD tab.

6.6.1 From Scan tab to CAD tab

- ① If there is no information input in the CAD Tab, it automatically copies the information when you go from the Scan Tab to the CAD Tab.
- ② If there is information in the CAD tab, you will have to copy it manually.
 - Right-click on the CAD tab.
 - Select “Convert Form Information from Scan” from the context menu.

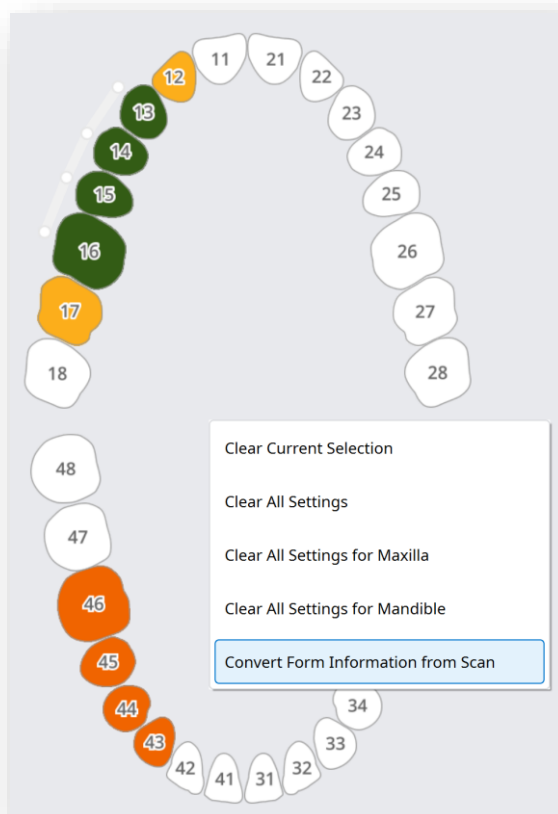


- The information in the Scan tab is copied to the CAD tab.

6.6.2 From CAD tab to Scan tab

To copy information, use the following methods:

- ① Right-click on the Scan tab.
- ② Select “Convert Form Information from CAD” from the context menu.

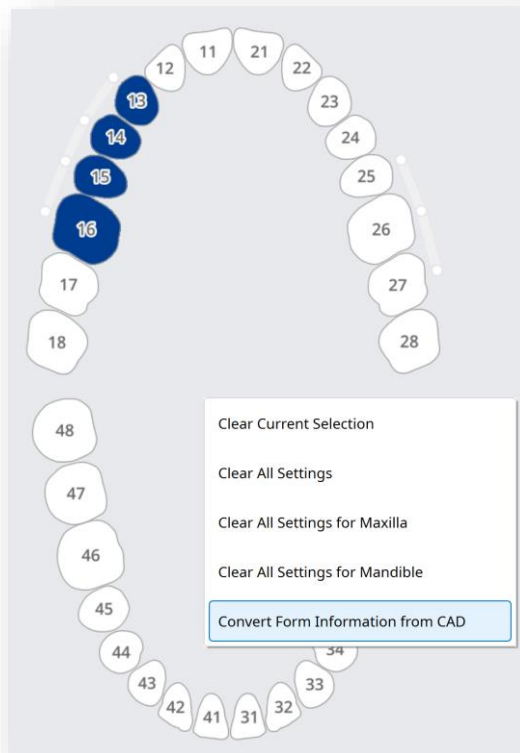

Note that only Type-Method information will be copied.

6.7 Scanning process

Based on the entered information, you proceed with scanning.

- ① Enter/save treatment information and select the necessary scan options.
- ② Click “**Scan**” button in the upper-right corner of the window.

6.8 CAD Process

Medit Link is integrated with exocad.

CAD	Registering treatment information	Confirmation of treatment information	Direct CAD Execution	CAD Results import back to Medit Link
exocad	0	0	0	0

- ① In Settings > Set Connected Program, you can link the CAD program of your choice from the list.
- ② Enter/save task information and design options in the CAD Tab.
- ③ Click “CAD” in the upper-right corner of the window.
- ④ The linked CAD program is run, and the data is sent automatically.
- ⑤ Save the results after performing CAD.

6.9 Exporting Result Files

In addition to sharing results with partners through orders, Medit Link allows you to store your results in the desired file format using the export function.

Export files in case

Options File List

Export Files To
C:/Users/Medit-Wonhoonchoi/Desktop/

Folder Name
2021-03-03-test

File Name
2021-03-03-test

CAD Type
exocad

Mesh Format
☐ Keep Original Mesh Format
☒ Mesh Format
☒ OBJ ☐ PLY ☐ STL

Image Format
☐ Keep Original Image Format
☒ Image Format
☒ JPG ☐ JPEG ☐ PNG ☐ BMP

☐ Compress Files
☐ Combine Individual Mesh
☒ Export All Occlusal Relationship

Close Export

Export Files To	Specifies the path where you want to store the results.
Folder Name	The system creates a new folder and stores all results in it. The default folder name is set by naming rule. But you can change the folder name directly. You can define the naming rule through Settings > Export Options > Rules for Folder Name
File Name	Sets the file name for the result. The default file name is set by naming rule. But you can change the file name directly. You can define the naming rule through Settings > Export Options > Rules for File Name
CAD Type	Select the CAD program to be used for your design program. Medit Link will adjust the files based on this selection, allowing to proceed with the design without any additional axis alignment.
Mesh Format	Specifies the format of the mesh file.
Image Format	Specifies the format of the exported image files.
Compress Files	Compresses all results into a single file.
Combine Individual Mesh	Combines the result of maxilla or mandible with individual tooth. Useful when importing several files into a program is not possible.
Export All Occlusal Relationship	When the multiple occlusal results generated from the scan app, you can export all occlusal results at once.

6.10 File Viewer

File Viewer groups and shows various data, including the data acquired from programs linked to Medit Link and data attached by users.



It provides basic management functions such as viewing/hiding data, adjusting transparency, and deleting data.

It supports previews of files in the following formats:

- 3D data: STL, OBJ, PLY
- 2D images: PNG, JPG, JPEG, BMP, GIF, TIF, TIFF
- Video files: MKV, MP4, AVI, WMV

6.10.1 Data Tree

It provides the ability to group and display data and perform actions such as viewing/hiding, transparency adjustment, and deletion through the data tree.

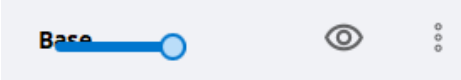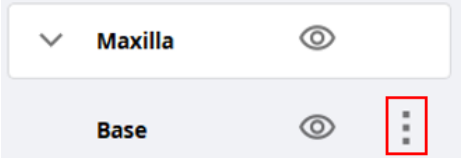	Show/Hide	You can adjust the view/hide status of your data.
	Transparency adjustment	You can adjust transparency for elements in the “View” state by moving the slider bar. 
	Extension Menu	Displays additional features available.  Show this only: Hides all other data and display only current data. Hide: Changes the current data to the Hide state. Duplicate: Copy the data and add it to the attachment group. Rename: Change the name of the data. Delete: Deletes the data.



Rename function is available only on online mode.

6.10.2 Model Control

3D Model Control

	Pan	Moves the model.
	Rotate	Rotates the model.
	Zoom in/out	Zooms in/out the model.
	Zoom Fit	Fits the model to the screen.
	Change Color Display Mode	Changes the color display mode for the model.
	Grid Settings	Show or hide the grid and control its position in relation to the model (overlay on/off).


2D Images Control

	Next/previous image	Moves among images.
	Edit	Edits the image.
	Pinch/Paste Window	You can detach the image window and adjust it to the desired position. You can resize the window. When the window is restored to its original location, it returns to its default size.

6.10.3 Capture and Memo


Capture the current screen as a 2D image, draw on it, or take notes.

You can save it as a file together with your notes and save them as a reference.

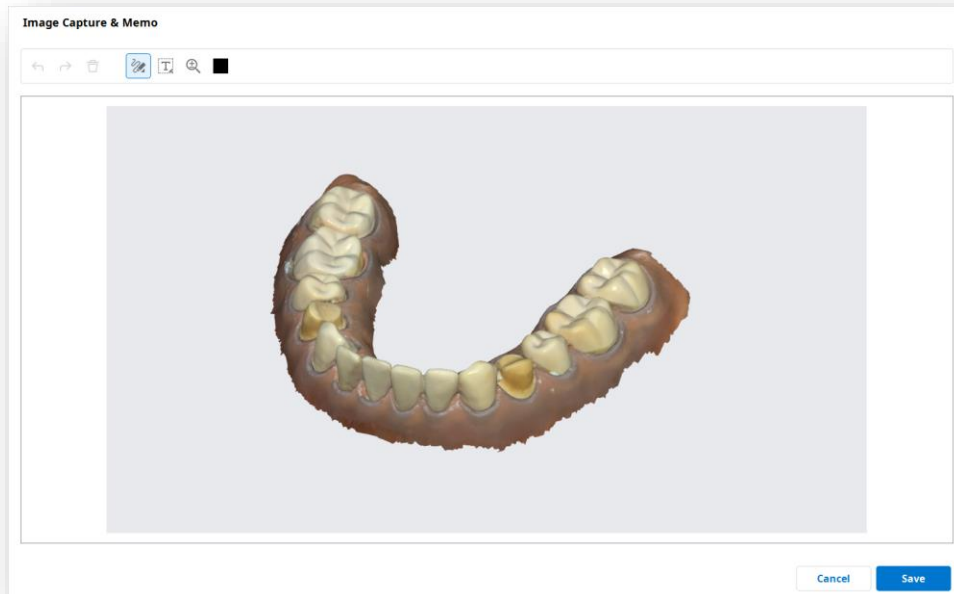


6.10.4 Orthodontic Capture

Captures the images in current viewing direction and top/bottom/left/right from the current viewing direction.

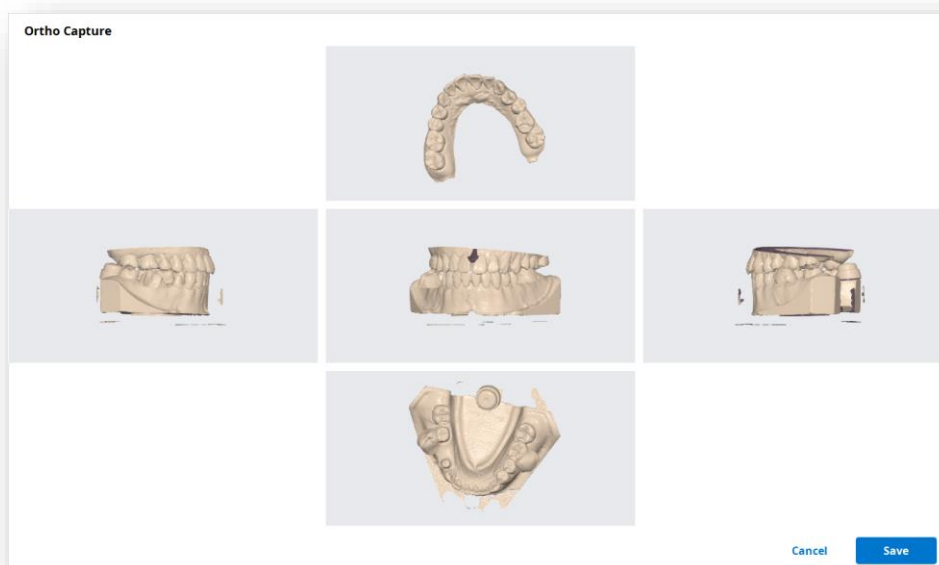

6.10.5 Fill Holes

Fill holes in 3D mesh data.

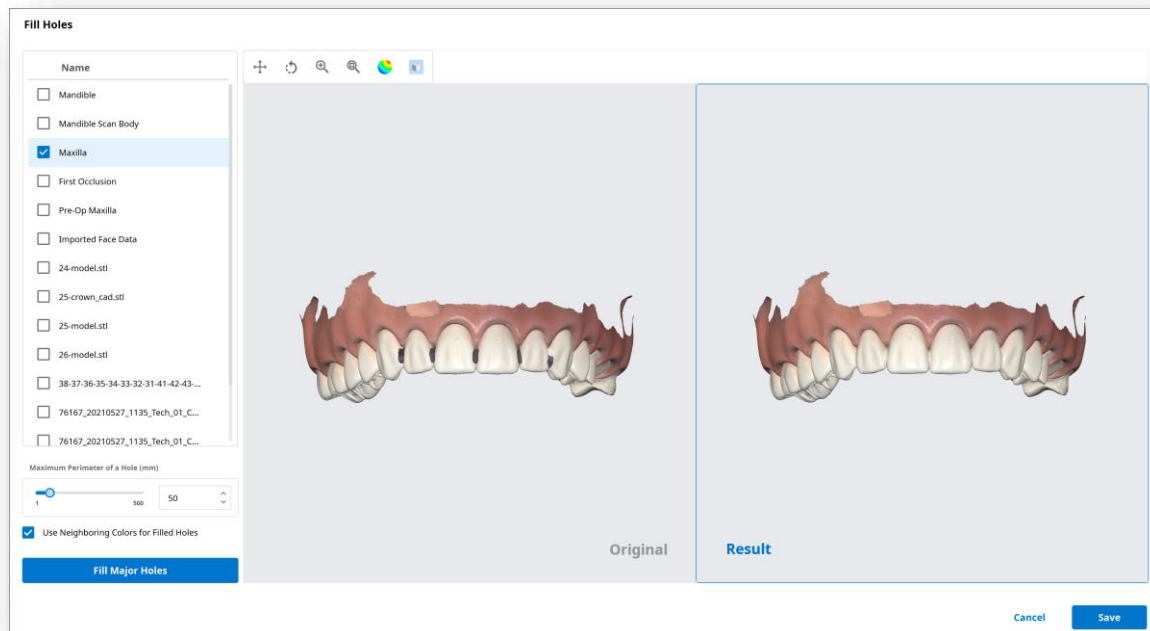

6.10.6 Trimming

Trimming tool provides various tools to edit data.

The selected 3D mesh data can be edited using different tools located on the bottom of the window.

Toolbox

Selection

	Smart Selection Teeth	 Automatically selects all teeth of the arch leaving out gingiva parts. This function is only available for the scan data that has been acquired by Medit Scan for Clinics with the “Use GPU” option on.
	Smart Single Tooth Selection	Automatically selects the area of a single tooth leaving out gingiva parts. Click, press, and drag the mouse on the tooth.
	Polyline Selection	Selects all entities within a polyline shape drawn on the screen.
	Brush Selection	Selects all entities on a freehand-drawn path on the screen. Only the front face will be selected. The brush comes in three different sizes.

After selecting an area

	Shrink Selected Area	Reduces the selected area each time you press the button.
	Expand Selected Area	Expands the selected area each time you press the button.
	Invert Selected Area	Inverts the selection.

Deselection



Selection/Deselection

When on, deselects the area using various tools.



Clear All Selection

Clears all selected areas.




Delete Selected Area

Deletes the data from selected area.


How to edit data

Using various tools to select the area you would like to edit.

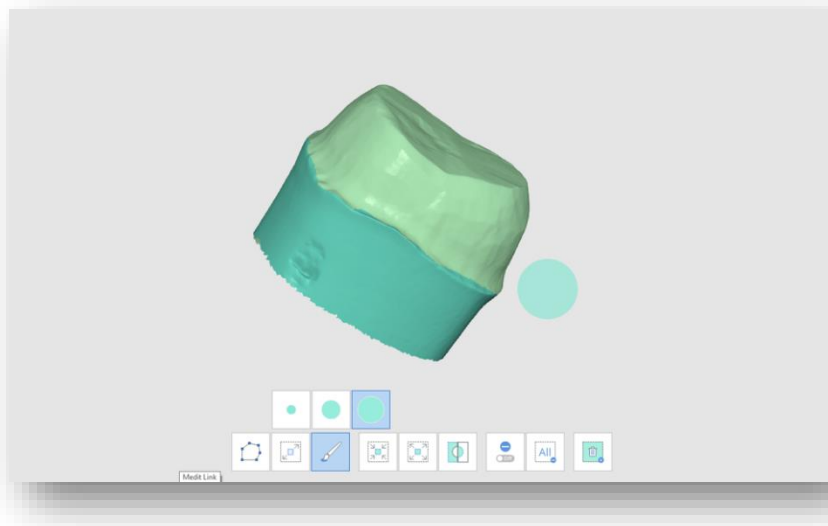


“Polyline” selects all entities within a polyline shape drawn on the screen, while “Brush” selects only the front face data.

① Revert the selection by pressing the “Invert Selected Area” button.

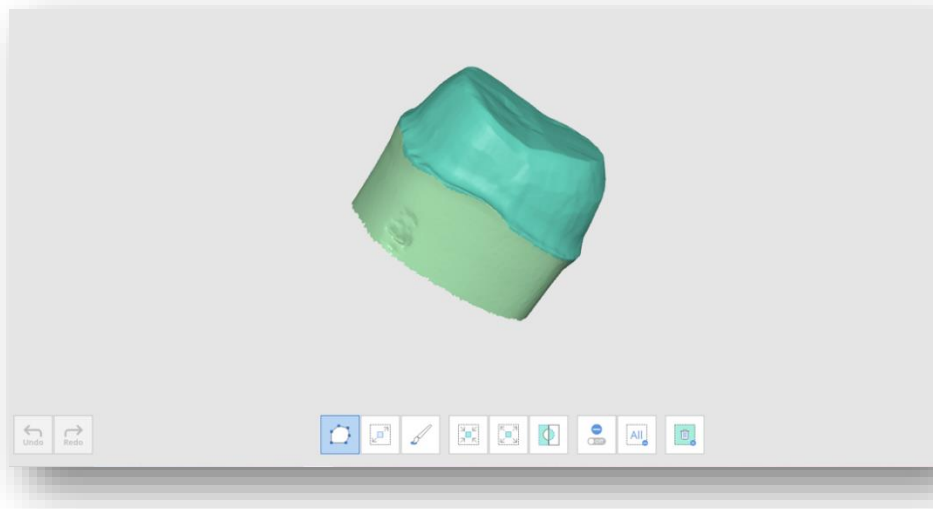


- Reduce/expand the selected areas using the “Shrink Selected Area” / “Expand Selected Area” tools.

Shrink the selection

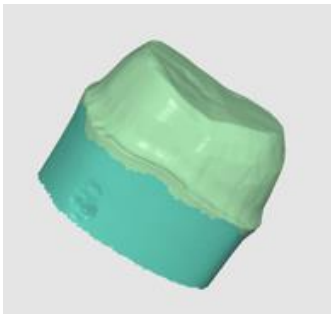

Expand the selection

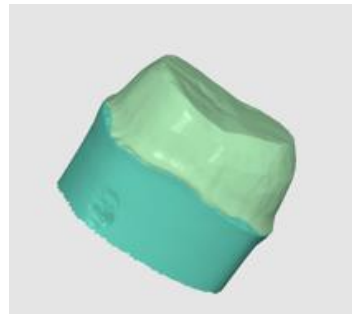

② Click the “Delete Selected Area” button.

6.10.7 Manage the Multi Occlusion Data in File Viewer


You can get multiple results of occlusion data from Medit Scan Application.

Medit Link provides the function to manage the results through the “File Viewer.”


How to manage multiple results

When the user gets multiple results of occlusion from Medit Scan for Clinics, Medit Link provides the option to select a result.

The list of results is listed by name, and you can select one item. The selected model will be shown in the “Data Tree” and “Model Display Area.”

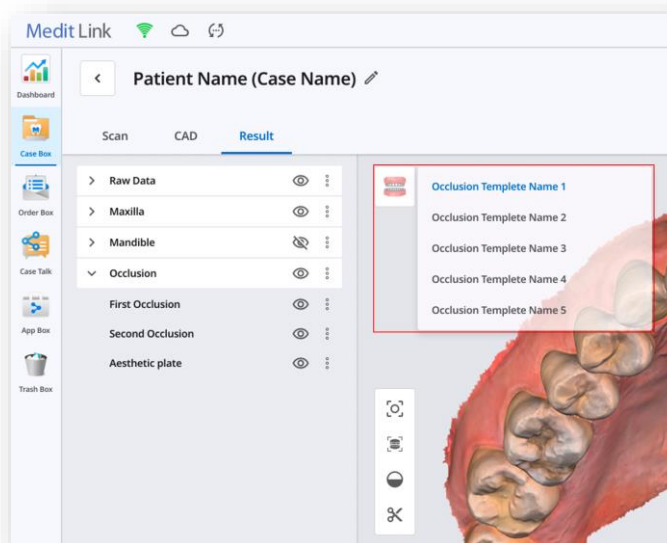

You can collapse the list and change the item using

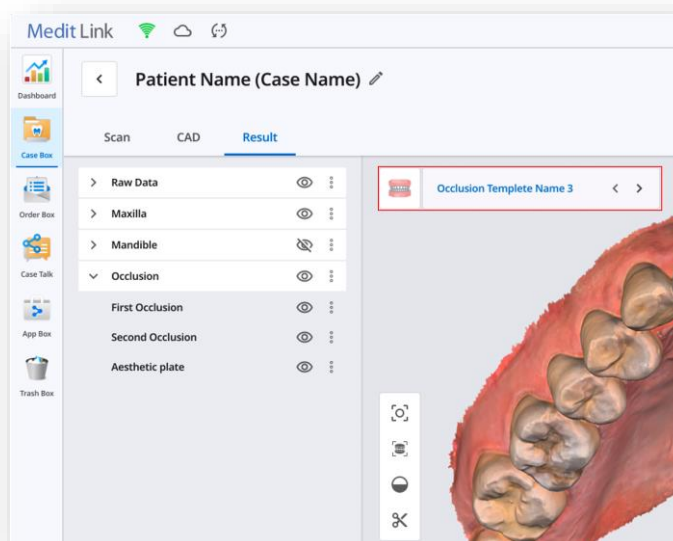


Utilizing the multi occlusion result

You can use the multi occlusion result from your CAD system. In this section, Introduce the simple work process to utilize the multi occlusion result with exocad.

① Acquire the multiple occlusion result from scan app and move to CAD Tab.


② Input the information for CAD work then click  button on the right top.

③ Select the occlusion result which you will use mainly for CAD work then click “Confirm.”

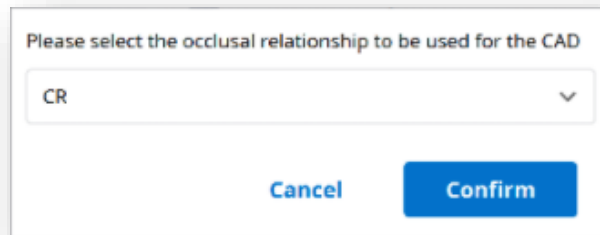


④ The selected occlusion group shown in the exocad.


⑤ Change to “Wizard Mode,” Click “Tools > Add/Remove Mesh.”

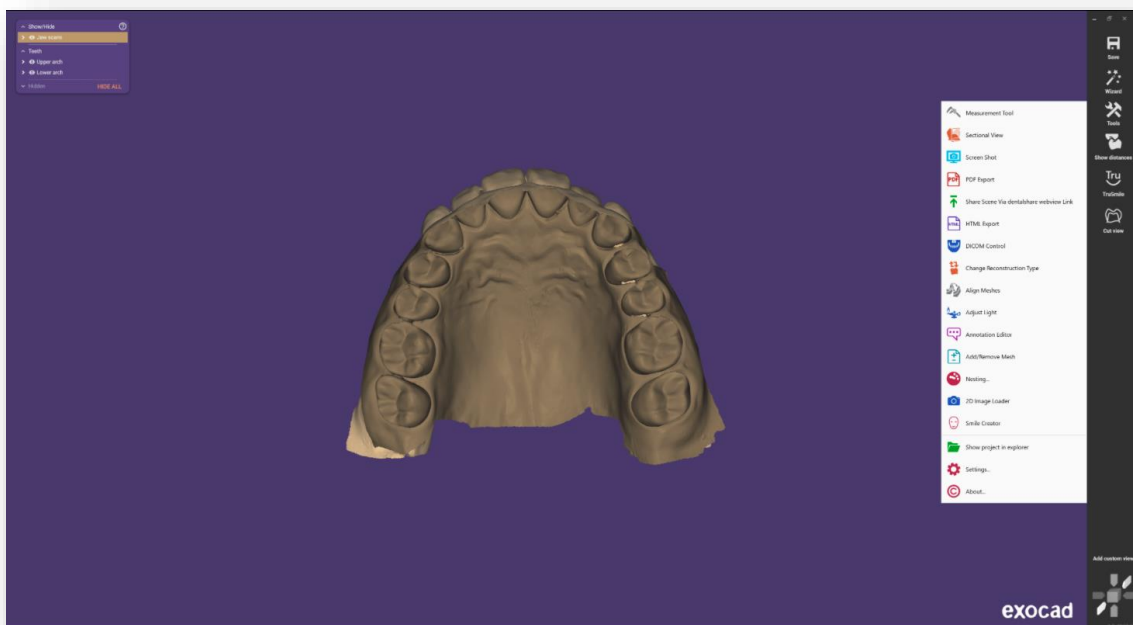

- ⑥ Select the mesh type to “Generic visualization mesh” and click “Load.”

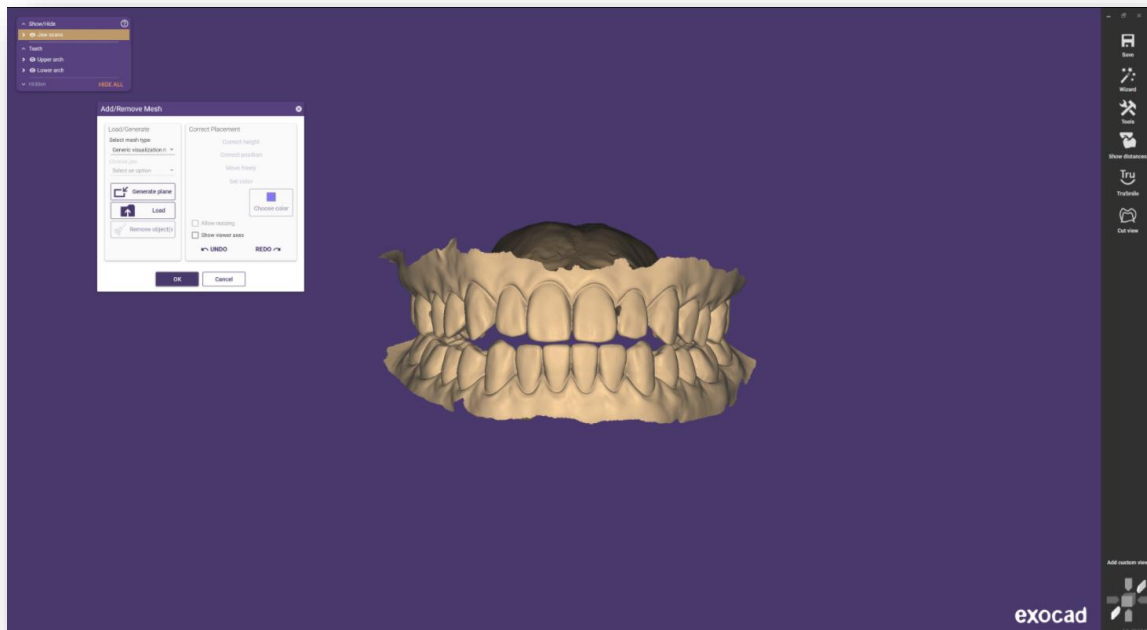

- ⑦ Select the additional occlusion data to reference for CAD work from the working folder.

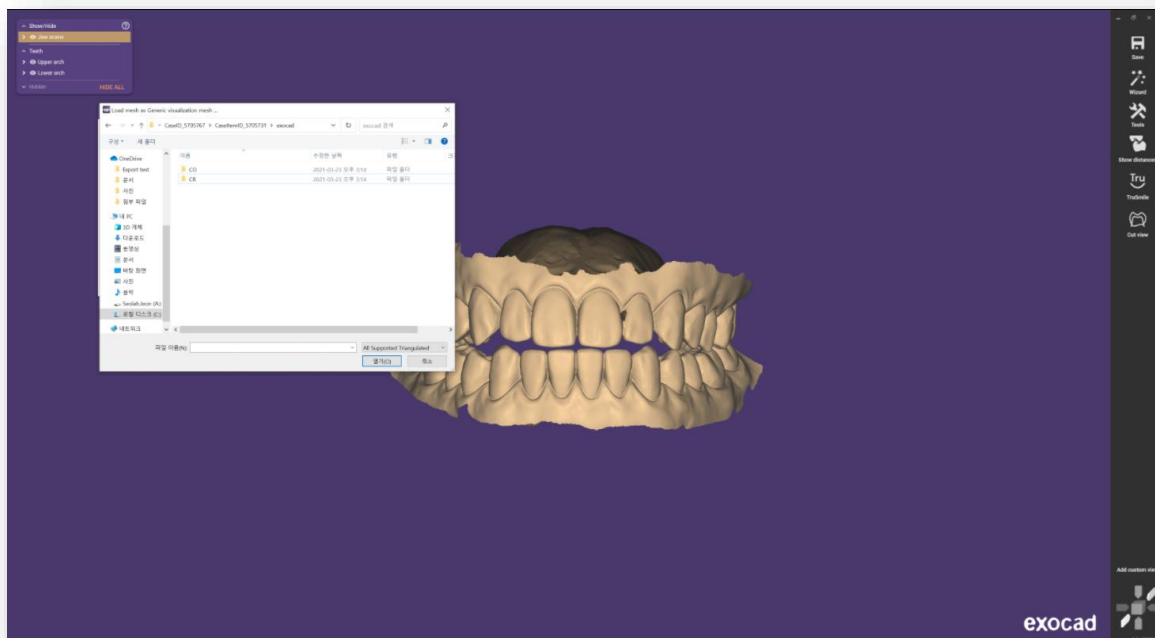

- ⑧ Proceed the CAD work by referencing imported data.

6.11 Order Process and Order Management

You can easily put both acquired and attached data into an order to your partner.

6.11.1 Ordering

You can order cases to your partners from Case Box.

- ① Click “Order” button.
- ② Check the information required for your order.

Order

Partner*
Wonhoonchoi_Lab_Dev Go to Partner

Payment Details **Files for Ordering**


No Price Information

Case Name*
Multi Occlusion Sample

Patient Name*
test

☒ Share Patient Name

Requested Delivery Date
2021-03-25 8:15

Memo

☒ Order All Occlusal Relationship

Cancel OK


- ③ Select the partner's account.
 - Press “Go to Partner” to view information about selected partners.
- ④ Check the product information and files.
 - You can decide whether to share patient name.
- ⑤ Enter the requested delivery date.
 - The default value is set to 7 days after the order date.
- ⑥ You can memos to your partner to outline some specific details.

- ⑦ When the multiple occlusal results generated from the scan app, you can send all occlusal results at once by checking “Order All Occlusal Relationship.”
- ⑧ Proceed with your order.
- **Skip payment:** Proceed with the order without using the payment system provided by Medit Link. Available when ordering from a partner with the “Contracted” status.
 - **Pay:** Utilize the payment system provided by Medit Link to proceed with the order. Available when ordering from a contracted partner.
 - **OK:** We proceed with the order without any conditions regarding the price of the product. Available when ordering from a partner with the “Partnership” status.

Confirm Order – In Box

① Confirm Order


- Accept the order upon receiving it.
- Accepted orders will be forwarded to the Work Box.

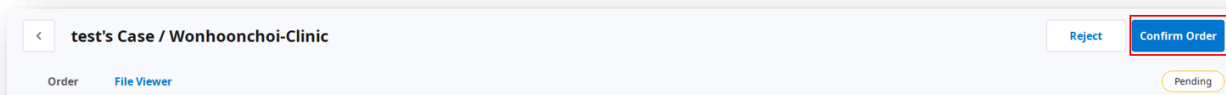

Cancel/Reject Order – Order Box/In Box


① Cancelling Order – Order Box

- You can cancel an order before your partner accepts it.
- Canceled cases can be re-ordered.

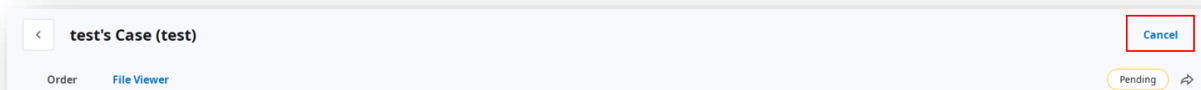


② Reject order – In Box

- You can reject a received order.
- Rejected cases can be re-ordered.

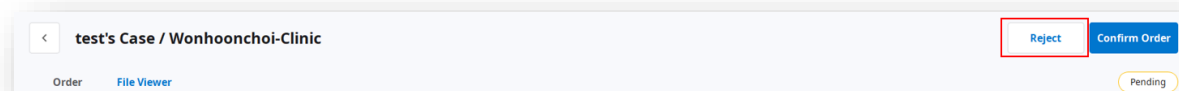

6.11.2 Job Completed – Work Box

Once all tasks are done, press “Done” to complete the case work.

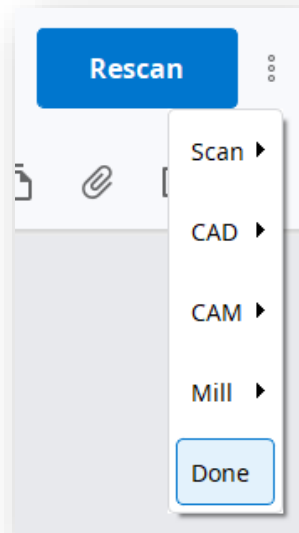

① Enter shipping information – In Box

- Enter the shipping information for the completed task and begin shipping.

② Enter shipping information

- Click on the “Shipment” tab.

MEDIT Link
Worhoochoi_Lab_Dev

2.3.6_Medit Scan for Clinics_Basic / Wonhoochoi_Clinic_Dev

Order
File Viewer
Ready

Order ID
3706292

Patient Name
2.3.6 Data

Partner Name
Wonhoochoi_Clinic_Dev

Date Scanned

Date Ordered
3/11/2021 8:39 AM

Requested Delivery Date
2021-03-18
AM 08:35


Expected Delivery Date
2021-03-18
AM 08:35

Memo
<Scanbody Library>
36 : nt-trading Preform L-Series 1.3 - RC - default
18 : nt-trading Preform L-Series 1.3 - RC - default
14 : DES - 3i OS Ti Base - Desktop Scan Abutment - RP - 1.5mm - Engaging

Order List
Shipment

No.	Type	Method	Material	Shade	Price
16	Implant Crown	Screw Type	Zirconia	A1	0.00
14	Custom Abutment	-	Zirconia	-	0.00
12	Veneer	Laminate	Press Ceramic	A1	0.00
11	Veneer	Laminate	Press Ceramic	A1	0.00
23	Custom Abutment	-	Zirconia	-	0.00
25	Coping	PFZ	Zirconia	A1	0.00
26	Pontic	Reduced	Zirconia	A1	0.00
27	Coping	PFZ	Zirconia	A1	0.00
36	Implant Crown	Screw Type	Zirconia	A1	0.00
35	Custom Abutment	-	Zirconia	-	0.00
46	Inlay	-	Zirconia	A1	0.00
Total					USD 0.00

③ Enter shipping information:

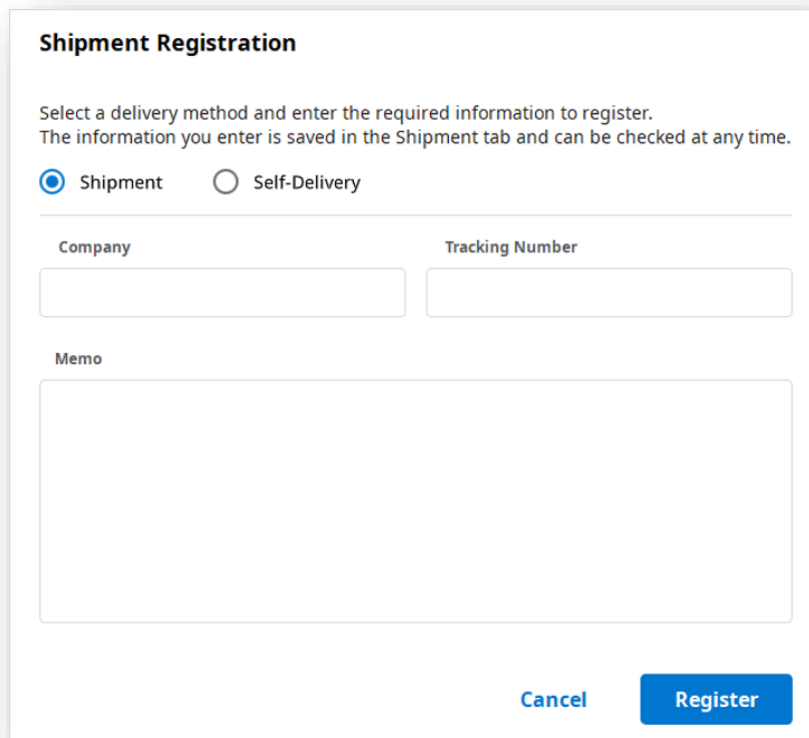

Shipment Registration

Select a delivery method and enter the required information to register.
The information you enter is saved in the Shipment tab and can be checked at any time.

☒ Shipment ☐ Self-Delivery

Company Tracking Number

Memo


Cancel Register

④ Click the “Register” button to complete the shipment information entry.

- The shipping information you entered will be forwarded to your partner.

⑤ Confirmation and Completion of Shipping – Order Box

- After receiving the product proceed with the delivery confirmation.
- The status of the case changes to Completed and the entire order process is completed.

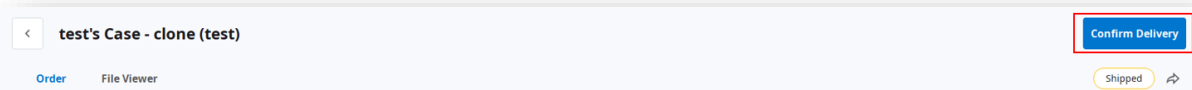

< test's Case - clone (test)

Order File Viewer

Confirm Delivery

Shipped ➔

6.11.3 Checking Order information

You can check the details of the case you ordered from Order Box in clinic account. You can also check the case you received from the In Box in lab account.

2.3.6_Medit Scan for Clinics_Basic / Wonhoonchoi_Clinic_Dev

Order

Order ID: 3706292

Patient Name: 2.3.6 Data

Partner Name: Wonhoonchoi_Clinic_Dev

Date Scanned:

Date Ordered: 3/11/2021 8:39 AM

Requested Delivery Date: 2021-03-18

Expected Delivery Date: 2021-03-18

AM 08:35

AM 08:35

Memo

<Scanbody Library>
36 : nt-trading Preform L-Series 1.3 - RC - default
16 : nt-trading Preform L-Series 1.3 - RC - default
14 : DESS - 3i OS TI Base - Desktop Scan Abutment - RP - 1.5mm - Engaging

Order List

No.	Type	Method	Material	Shade	Price
16	Implant Crown	Screw Type	Zirconia	A1	0.00
14	Custom Abutment	-	Zirconia	-	0.00
12	Veneer	Laminate	Press Ceramic	A1	0.00
11	Veneer	Laminate	Press Ceramic	A1	0.00
23	Custom Abutment	-	Zirconia	-	0.00
25	Coping	PFZ	Zirconia	A1	0.00
26	Pontic	Reduced	Zirconia	A1	0.00
27	Coping	PFZ	Zirconia	A1	0.00
36	Implant Crown	Screw Type	Zirconia	A1	0.00
35	Custom Abutment	-	Zirconia	-	0.00
46	Inlay	-	Zirconia	A1	0.00

Total USD 0.00

① Order information

The information that can be included in the order is as following: Order ID/patient name/account name/scanning date/order date/delivery date/delivery due date/memo

You can check the details of the ordered product and verifying the delivery information.

② Order List

Order list displays the details of the ordered product.

③ Shipping

Enter/check the shipping information.

6.11.4 Managing Expected/Requested Delivery Dates

You can enter the expected delivery date when placing an order, and this information will be forwarded to your partner.

① Enter requested delivery date – Clinic Account

- When ordering from Case Box, you can enter the requested delivery date.

The screenshot shows a web form titled "Order Case". At the top, there is a "Partner*" dropdown menu with "Wonhoonchoi_Lab_Dev" selected and a "Go To Partner" button. Below this are two tabs: "Payment Details" (active) and "File For Order". The main content area has a graphic of three overlapping document icons and the text "No Price Information". Below this are input fields for "Case Name*" (containing "test's Case") and "Patient Name*" (containing "test"). A checkbox labeled "Share Patient Name" is checked. Below the checkbox is the "Delivery Desired Day" section, which is highlighted with a red box. It contains a date picker showing "2020-09-08" and a time dropdown showing "오전 12:00". At the bottom of the form is a "Memo" text area and two buttons: "OK" and "Close".

- The default value is 7 days after the date of the order.
- Enter the date/time you wish to receive the products and proceed with the order.

② Enter expected delivery date – Lab Account

- You can enter the expected delivery date for the case ordered from In Box.

MEDIT Link

test's Case - Clone / Wonhoonchoi_Clinic_Dev

Order ID: 3706004

Patient Name: test

Partner Name: Wonhoonchoi_Clinic_Dev

Date Scanned: 3/23/2021 6:33 PM

Date Ordered: 3/23/2021 6:33 PM

Requested Delivery Date: 2021-03-30 PM 06:30

Expected Delivery Date: 2021-03-30 PM 06:30

Memo

No.	Type	Method	Material	Shade	Price
23	Inlay	-	Zirconia	None Shade	0.00
24	Inlay	-	Zirconia	None Shade	0.00
25	Inlay	-	Zirconia	None Shade	0.00
26	Inlay	-	Zirconia	None Shade	0.00

Total USD 0.00

③ Check the requested delivery date with the partner and make adjustments using .

④ Enter the expected delivery date considering the actual working day. If required, you can leave a comment also.

⑤ Click "Apply."

March2021>

S	M	T	W	T	F	S
28	1	2	3	4	5	6
7	8	9	10	11	12	13
14	15	16	17	18	19	20
21	22	23	24	25	26	27
28	29	30	31	1	2	3
4	5	6	7	8	9	10

Delivery Time

🕒

6:30 PM

Comment

You can leave the comment to your partner

Cancel

Apply

The expected delivery date and memo will be forwarded together with the order.

- You can check the date and memos in the Order Box.
- You can change the expected delivery date, and the changes will be forwarded to your partner.

- You can check the case history using  .

6.12 Patient Management

Patient management tool provides functions for adding/editing/merging patients.

Patient Management

Allows to select patients to delete or merge their information from different cases.
Allows to merge together information of a single patient spread across multiple cases. After merging, all cases of the patient are combined together under the patient name you selected.
Press 'Delete' to remove patients from the system.

All Search / Contains + New Patient

Delete Merge Patients Patients with No Cases

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Date Registered	Recent Visit Date	
☒	5-9_Scan_CAMOrder_Teeth2				3	2020-09-18	2020-09-18	
☒	5-9_Scan_CAMOrder_Teeth2				1	2020-09-18	2020-09-18	
☒	test				2	2020-09-17	2020-09-17	
☒	test				22	2019-07-26	2020-09-17	
☒	test				1	2020-09-17	2020-09-17	
☒	ASD Demo				1	2020-09-16	2020-09-16	
☒	ASD Demo				1	2020-09-16	2020-09-16	
☐	gitt				1	2020-09-15	2020-09-15	
☐	Medit Link 2.3.0			2020-9-24	2	2020-09-14	2020-09-14	
☐	test				1	2020-09-13	2020-09-13	
☐	목/치아스캔				2	2020-09-13	2020-09-13	
☐	test3				1	2020-09-13	2020-09-13	
☐	test1				1	2020-09-13	2020-09-13	
☐	test1231				26	2019-07-25	2020-09-11	

Close

6.12.1 Adding Patients

- ① Click “New Patient.”

New Patient

Patient Name*

Add new patient

General

☐ Female ☐ Male

Birth Date

yyyy mm dd

Patient ID

Phone Number

Email Address

Memo

Creating patient information involves consent from the patient.
Patient information can be shared with 3-parties on order.

[Terms & Conditions](#)

Cancel Save

- ② Enter the required information in the pop-up window, then click “Save” to save the information.

6.12.2 Deleting Patients

- ① Select the patient you want to delete from the list.
 - Patient deletion is possible if there no cases to the patient.

Patient Management

Allows to select patients to delete or merge their information from different cases.
Allows to merge together information of a single patient spread across multiple cases. After merging, all cases of the patient are combined together under the patient name you selected.
Press Delete to remove patients from the system.

All Search / Contains + New Patient

Delete Merge Patients Patients with No Cases

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Date Registered	Recent Visit Date		
☐	test				0	2019-08-28	2019-09-02		
☐	test				0	2019-08-28	2019-08-28		
☐	test				0	2019-08-28	2019-08-28		
☐	A.I Abutment2				0	2019-08-28	2019-08-28		
☐	Scanbody				0	2019-08-28	2019-08-26		
☐	Margin Line Creation				0	2019-08-28	2019-08-22		
☐	Abutment Library				0	2019-08-28	2019-08-22		
☐	GuideMessage				0	2019-08-28	2019-08-16		
☒	margin				0	2019-08-28	2019-08-07		
☒	test2				0	2019-08-28	2019-08-07		
☒	saligadiskf				0	2019-08-28	2019-08-07		
☒	Test				0	2019-08-28	2019-08-07		
☒	test				0	2019-08-28	2019-07-26		
☒	test				0	2019-08-28	2019-07-25		

Close

- ② Click "Delete."

6.12.3 Merging Patients

Manage all cases belonging to one patient all together by merging patient information.

- ① Select the patient you want to merge from the list.

Patient Management

Allows to select patients to delete or merge their information from different cases.
Allows to merge together information of a single patient spread across multiple cases. After merging, all cases of the patient are combined together under the patient name you selected.
Press Delete to remove patients from the system.

Buttons: Delete, Merge Patients, New Patient

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Date Registered	Recent Visit Date	
☐	5-9_iScan_CAMOrder_test12				3	2020-09-18	2020-09-18	
☐	5-9_iScan_CAMOrder_test12				1	2020-09-18	2020-09-18	
☒	test				2	2020-09-17	2020-09-17	
☒	test				22	2019-07-26	2020-09-17	
☒	test				1	2020-09-17	2020-09-17	
☒	ASD Demo				1	2020-09-16	2020-09-16	
☒	ASD Demo				1	2020-09-16	2020-09-16	
☒	gfff				1	2020-09-15	2020-09-15	
☐	Mredit Link 2.3.0			2020-09-24	2	2020-09-14	2020-09-14	
☐	test				1	2020-09-13	2020-09-13	
☐	测试测试				2	2020-09-13	2020-09-13	
☐	test12				1	2020-09-13	2020-09-13	
☐	test1				1	2020-09-13	2020-09-13	
☐	test1231				26	2019-07-25	2020-09-11	

Buttons: Close

- ② Click “Merge Patients.”
- ③ Select the patient you would like to merge all cases to.
- ④ All cases are merged based on the selected information.

Merge all patient information to

Please select the patient you'd like to merge all data to.
After merging all data together, all information in other patient files will be deleted.

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Date Registered	Recent Visit Date
☒	5-9_iScan_CAMOrd...				1	2020-09-18	2020-09-18
☐	test				2	2020-09-17	2020-09-17
☐	test				22	2019-07-26	2020-09-17

Buttons: Cancel, Apply


- ⑤ Click “Apply” to proceed with the merge.

6.13 Case filtering

Use various filtering options to quickly identify the desired patient and case.

6.13.1 Filtering by Status

Provides the ability to filter cases according to their status.

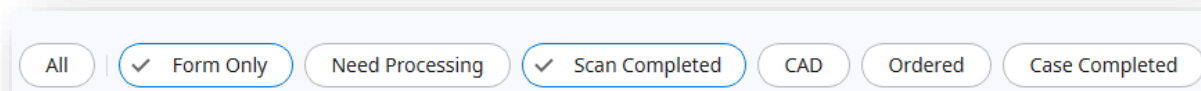

- ① You can click the status you want to view to display only cases in that state.
- ② You can select multiple states at a time.

Different filtering options for each box:


Case Box	All/Form Only/Need Processing/Scan Completed/CAD/Ordered/Case Completed
Order Box	All/Pending/Canceled/Rejected/Accepted/Shipped/Case Completed
In Box	All/Pending/Canceled/Rejected/Accepted/Ready/Shipped/Case Completed
Work Box	All/Form Only/Need Processing/Scan Completed/CAD/Done
Trash Box	Not supported

6.13.2 Filtering Using Period Settings

Filters and displays cases within the period of time you set.

① Click  2020-08-10 ~ 2020-09-09 .

② Select the date range you want to view.

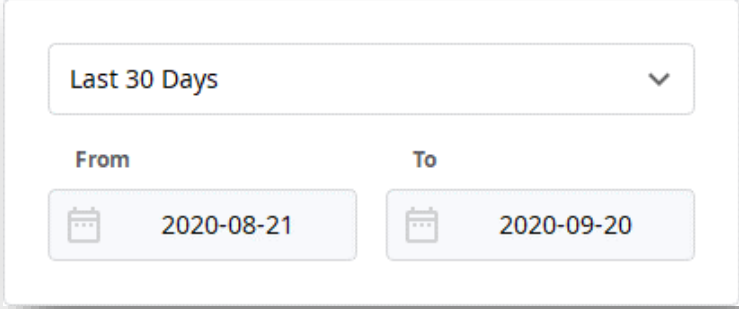


The screenshot shows a date range selection interface. At the top, there is a dropdown menu with the text "Last 30 Days" and a downward arrow. Below this, there are two date pickers. The first is labeled "From" and shows the date "2020-08-21". The second is labeled "To" and shows the date "2020-09-20". Each date picker has a calendar icon to its left.

- All: Shows all cases.
- Custom Range: Specifies the start and end dates to display cases within this period.
- Today: Shows the cases for today's date.
- Yesterday: Shows the cases with yesterday's date.
- Past 7 days: Shows cases before 7 days except today.
- Past 30 days: Shows cases before 30 days except today.
- This Month: Show this month's cases.

The case corresponding to the date range chosen will be updated in real time.

6.13.3 Grouping Cases

Use various grouping and sorting methods for your cases.

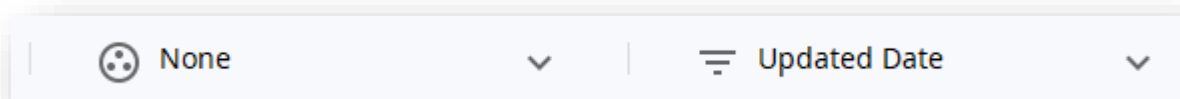

- Grouping Method: Displays the entire case grouped by the selected method.
- Sort by: Displays cases within a group sorted by the selected method.

For each box, the following grouping and sorting methods can be applied:


Box	Grouping	Sorting
Case Box	None	Patient Name / Last Modified Date / Scanning Date/ Order Date / Lab Name
	Patient Name	Last Modified Date / Scanning Date/ Order Date / Lab Name
	Scanning Date	Patient Name / Last Modified Date / Order Date / Lab Name
	Order Date	Patient Name / Last Modified Date / Scanning Date/ Lab Name
Order Box	None	Patient Name / Scanning Date / Order Date/ Requested Delivery Date / Lab Name
	Scanning Date	Patient Name / Order Date / Requested Delivery Date / Lab Name
	Order Date	Patient Name / Scanning Date / Requested Delivery Date / Lab Name
	Requested Delivery Date	Patient Name / Scanning Date / Order Date / Lab Name
In Box	None	Patient Name / Order Date / Requested Delivery Date / Clinic Name
	Order Date	Patient Name / Scanning Date / Requested Delivery Date / Clinic Name
	Requested Delivery Date	Patient Name/ Order Date / Clinic Name
Work Box	None	Patient Name / Last Modified Date / Requested Delivery Date / Clinic Name
	Requested Delivery Date	Patient Name / Last Modified Date / Clinic Name



The Trash Box and Case Talk does not provide a grouping and sorting function.

6.13.4 Filtering with Search Capabilities

Leverage search capabilities to quickly display the information you want.

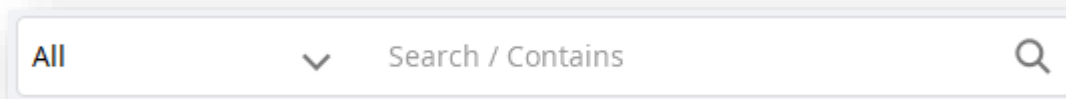

Each box provides the following search options:

Case Box	All, Case Name, Patient Name, Patient ID, Form Information, Lab Name / Contains, Equal
Order Box	All, Case Name, Patient Name, Patient ID, Lab Name, Order ID / Contains, Equal
In Box	All, Case Name, Patient Name, Clinic Name, Order ID / Contains, Equal
Work Box	All, Case Name, Patient Name, Form Information, Clinic Name, Order ID / Contains, Equal
Case Talk	Case Talk Name (Contains), Patient Name (Equals)
Trash Box – Clinic	All, Case Name, Patient Name, Lab Name, Deleted By / Contains, Equal
Trash Box – Lab	All, Case Name, Patient Name, Clinic Name, Deleted By / Contains, Equal